

Alexa Fluor® 488 anti-human CD27 Antibody

Catalog# / Size	302871 / 25 tests 302872 / 100 tests
Clone	O323
Regulatory Status	RUO
Workshop	IV T-186
Other Names	S152, T14, TNFRSF7
Isotype	Mouse IgG1, κ
Description	CD27 is a 50-55 kD type I membrane protein also known as S152 and T14. It is a lymphocyte-specific member of the TNF-receptor superfamily. CD27 is expressed on medullary thymocytes, virtually all mature T cells, some B cells, and NK cells. CD27 binds to CD70 and plays an important role in costimulation of T cell activation, and regulation of B cell differentiation and proliferation. The cytoplasmic domains of CD27 have also been shown to interact with TRAF2 and TRAF5 to elicit NF-κB and SAPK/JNK activation.

Product Details

Verified Reactivity	Human, Cynomolgus, Rhesus
Reported Reactivity	African Green, Baboon, Squirrel Monkey
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Alexa Fluor® 488 under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * Alexa Fluor® 488 has a maximum emission of 519 nm when it is excited at 488 nm. Alexa Fluor® and Pacific Blue™ are trademarks of Life Technologies Corporation. View full statement regarding label licenses
Excitation Laser	Blue Laser (488 nm)
Application References	1. Knapp W, <i>et al.</i> Eds. 1989. Leucocyte Typing IV. Oxford University Press. New York. 2. Correia DV, <i>et al.</i> 2011. <i>Blood</i> 118:992. (FC) PubMed
(PubMed link indicates BioLegend citation)	
RRID	AB_3683308 (BioLegend Cat. No. 302871) AB_3683308 (BioLegend Cat. No. 302872)

Antigen Details

Structure	TNF-R superfamily, type I transmembrane glycoprotein, 50-55 kD
Distribution	Medullary thymocytes, T and B cell subsets, NK cells
Function	T cell costimulation
Ligand/Receptor	CD70
Cell Type	B cells, NK cells, T cells, Thymocytes, Tregs
Biology Area	Costimulatory Molecules, Immunology
Molecular Family	CD Molecules
Antigen References	1. Hintzen R, <i>et al.</i> 1994. <i>Immunol. Today</i> 15:307. 2. Agematsu K, <i>et al.</i> 1995. <i>J. Immunol.</i> 154:3627.
Gene ID	939


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-human CD27, Biotin anti-human CD27, FITC anti-human CD27, PE anti-human CD27, Purified anti-human CD27, Alexa Fluor® 647 anti-human CD27, Alexa Fluor® 700 anti-human CD27, APC/Cyanine7 anti-human CD27, PerCP anti-human CD27, PerCP/Cyanine5.5 anti-human CD27, Pacific Blue™ anti-human CD27, Brilliant Violet 421™ anti-human CD27, Brilliant Violet 570™ anti-human CD27, Brilliant Violet 650™ anti-human CD27, Brilliant Violet 605™ anti-human CD27, Brilliant Violet 711™ anti-human CD27, Brilliant Violet 785™ anti-human CD27, Brilliant Violet 510™ anti-human CD27, PE/Cyanine7 anti-human CD27, Purified anti-human CD27 (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD27, APC/Fire™ 750 anti-human CD27, TotalSeq™-A0154 anti-human CD27, Brilliant Violet 750™ anti-human CD27, TotalSeq™-B0154 anti-human CD27, TotalSeq™-C0154 anti-human CD27, Spark NIR™ 685 anti-human CD27, PE/Fire™ 810 anti-human CD27, TotalSeq™-D0154 anti-human CD27, APC/Fire™ 810 anti-human CD27, PE/Cyanine5 anti-human CD27 Antibody, Spark UV™ 387 anti-human CD27, Spark Blue™ 515 anti-human CD27, Spark Red™ 718 anti-human CD27, Alexa Fluor® 488 anti-human CD27, PerCP/Fire™ 806 anti-human CD27 Antibody, Spark YG™ 581 anti-human CD27 (Flexi-Fluor™), Spark YG™ 593 anti-human CD27 (Flexi-Fluor™) Antibody

Product Data

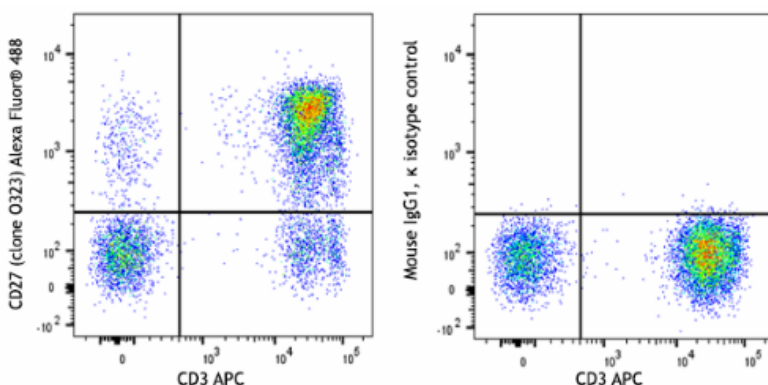

Human peripheral blood lymphocytes were stained with anti-human CD3 (clone UCHT1) APC and anti-human CD27 (clone O323) Alexa Fluor® 488 (left) or mouse IgG1, κ Alexa Fluor® 488 isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com

Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587