

Alexa Fluor® 555 Mouse IgG2a, κ Isotype Ctrl Antibody

Catalog# / Size	402543 / 100 µg
Clone	MOPC-173
Regulatory Status	RUO
Isotype	Mouse IgG2a, κ
Description	The MOPC-173 immunoglobulin has unknown specificity. The isotype of this antibody is mouse IgG2a, κ. This antibody was chosen as an isotype control after screening on a variety of resting, activated, live, and fixed mouse, rat, and human tissues.

Product Details

Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Preparation	The antibody was purified by affinity chromatography and conjugated with Alexa Fluor® 555 under optimal conditions.
Concentration	0.5 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	ICC - Quality tested IHC-P - Verified
Recommended Usage	Each lot of this antibody is quality control tested by immunocytochemistry (ICC). For ICC, a concentration range of 2.5 - 10 µg/mL is recommended. For immunohistochemistry on formalin-fixed paraffin-embedded tissues (IHC-P), a concentration range of 5.0 -10.0 µg/mL is recommended. It is recommended that the reagent be titrated for optimal performance in your experimental conditions. * Alexa Fluor® 555 has an excitation maximum of 555 nm, and a maximum emission of 565 nm. Alexa Fluor® and Pacific Blue™ are trademarks of Life Technologies Corporation. View full statement regarding label licenses
Excitation Laser	Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm)

Antigen Details

Gene ID	NA
----------------	----

Related Protocols


- [Immunocytochemistry Staining Protocol](#)
- [Immunohistochemistry Protocol for Paraffin-Embedded Sections](#)

Spark YG™ 581 Mouse IgG2a, κ Isotype Ctrl, Spark Violet™ 538 Mouse IgG2a, κ Isotype Ctrl, TotalSeq™-D0091 Mouse IgG2a, κ Isotype Ctrl, Spark YG™ 593 Mouse IgG2a, κ Isotype Ctrl, GMP PE Mouse IgG2a, κ Isotype Ctrl, GMP FITC Mouse IgG2a, κ Isotype


Other Formats

Ctrl, GMP APC Mouse IgG2a, κ Isotype Ctrl, GMP PE/Dazzle™ 594 Mouse IgG2a, κ Isotype Ctrl, GMP PerCP/Cyanine5.5 Mouse IgG2a, κ Isotype Ctrl, GMP APC/Fire™ 750 Mouse IgG2a, κ Isotype Ctrl, PE/Fire™ 700 Mouse IgG2a, κ Isotype Ctrl, GMP Pacific Blue™ Mouse IgG2a, κ Isotype Ctrl, Spark YG™ 570 Mouse IgG2a, κ Isotype Ctrl, GMP PE/Cyanine7 Mouse IgG2a, κ Isotype Ctrl, PE/Fire™ 640 Mouse IgG2a, κ Isotype Ctrl, PerCP/Fire™ 806 Mouse IgG2a, κ Isotype Ctrl, PerCP/Fire™ 780 Mouse IgG2a, κ Isotype Ctrl, Spark Violet™ 423 Mouse IgG2a, κ Isotype Ctrl, APC/Fire™ 810 Mouse IgG2a, κ Isotype Ctrl, Spark Blue™ 515 Mouse IgG2a, κ Isotype Ctrl, Spark Blue™ 574 Mouse IgG2a, κ Isotype Ctrl, Alexa Fluor® 660 Mouse IgG2a, κ Isotype Ctrl, Spark UV™ 387 Mouse IgG2a, κ Isotype Ctrl, Alexa Fluor® 594 Mouse IgG2a, κ Isotype Ctrl, APC/Fire™ 750 Mouse IgG2a, κ Isotype Ctrl, PE Mouse IgG2a, κ Isotype Ctrl, Purified Mouse IgG2a, κ Isotype Ctrl, PE/Dazzle™ 594 Mouse IgG2a, κ Isotype Ctrl, PE/Cyanine5 Mouse IgG2a, κ Isotype Ctrl, APC Mouse IgG2a, κ Isotype Ctrl, Brilliant Violet 570™ Mouse IgG2a, κ Isotype Ctrl, Brilliant Violet 421™ Mouse IgG2a, κ Isotype Ctrl, Ultra-LEAF™ Purified Mouse IgG2a, κ Isotype Ctrl, Brilliant Violet 510™ Mouse IgG2a, κ Isotype Ctrl, Brilliant Violet 785™ Mouse IgG2a, κ Isotype Ctrl, PE/Fire™ 810 Mouse IgG2a, κ Isotype Ctrl, Spark PLUS UV395™ Mouse IgG2a, κ Isotype Ctrl, TotalSeq™-B1472 Mouse IgG2a, κ Isotype Ctrl Antibody, Spark PLUS B550™ Mouse IgG2a, κ Isotype Ctrl, Spark Red™ 718 Mouse IgG2a, κ Isotype Ctrl, PE/Fire™ 744 Mouse IgG2a, κ Isotype Ctrl, Alexa Fluor® 488 Mouse IgG2a, κ Isotype Ctrl, PE Mouse IgG2a, κ Isotype Ctrl (ICFC), PerCP Mouse IgG2a, κ Isotype Ctrl, FITC Mouse IgG2a, κ Isotype Ctrl, TotalSeq™-A0091 Mouse IgG2a, κ Isotype Ctrl, PE/Cyanine7 Mouse IgG2a, κ Isotype Ctrl, Pacific Blue™ Mouse IgG2a, κ Isotype Ctrl, Biotin Mouse IgG2a, κ Isotype Ctrl, FITC Mouse IgG2a, κ Isotype Ctrl, Spark Blue™ 550 Mouse IgG2a, κ Isotype Ctrl, APC/Cyanine7 Mouse IgG2a, κ Isotype Ctrl, Brilliant Violet 711™ Mouse IgG2a, κ Isotype Ctrl, Alexa Fluor® 647 Mouse IgG2a, κ Isotype Ctrl, PerCP/Cyanine5.5 Mouse IgG2a, κ Isotype Ctrl, TotalSeq™-B0091 Mouse IgG2a, κ Isotype Ctrl, Brilliant Violet 650™ Mouse IgG2a, κ Isotype Ctrl, Brilliant Violet 605™ Mouse IgG2a, κ Isotype Ctrl, TotalSeq™-C1472 Mouse IgG2a, κ Isotype Ctrl Antibody, Alexa Fluor® 555 Mouse IgG2a, κ Isotype Ctrl, FITC Mouse IgG2a, κ Isotype Ctrl (FC), PE Mouse IgG2a, κ Isotype Ctrl (FC), APC Mouse IgG2a, κ Isotype Ctrl (FC), Alexa Fluor® 700 Mouse IgG2a, κ Isotype Ctrl, TotalSeq™-C0091 Mouse IgG2a, κ Isotype Ctrl, KIRAVIA Blue 520™ Mouse IgG2a, κ Isotype Ctrl, Spark NIR™ 685 Mouse IgG2a, κ Isotype Ctrl

Product Data

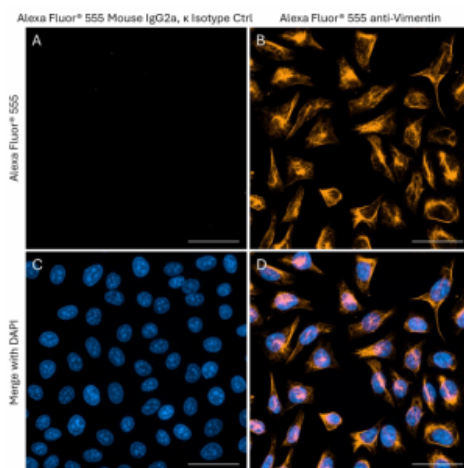

HeLa cells were fixed and permeabilized with 100% ice-cold methanol for 10 minutes, and blocked with 5% FBS for 1 hour at room temperature. Cells were then stained with Alexa Fluor® 555 Mouse IgG2a, κ Isotype Ctrl (clone MOPC-173) (panel A) or Alexa Fluor® 555 anti-Vimentin (clone O91D3) (panel B). Nuclei were counterstained with DAPI (Cat. No. 422801). The images were captured on a Revvity Operetta CLS™ High Content Analysis System with a 63X objective and merged with DAPI (panels C and D). Scale bar: 50 μm

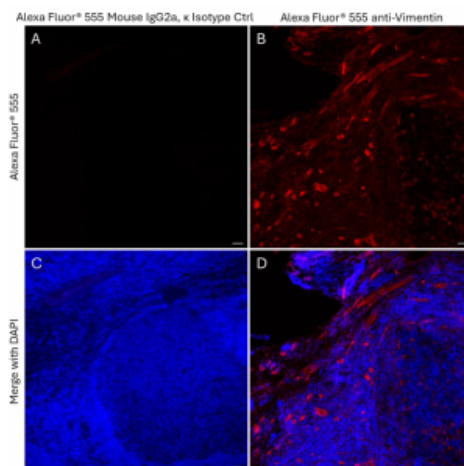

IHC staining on formalin-fixed paraffin-embedded human tonsil tissue. Following antigen retrieval using 1X Tris-EDTA pH 9.0 Antigen Retrieval Buffer (Cat. No. 422704), the tissue was incubated with Alexa Fluor® 555 Mouse IgG2a, κ Isotype Ctrl (clone MOPC-173) (panel A) or Alexa Fluor® 555 anti-Vimentin (clone O91D3) (panel B). Nuclei were counterstained with DAPI (Cat. No. 422801). Images were captured with a 10X objective and merged (panels C and D). Scale bar: 50 μm

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biollegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587