

APC/Fire™ 810 anti-human CD45RA Antibody

Catalog# / Size	304191 / 25 tests 304192 / 100 tests
Clone	HI100
Regulatory Status	RUO
Workshop	IV N906
Other Names	GP180, L-CA, LCA, LY5, T200, PTPRC
Isotype	Mouse IgG2b, κ
Description	CD45RA is a 205-220 kD single chain type I glycoprotein. It is an exon 4 splice variant of the tyrosine phosphatase CD45. The CD45RA isoform is expressed on resting/naïve T cells, medullary thymocytes, B cells and monocytes. CD45RA enhances both T cell receptor and B cell receptor signaling. CD45 non-covalently associates with lymphocyte phosphatase-associated phosphoprotein (LPAP) on T and B lymphocytes. CD45 has been reported to be associated with several other cell surface antigens including CD1, CD2, CD3, and CD4. CD45 has also been reported to bind galectin-1. CD45 isoform expression can change in response to cytokines.

Product Details

Verified Reactivity	Human
Reported Reactivity	Chimpanzee
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with APC/Fire™ 810 under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>* APC/Fire™ 810 has a maximum excitation of 650 nm and a maximum emission of 810 nm.</p> <p>Excessive exposure to light, and commonly used fixation, permeabilization buffers can affect APC/Fire™ 810 fluorescence signal intensity and spread. Please keep conjugates protected from light exposure. For more information and representative data, visit our Fire Dyes page.</p>
Excitation Laser	Red Laser (633 nm)
Application Notes	Additional reported applications (for relevant formats of this clone) include: inhibition of CD45 functions ² , immunohistochemical staining of frozen tissue sections ³ and formalin-fixed paraffin-embedded tissue sections ⁴ , and immunocytochemistry ^{15,16} .

Application References

(PubMed link indicates
BioLegend citation)

1. Knapp W, *et al.* 1989. Leucocyte Typing IV. Oxford University Press. New York.
2. Yamada T, *et al.* 2002. *J. Biol. Chem.* 277:28830. (WB, Block)
3. Weninger W, *et al.* 2003 *J. Immunol.* 170:4638. (IHC-F)
4. Imaniguli MM, *et al.* 2009. *Blood.* 113:3620 (IHC-P)
5. Roque S, *et al.* 2007. *J. Immunol.* 178:8028. (FC) [PubMed](#)
6. Smeltz RB. 2007. *J. Immunol.* 178:4786. (FC) [PubMed](#)
7. Palendira U, *et al.* 2008. *Blood* (FC) [PubMed](#)
8. Kuttruff S, *et al.* 2009. *Blood* 113:358. (FC) [PubMed](#)
9. Thakral D, *et al.* 2008. *J. Immunol.* 180:7431. (FC) [PubMed](#)
10. Alanio C, *et al.* 2010. *Blood* 115:3718. (FC) [PubMed](#)
11. Iannello A, *et al.* 2010. *J. Immunol.* 184:114. (FC) [PubMed](#)
12. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)

[See More](#)

RRID

AB_3675028 (BioLegend Cat. No. 304191)
AB_3675028 (BioLegend Cat. No. 304192)

Antigen Details

Structure	Tyrosine phosphatases, type I transmembrane (exon 4 splicing of CD45 gene), 205-220 kD
Distribution	B cells, naïve T cells, monocytes
Function	Enhances TCR and BCR signaling
Ligand/Receptor	Galectin-1, CD2, CD3, CD4
Cell Type	B cells, Monocytes, T cells, Tregs
Biology Area	Cell Biology, Immunology, Inhibitory Molecules, Neuroscience, Neuroscience Cell Markers
Molecular Family	CD Molecules
Antigen References	1. Thomas M. 1989. <i>Annu. Rev. Immunol.</i> 7:339. 2. Trowbridge I, <i>et al.</i> 1994. <i>Annu. Rev. Immunol.</i> 12:85.
Gene ID	5788


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-human CD45RA, Biotin anti-human CD45RA, FITC anti-human CD45RA, PE anti-human CD45RA, PE/Cyanine5 anti-human CD45RA, Purified anti-human CD45RA, Alexa Fluor® 488 anti-human CD45RA, Alexa Fluor® 647 anti-human CD45RA, Pacific Blue™ anti-human CD45RA, Alexa Fluor® 700 anti-human CD45RA, PerCP/Cyanine5.5 anti-human CD45RA, PE/Cyanine7 anti-human CD45RA, APC/Cyanine7 anti-human CD45RA, Brilliant Violet 421™ anti-human CD45RA, Brilliant Violet 570™ anti-human CD45RA, Brilliant Violet 605™ anti-human CD45RA, Brilliant Violet 650™ anti-human CD45RA, Brilliant Violet 711™ anti-human CD45RA, Brilliant Violet 785™ anti-human CD45RA, Brilliant Violet 510™ anti-human CD45RA, Purified anti-human CD45RA (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD45RA, APC/Fire™ 750 anti-human CD45RA, PerCP anti-human CD45RA, TotalSeq™-A0063 anti-human CD45RA, Alexa Fluor® 594 anti-human CD45RA, TotalSeq™-B0063 anti-human CD45RA, TotalSeq™-C0063 anti-human CD45RA, Brilliant Violet 750™ anti-human CD45RA, Spark NIR™ 685 anti-human CD45RA, PE/Fire™ 640 anti-human CD45RA, PE/Fire™ 700 anti-human CD45RA Antibody, Spark YG™ 581 anti-human CD45RA, TotalSeq™-D0063 anti-human CD45RA, Spark Violet™ 423 anti-human CD45RA, GMP FITC anti-human CD45RA, Spark UV™ 387 anti-human CD45RA, GMP APC anti-human CD45RA, TotalSeq™-Bn0063 anti-human CD45RA, Spark Blue™ 550 anti-human CD45RA, GMP PE/Dazzle™ 594 anti-human CD45RA, GMP APC/Fire™ 750 anti-human CD45RA, Spark PLUS UV395™ anti-human CD45RA, Spark Red™ 718 anti-human CD45RA, APC/Fire™ 810 anti-human CD45RA, PE/Fire™ 810 anti-human CD45RA, Spark Blue™ 574 anti-human CD45RA (Flexi-Fluor™), PerCP/Fire™ 780 anti-human CD45RA Antibody, GMP PE/Cyanine7 anti-human CD45RA, StarBright UltraViolet 795 anti-human CD45RA, Spark YG™ 593 anti-human CD45RA (Flexi-Fluor™) Antibody

Product Data

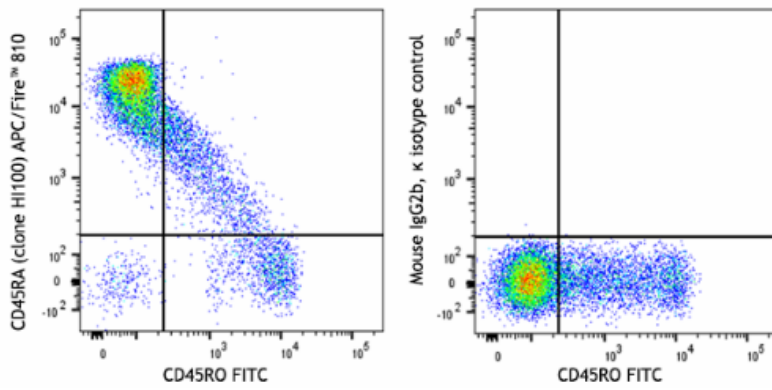

Human peripheral blood lymphocytes were stained with anti-human CD45RO (clone UCHL1) FITC and anti-human CD45RA (clone HI100) APC/Fire™ 810 (left) or mouse IgG2b, κ APC/Fire™ 810 isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587