

APC/Fire™ 810 anti-mouse CD4 Antibody

Catalog# / Size	100581 / 25 µg 100582 / 100 µg
Clone	RM4-5
Regulatory Status	RUO
Other Names	L3T4, T4
Isotype	Rat IgG2a, κ
Description	CD4 is a 55 kD protein also known as L3T4 or T4. It is a member of the Ig superfamily, primarily expressed on most thymocytes and a subset of T cells, and weakly on macrophages and dendritic cells. It acts as a co-receptor with the TCR during T cell activation and thymic differentiation by binding MHC class II and associating with the protein tyrosine kinase Lck.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Rat
Immunogen	BALB/c mouse thymocytes
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Preparation	The antibody was purified by affinity chromatography and conjugated with APC/Fire™ 810 under optimal conditions.
Concentration	0.2 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is ≤ 1.0 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>* APC/Fire™ 810 has a maximum excitation of 650 nm and a maximum emission of 810 nm.</p> <p>Excessive exposure to light, and commonly used fixation, permeabilization buffers can affect APC/Fire™ 810 fluorescence signal intensity and spread. Please keep conjugates protected from light exposure. For more information and representative data, visit our Fire Dyes page.</p>
Excitation Laser	Red Laser (633 nm)
Application Notes	The RM4-5 antibody blocks the binding of GK1.5 antibody and H129.19 antibody to CD4 ⁺ T cells, but not RM4-4 antibody. Additional reported applications (for the relevant formats) include: blocking of ligand binding, <i>in vivo</i> depletion of CD4 ⁺ cells ¹ , and immunohistochemistry of acetone-fixed frozen tissue sections ^{2,3,11} and paraffin-embedded sections ¹¹ . Clone RM4-5 is not recommended for immunohistochemistry of formalin-fixed paraffin sections. Instead, acetone frozen or zinc-fixed paraffin sections are recommended. The Ultra-LEAF™ Purified antibody (Endotoxin < 0.01 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for functional assays (Cat. No. 100575 and 100576).

Application References

(PubMed link indicates
BioLegend citation)

1. Kruisbeek AM. 1991. *In Curr. Protocols Immunol.* pp. 4.1.1-4.1.5. (Block, Deplete)
2. Nitta H, *et al.* 1997. *Cell Vision* 4:73. (IHC)
3. Fan WY, *et al.* 2001. *Exp. Biol. Med.* 226:1045.
4. Muraille E, *et al.* 2003. *Infect. Immun.* 71:2704. (IHC)
5. León-Ponte M, *et al.* 2007. *Blood* 109:3139. (FC)
6. Bourdeau A, *et al.* 2007. *Blood* doi:10.1182/blood-2006-08-044370. (FC)
7. Matsumoto M, *et al.* 2007. *J. Immunol.* 178:2499. [PubMed](#)
8. Shigeta A, *et al.* 2008. *Blood* 112:4915. [PubMed](#)
9. Zaborsky N, *et al.* 2010. *J. Immunol.* 184:725. [PubMed](#)
10. Rodrigues-Manzanet R, *et al.* 2010. *P. Natl Acad Sci USA* 107:8706. [PubMed](#)
11. Whiteland JL, *et al.* 1995. *J. Histochem. Cytochem.* 43:313. (IHC)

RRID AB_3683125 (BioLegend Cat. No. 100581)
AB_3683125 (BioLegend Cat. No. 100582)

Antigen Details

Structure	Ig superfamily, 55 kD
Distribution	Majority of thymocytes, T cell subset
Function	TCR co-receptor, T cell activation
Ligand/Receptor	MHC class II molecule
Cell Type	Dendritic cells, T cells, Thymocytes, Tregs
Biology Area	Immunology
Molecular Family	CD Molecules
Antigen References	<ol style="list-style-type: none">1. Barclay A, <i>et al.</i> 1997. <i>The Leukocyte Antigen FactsBook</i> Academic Press.2. Bierer BE, <i>et al.</i> 1989. <i>Annu. Rev. Immunol.</i> 7:579.3. Janeway CA. 1992. <i>Annu. Rev. Immunol.</i> 10:645.
Gene ID	12504


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-mouse CD4, Biotin anti-mouse CD4, FITC anti-mouse CD4, PE anti-mouse CD4, PE/Cyanine5 anti-mouse CD4, Purified anti-mouse CD4, PE/Cyanine7 anti-mouse CD4, APC/Cyanine7 anti-mouse CD4, Alexa Fluor® 647 anti-mouse CD4, Alexa Fluor® 488 anti-mouse CD4, Pacific Blue™ anti-mouse CD4, Alexa Fluor® 700 anti-mouse CD4, PerCP anti-mouse CD4, PerCP/Cyanine5.5 anti-mouse CD4, Brilliant Violet 421™ anti-mouse CD4, APC/Fire™ 750 anti-mouse CD4, Brilliant Violet 570™ anti-mouse CD4, Brilliant Violet 605™ anti-mouse CD4, Brilliant Violet 650™ anti-mouse CD4, Brilliant Violet 711™ anti-mouse CD4, Brilliant Violet 785™ anti-mouse CD4, Brilliant Violet 510™ anti-mouse CD4, Purified anti-mouse CD4 (Maxpar® Ready), PE/Dazzle™ 594 anti-mouse CD4, TotalSeq™-A0001 anti-mouse CD4, TotalSeq™-B0001 anti-mouse CD4, TotalSeq™-C0001 anti-mouse CD4, Ultra-LEAF™ Purified anti-mouse CD4, Spark Violet™ 423 anti-mouse CD4, Spark Red™ 718 anti-mouse CD4 (Flexi-Fluor™), APC/Fire™ 810 anti-mouse CD4, PE/Fire™ 640 anti-mouse CD4 Antibody, Brilliant Violet 750™ anti-mouse CD4 Antibody, Spark YG™ 581 anti-mouse CD4 (Flexi-Fluor™), StarBright UltraViolet 740 anti-mouse CD4

Product Data

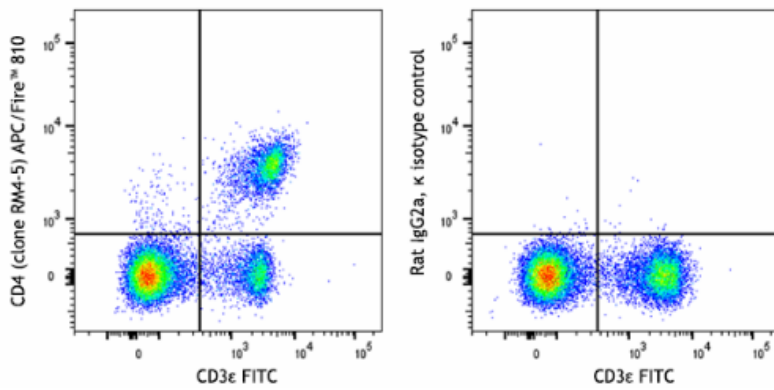

C57BL/6 mouse splenocytes were stained with anti-mouse CD3ε (clone 145-2C11) FITC and anti-mouse CD4 (clone RM4-5) APC/Fire™ 810 (left) or rat IgG2a, κ APC/Fire™ 810 isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587