

## APC/Fire™ 810 anti-mouse F4/80 Recombinant Antibody

<b>Catalog# / Size</b>	157323 / 25 µg 157324 / 100 µg
<b>Clone</b>	QA17A29
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	EMR1, Ly71
<b>Isotype</b>	Mouse IgG1, κ
<b>Description</b>	F4/80 is a 160 kD glycoprotein. It is characterized as a member of the epidermal growth factor (EGF)-transmembrane 7 (TM7) family. F4/80, also known as EMR1 or Ly71, has been widely used as a murine macrophage marker, which is expressed on the majority of tissue macrophages including peritoneal macrophages, macrophages in lung, gut, thymus and red pulp of spleen (but not on the macrophages located in T cell areas of the spleen, lymph node and Peyer's patch), Kuffer cells, Langerhans cells, and bone marrow stromal cells. F4/80 has also been shown on a subset of dendritic cells. The biological ligand of F4/80 has not been identified, but it has been reported that F4/80 is required for induction of CD8 <sup>+</sup> T cells-mediated peripheral tolerance.

### Product Details

---

<b>Verified Reactivity</b>	Mouse
<b>Antibody Type</b>	Recombinant
<b>Host Species</b>	Mouse
<b>Immunogen</b>	Murine macrophages
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with APC/Fire™ 810 under optimal conditions.
<b>Concentration</b>	0.2 mg/mL
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	<p>Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a>. For flow cytometric staining, the suggested use of this reagent is ≤ 0.25 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>* APC/Fire™ 810 has a maximum excitation of 650 nm and a maximum emission of 810 nm.</p> <p>Excessive exposure to light, and commonly used fixation, permeabilization buffers can affect APC/Fire™ 810 fluorescence signal intensity and spread. Please keep conjugates protected from light exposure. For more information and representative data, visit our <a href="#">Fire Dyes</a> page.</p>
<b>Excitation Laser</b>	Red Laser (633 nm)
<b>RRID</b>	AB_3683173 (BioLegend Cat. No. 157323) AB_3683173 (BioLegend Cat. No. 157324)

### Antigen Details

---

<b>Structure</b>	A 160 kD glycoprotein and a member of the epidermal growth factor (EGF)-transmembrane 7 (TM7) family
<b>Distribution</b>	Expressed on the majority of tissue macrophages
<b>Function</b>	Required for induction of CD8 <sup>+</sup> T cells-mediated peripheral tolerance
<b>Cell Type</b>	Macrophages, Monocytes
<b>Biology Area</b>	Immunology, Innate Immunity

#### Antigen References

- Schaller E, *et al.* 2002. *Mol. Cell. Biol.* 22:8035.
- Stevceva L, *et al.* 2001. *BMC Clin Pathol.* 1:3.
- Kobayashi M, *et al.* 2008. *J. Leukoc. Biol.* 83:1354.
- Poeckel D, *et al.* 2009. *J. Biol Chem.* 284:21077.
- Glass AM, *et al.* 2013. *J. Immunol.* 190:4830.
- Koehm S, *et al.* 2007. *J. Allergy Clin. Immunol.* 120:570.
- Rankin AL, *et al.* 2010. *J. Immunol.* 184:1526.
- Sasi SP, *et al.* 2014. *J Biol Chem.* 289:14178.
- Thakus VS, *et al.* 2014. *Toxicol Lett.* 230:322.
- Ohnishi A, *et al.* 2015. *PLoS One.* 10:117616.

**Gene ID** [13733](#)


## Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

## Other Formats

Purified anti-mouse F4/80 Recombinant Antibody, PE anti-mouse F4/80 Recombinant Antibody, APC anti-mouse F4/80 Recombinant Antibody, PE/Cyanine7 anti-mouse F4/80 Recombinant Antibody, FITC anti-mouse F4/80 Recombinant Antibody, Spark YG™ 593 anti-mouse F4/80 Recombinant Antibody, APC/Cyanine7 anti-mouse F4/80 Recombinant Antibody, PerCP/Cyanine5.5 anti-mouse F4/80 Recombinant Antibody, Alexa Fluor® 647 anti-mouse F4/80 Recombinant Antibody, PE/Fire™ 640 anti-mouse F4/80 Recombinant Antibody, Spark Red™ 718 anti-mouse F4/80 Recombinant (Flexi-Fluor™) Antibody, PE/Fire™ 810 anti-mouse F4/80 Recombinant Antibody, APC/Fire™ 810 anti-mouse F4/80 Recombinant Antibody

## Product Data

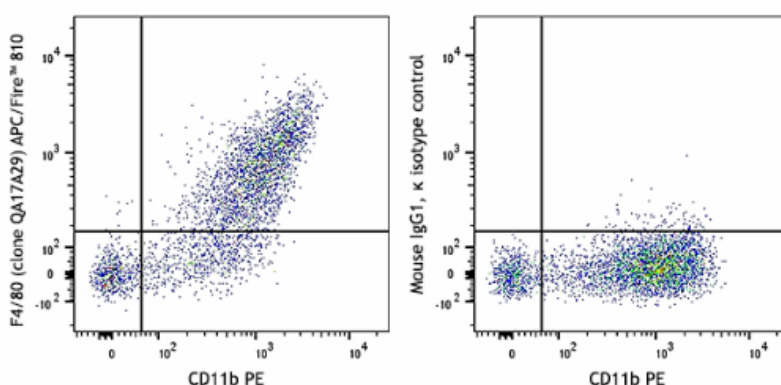

Thioglycollate-elicited BALB/c mouse peritoneal macrophages were stained with anti-mouse CD11b (clone M1/70) PE and anti-mouse F4/80 recombinant (clone QA17A29) APC/Fire™ 810 (left) or mouse IgG1, κ APC/Fire™ 810 isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at [www.biolegend.com/terms](http://www.biolegend.com/terms) ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587