

Brilliant Violet 421™ anti-human TCR Vβ13.1 Antibody

Catalog# / Size	362419 / 25 tests 362420 / 100 tests
Clone	H131
Regulatory Status	RUO
Other Names	T cell receptor Vβ13.1 chain, TCRBV13.1, TCRBV13
Isotype	Mouse IgG2b, κ
Description	TCR Vβ13.1 is a variant of the TCR β chain. The receptor is complexed with the TCR α chain and belongs to the immunoglobulin superfamily. It is expressed on a subset of T cells and some T cell clones. Variability in the β chain is generated by Vβ, Dβ, and Jβ gene rearrangement, while variability in the α chain is generated by Vα and Jα rearrangement. TCR Vβ13.1 has been shown to be related to antigen recognition and inflammation as well as diseases, such as HIV and multiple sclerosis.

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Mouse T cell hybridoma
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 421™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. Brilliant Violet 421™ excites at 405 nm and emits at 421 nm. The standard bandpass filter 450/50 nm is recommended for detection. Brilliant Violet 421™ is a trademark of Sirigen Group Ltd. Learn more about Brilliant Violet™. This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.
Excitation Laser	Violet Laser (405 nm)
Application Notes	Additional reported applications (for the relevant formats of this clone) include: activation of TCR Vβ13.1 expressing cells ¹ .

Application References

1. Choi YW, *et al.* 1991. *Proc. Natl. Acad. Sci. USA* 88:8357.

(PubMed link indicates
BioLegend citation)

Antigen Details

Structure	Immunoglobulin superfamily
Distribution	Subset of T cells
Function	Antigen recognition, activation of TCR V β 13.1 bearing T cells
Ligand/Receptor	Peptide bound MHC molecules
Cell Type	T cells
Biology Area	Adaptive Immunity, Immunology
Molecular Family	TCRs

Antigen References	1. Hong J, <i>et al.</i> 1999. <i>J. Immunol.</i> 163:3530. 2. Jason J, <i>et al.</i> 1997. <i>Scand. J. Immunol.</i> 45:81. 3. Chang JC, <i>et al.</i> 1995. <i>Ann. NY. Acad. Sci.</i> 756:370. 4. Miyahara E, <i>et al.</i> 1999. <i>Anticancer Res.</i> 19:2057. 5. Usuku K, <i>et al.</i> 1993. <i>Immunogenetics.</i> 38:193. 6. Bakakos P, <i>et al.</i> 2002. <i>Clin. Exp. Immunol.</i> 129:370.
--------------------	--

Gene ID [6957](#)


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

FITC anti-human TCR V β 13.1, PE/Cyanine7 anti-human TCR V β 13.1, Purified anti-human TCR V β 13.1, Brilliant Violet 605™ anti-human TCR V β 13.1 Antibody, PE anti-human TCR V β 13.1, TotalSeq™-A0908 anti-human TCR V β 13.1, TotalSeq™-C0908 anti-human TCR V β 13.1, Brilliant Violet 421™ anti-human TCR V β 13.1, APC anti-human TCR V β 13.1, TotalSeq™-B0908 anti-human TCR V β 13.1 Antibody

Product Data

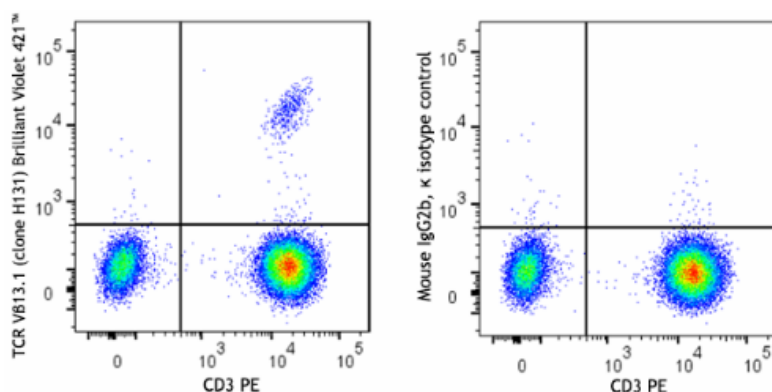

Human peripheral blood lymphocytes were stained with anti-human CD3 (clone UCHT1) PE and with anti-human TCR V β 13.1 (clone H131) Brilliant Violet 421™ (left) or mouse IgG2b, κ Brilliant Violet 421™ isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of

BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587