

Brilliant Violet 510™ anti-human CD11a Antibody

Catalog# / Size	301241 / 25 tests 301242 / 100 tests
Clone	HI111
Regulatory Status	RUO
Workshop	IV N231
Other Names	LFA-1 α chain, Lymphocyte function-associated antigen 1, α L Integrin, ITGAL
Isotype	Mouse IgG1, κ
Description	CD11a is a 170-180 kD type I transmembrane glycoprotein also known as LFA-1 α chain and integrin α_L subunit. CD11a non-covalently associates with integrin β_2 (CD18) to form LFA-1. It is expressed on all leukocytes, including B and T lymphocytes, monocytes, macrophages, neutrophils, basophils and eosinophils. It is absent on non-hematopoietic tissues and platelets. CD11a plays a central role in leukocyte cell-cell interactions and is important in lymphocyte costimulation. CD11a/CD18 binds to ICAM-1 (CD54), ICAM-2 (CD102), and ICAM-3 (CD50).

Product Details

Verified Reactivity	Human
Reported Reactivity	African Green, Baboon, Cow, Capuchin monkey, Chimpanzee, Cynomolgus, Dog, Horse, Rabbit, Sheep
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 510™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 μ L per million cells in 100 μ L staining volume or 5 μ L per 100 μ L of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.

Brilliant Violet 510™ excites at 405 nm and emits at 510 nm. The bandpass filter 510/50 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 510™ is a trademark of Sirigen Group Ltd.

[Learn more about Brilliant Violet™.](#)

This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.

Excitation Laser	Violet Laser (405 nm)
Application Notes	Clone HI111 epitope maps to the top region of the I domain that is close to the putative ligand-binding site surrounding the MIDAS (metal ion-dependent adhesion site). HI111 is specific for the closed confirmation of the integrin ⁸ . Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen sections, Western blotting ² , and blocking of cell-cell interaction and inhibition the binding of ICAM-1 ⁴ . This clone was tested in-house and does not work on formalin fixed paraffin-embedded (FFPE) tissue. The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/μg, Azide-Free, 0.2 μm filtered) is recommended for functional assays (Cat. No. 301233 & 301234).
Application References	1. Knapp W, <i>et al.</i> 1989. Leucocyte Typing IV. Oxford University Press New York.
(PubMed link indicates BioLegend citation)	2. Leite F, <i>et al.</i> 2002. <i>Infec. Immun.</i> 70:4336. 3. Jiang Y, <i>et al.</i> 2005. <i>Clin. Hemorheol. Microcircul.</i> 32:261. 4. BTchard D, <i>et al.</i> 2001. <i>J. Immunol.</i> 167:3099. 5. Sithu SD, <i>et al.</i> 2007. <i>J. Biol. Chem.</i> doi:10.1074/jbc.M611273200. 6. Choi EY, <i>et al.</i> 2008. <i>Blood</i> 111:3607. PubMed 7. Yoshino N, <i>et al.</i> 2000. <i>Exp. Anim. (Tokyo)</i> 49:97. (FC) 8. Ma Q, <i>et al.</i> 2002. <i>J. Biol. Chem.</i> 277:10638.
RRID	AB_3675024 (BioLegend Cat. No. 301241) AB_3675024 (BioLegend Cat. No. 301242)

Antigen Details

Structure	Integrin, type I transmembrane glycoprotein, noncovalently linked with integrin β ₂ (CD18) forms LFA-1, 170-180 kD
Distribution	Leukocytes
Function	Adhesion, costimulation
Ligand/Receptor	ICAM-1(CD54), ICAM-2(CD102), ICAM-3(CD50)
Cell Type	Leukocytes, Tregs
Biology Area	Cell Adhesion, Cell Biology, Costimulatory Molecules, Immunology, Innate Immunity, Neuroinflammation, Neuroscience
Molecular Family	Adhesion Molecules, CD Molecules
Antigen References	1. Lub M, <i>et al.</i> 1995. <i>Immunol. Today</i> 16:479. 2. Parsons J. 1996. <i>Curr. Opin. Cell Biol.</i> 8:146.
Gene ID	3683


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-human CD11a, Biotin anti-human CD11a, FITC anti-human CD11a, PE anti-human CD11a, PE/Cyanine5 anti-human CD11a, Purified anti-human CD11a, Alexa Fluor® 488 anti-human CD11a, Alexa Fluor® 647 anti-human CD11a, PE/Cyanine7 anti-human CD11a, Alexa Fluor® 594 anti-human CD11a, Purified anti-human CD11a (Maxpar® Ready), APC/Fire™ 750 anti-human CD11a, Alexa Fluor® 700 anti-human CD11a, PerCP/Cyanine5.5 anti-human CD11a, PE/Dazzle™ 594 anti-human CD11a, Ultra-LEAF™ Purified anti-human CD11a, APC/Cyanine7 anti-human CD11a, Brilliant Violet 711™ anti-human CD11a, Brilliant Violet 650™ anti-human CD11a, Spark Red™ 718 anti-human CD11a (Flexi-Fluor™), Brilliant Violet 510™ anti-human CD11a

Product Data

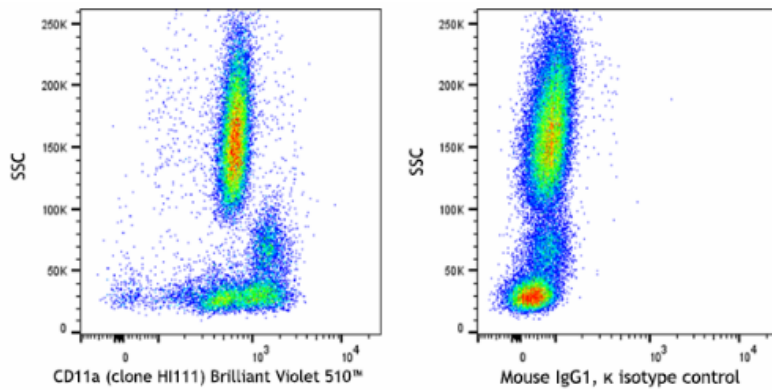

Human lysed whole blood was stained with anti-human CD11a (clone HI111) Brilliant Violet 510™ (left) or mouse IgG1, κ isotype control Brilliant Violet 510™ (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587