

Brilliant Violet 570™ anti-human CD4 Antibody

Catalog# / Size	344667 / 25 tests 344668 / 100 tests
Clone	SK3
Regulatory Status	RUO
Other Names	T4, Leu3a
Isotype	Mouse IgG1, κ
Description	CD4, also known as T4, is a 55 kD single-chain type I transmembrane glycoprotein expressed on most thymocytes, a subset of T cells, and monocytes/macrophages. CD4, a member of the Ig superfamily, recognizes antigens associated with MHC class II molecules and participates in cell-cell interactions, thymic differentiation, and signal transduction. CD4 acts as a primary receptor for HIV, binding to HIV gp120. CD4 has also been shown to interact with IL-16.

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 570™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>Brilliant Violet 570™ excites at 405 nm and emits at 570 nm. The bandpass filter 585/42 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 570™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)

Application References

(PubMed link indicates
BioLegend citation)

1. Evans RL, *et al.* 1981. *Immunol.* 78:544
2. Arno A *et al.* 1999. *J. Infect. Dis.* 180:56
3. Muech M, *et al.* 1997. *Blood* 89:1364
4. Wang L, *et al.* 2012. *Cytometry A.* 81:567. [PubMed](#)

RRID

AB_3683365 (BioLegend Cat. No. 344667)
AB_3683365 (BioLegend Cat. No. 344668)

Antigen Details

Structure	Ig superfamily, type I transmembrane glycoprotein, 55 kD
Distribution	T cell subset, majority of thymocytes, monocytes/macrophages
Function	MHC class II co-receptor, lymphocyte adhesion, thymic differentiation, HIV receptor
Ligand/Receptor	MHC class II molecules, HIV gp120, IL-16
Cell Type	Macrophages, Monocytes, T cells, Thymocytes, Tregs
Biology Area	Immunology
Molecular Family	CD Molecules
Antigen References	1. Center D <i>et al.</i> 1996. <i>Immunol. Today</i> 17:476. 2. Gaubin M <i>et al.</i> 1996. <i>Eur. J. Clin. Chem. Clin. Biochem.</i> 34:723.
Gene ID	920


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

PerCP anti-human CD4, Purified anti-human CD4, FITC anti-human CD4, PE anti-human CD4, PerCP/Cyanine5.5 anti-human CD4, Biotin anti-human CD4, PE/Cyanine7 anti-human CD4, APC anti-human CD4, APC/Cyanine7 anti-human CD4, Alexa Fluor® 488 anti-human CD4, Pacific Blue™ anti-human CD4, Alexa Fluor® 700 anti-human CD4, Purified anti-human CD4 (Maxpar® Ready), Brilliant Violet 421™ anti-human CD4, Brilliant Violet 510™ anti-human CD4, Alexa Fluor® 647 anti-human CD4, APC/Fire™ 750 anti-human CD4, PE/Dazzle™ 594 anti-human CD4, Brilliant Violet 785™ anti-human CD4, Brilliant Violet 750™ anti-human CD4, Brilliant Violet 605™ anti-human CD4, Brilliant Violet 711™ anti-human CD4, TotalSeq™-A0045 anti-human CD4, TotalSeq™-C0045 anti-human CD4, PE/Cyanine5 anti-human CD4, Spark Blue™ 550 anti-human CD4, Spark NIR™ 685 anti-human CD4, KIRAVIA Blue 520™ anti-human CD4, APC/Fire™ 810 anti-human CD4, PE/Fire™ 640 anti-human CD4, PE/Fire™ 700 anti-human CD4, Spark Violet™ 538 anti-human CD4 Antibody, Spark YG™ 581 anti-human CD4, Alexa Fluor® 660 anti-human CD4, Spark YG™ 593 anti-human CD4, PE/Fire™ 810 anti-human CD4 Antibody, APC/Fire™ 750 anti-human CD4, Spark Blue™ 574 anti-human CD4, GMP APC anti-human CD4, TotalSeq™-B0045 anti-human CD4 Antibody, Spark Violet™ 423 anti-human CD4 Antibody, GMP Pacific Blue™ anti-human CD4, Spark UV™ 387 anti-human CD4, GMP PE anti-human CD4, GMP FITC anti-human CD4, Spark Red™ 718 anti-human CD4, GMP PerCP anti-human CD4, GMP PE/Cyanine7 anti-human CD4, Spark Violet™ 500 anti-human CD4, GMP PerCP/Cyanine5.5 anti-human CD4, Brilliant Violet 650™ anti-human CD4, PerCP/Fire™ 806 anti-human CD4, Spark Blue™ 515 anti-human CD4, GMP APC/Fire™ 750 anti-human CD4, Spark Violet™ 423 anti-human CD4, PerCP/Fire™ 780 anti-human CD4, PE/Fire™ 744 anti-human CD4, Spark PLUS UV395™ anti-human CD4, Spark PLUS B550™ anti-human CD4, Brilliant Violet 570™ anti-human CD4, GMP Spark Violet™ 423 anti-human CD4, StarBright UltraViolet 575 anti-human CD4, StarBright UltraViolet 740 anti-human CD4, StarBright UltraViolet 795 anti-human CD4, Spark PLUS V475™ anti-human CD4 Antibody

Product Data

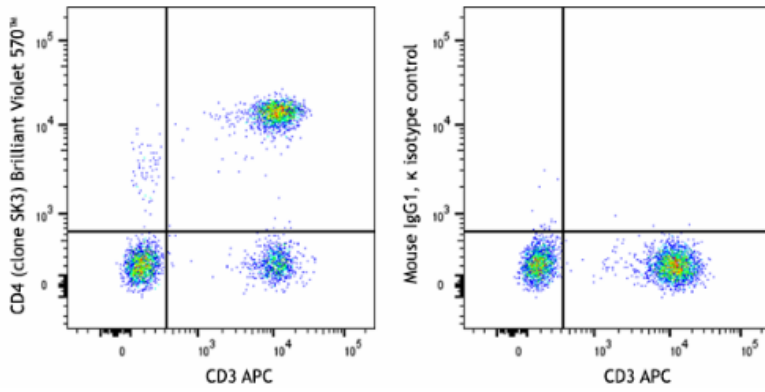

Human peripheral blood lymphocytes were stained anti-human CD3 (clone UCHT1) APC and anti-human CD4 (clone SK3) Brilliant Violet 570™ (left) or mouse IgG1, κ Brilliant Violet 570™ isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587