

Brilliant Violet 570™ anti-human CD45 Antibody

Catalog# / Size	368565 / 25 tests 368566 / 100 tests
Clone	2D1
Regulatory Status	RUO
Other Names	Leukocyte Common Antigen (LCA), T200
Isotype	Mouse IgG1, κ
Description	CD45 is a 180 - 240 kD single chain type I membrane glycoprotein also known as leukocyte common antigen (LCA) and T200. It is a tyrosine phosphatase expressed on the plasma membrane of all hematopoietic cells, except erythrocytes or platelets. CD45 is a signaling molecule that regulates a variety of cellular processes including cell growth, differentiation, cell cycle, and oncogenic transformation. CD45 plays a critical role in T and B cell antigen receptor-mediated activation by dephosphorylating substrates including p56Lck, p59Fyn, and other Src family kinases. CD45 non-covalently associates with lymphocyte phosphatase-associated phosphoprotein (LPAP) on T and B lymphocytes. CD45 has been reported to bind galectin-1 and to be associated with several other cell surface antigens including CD1, CD2, CD3, and CD4.

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Human PBMC
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 570™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>Brilliant Violet 570™ excites at 405 nm and emits at 570 nm. The bandpass filter 585/42 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 570™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)
Application Notes	It was found that the HI30 clone and the 2D1 clone can cross block each other's binding.

Application References

(PubMed link indicates BioLegend citation)

1. Bradstock KF, *et al.* 1980. *J. Natl. Cancer Inst.* 65:33.
2. Csiba A, *et al.* 1984. *Br. J. Cancer* 50:699.
3. Tchilian EZ, *et al.* 2001. *J. Immunol.* 166:1308.
4. Lee MS, *et al.* 2004. *Int. Immunol.* 16:1109.

RRID AB_3699010 (BioLegend Cat. No. 368565)
AB_3699010 (BioLegend Cat. No. 368566)

Antigen Details

Structure	Type I transmembrane protein, 180 - 240 kD (multiple isoforms)
Distribution	Hematopoietic cells, except circulating erythrocytes or platelets.
Function	Tyrosine phosphatases, signaling, co-stimulation (co-inhibition), TCR and BCR mediated activation.
Ligand/Receptor	Galectin-1, CD2, CD3, and CD4.
Cell Type	B cells, Dendritic cells, Neutrophils
Biology Area	Cell Biology, Immunology, Inhibitory Molecules, Neuroscience, Neuroscience Cell Markers
Molecular Family	CD Molecules, Protein Kinases/Phosphatase, TCRs
Antigen References	1. Thomas M. 1989. <i>Annu. Rev. Immunol.</i> 7:339. 2. Trowbridge I, <i>et al.</i> 1994. <i>Annu. Rev. Immunol.</i> 12:85.
Gene ID	5788


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)


Other Formats

Purified anti-human CD45, PerCP anti-human CD45, FITC anti-human CD45, PE anti-human CD45, APC anti-human CD45, Alexa Fluor® 700 anti-human CD45, APC/Cyanine7 anti-human CD45, PerCP/Cyanine5.5 anti-human CD45, APC/Fire™ 750 anti-human CD45, Alexa Fluor® 594 anti-human CD45, Brilliant Violet 510™ anti-human CD45, Brilliant Violet 421™ anti-human CD45, PE/Dazzle™ 594 anti-human CD45, Brilliant Violet 605™ anti-human CD45, PE/Cyanine7 anti-human CD45, Pacific Blue™ anti-human CD45, Alexa Fluor® 647 anti-human CD45, Biotin anti-human CD45, Alexa Fluor® 488 anti-human CD45, TotalSeq™-A0048 anti-human CD45, Brilliant Violet 750™ anti-human CD45, TotalSeq™-C0048 anti-human CD45, TotalSeq™-B0048 anti-human CD45, Spark Blue™ 550 anti-human CD45, Spark NIR™ 685 anti-human CD45, TotalSeq™-D0048 anti-human CD45, PE/Cyanine5 anti-human CD45, Spark Blue™ 574 anti-human CD45 Antibody, Spark Violet™ 500 anti-human CD45, Spark Red™ 718 anti-human CD45, APC/Fire™ 810 anti-human CD45 Antibody, Spark YG™ 581 anti-human CD45 (Flexi-Fluor™), Brilliant Violet 570™ anti-human CD45

Product Data

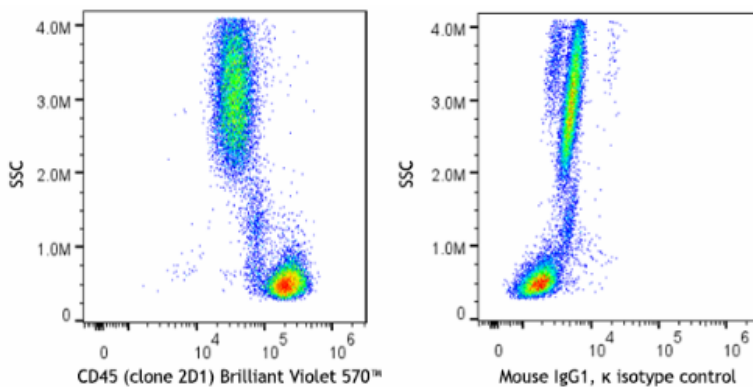

Human lysed whole blood was stained with anti-human CD45 (clone 2D1) Brilliant Violet 570™ (left) or mouse IgG1, κ Brilliant Violet 570™ isotype control (right).

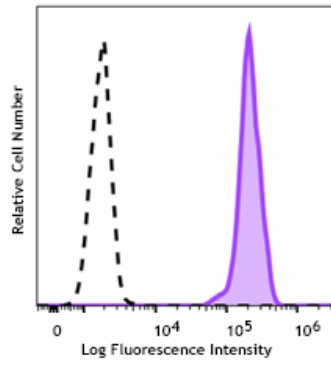

Human peripheral blood lymphocytes were stained with anti-human CD45 (clone 2D1) Brilliant Violet 570™ (filled histogram) or mouse IgG1, κ Brilliant Violet 570™ isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587