

Brilliant Violet 605™ anti-human IFN-γ Antibody

Catalog# / Size	383409 / 25 tests 383410 / 100 tests
Clone	W19227C
Regulatory Status	RUO
Other Names	Interferon-γ, Immune interferon, Type II interferon, T cell interferon, Macrophage-activating factor (MAF), IFN-g, IFN-gamma
Isotype	Rat IgG2b, λ
Description	Human interferon-gamma (IFN-γ, type II interferon) is a potent multifunctional cytokine which is secreted primarily by activated NK cells and T cells in response to antigens and other cytokines. Originally characterized based on anti-viral activities. IFN-γ also exerts tumor anti-proliferative, immunoregulatory, and proinflammatory activities. IFN-γ can upregulate MHC class I and II, and macrophage activation. This human IFN-γ antibody does not show cross-reactivity with mouse.

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Rat
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 605™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool).
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	ICFC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by intracellular flow cytometry using our Cyto-Fast™ Fix/Perm Buffer Set. For flow cytometric staining, the suggested use of this reagent is 5 µL per 10⁶ cells in 100 µL volume. It is highly recommended that the reagent be titrated for optimal performance for each application. Brilliant Violet 605™ excites at 405 nm and emits at 603 nm. The bandpass filter 610/20 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 605™ is a trademark of Sirigen Group Ltd. Learn more about Brilliant Violet™. This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.
Excitation Laser	Violet Laser (405 nm)

RRID AB_3720571 (BioLegend Cat. No. 383409)
AB_3720571 (BioLegend Cat. No. 383410)

Antigen Details

Structure	Cytokine; dimer; 20-25 kD (Mammalian)
Distribution	CD8 ⁺ and CD4 ⁺ T cells, NK cells
Function	Upregulated by IL-2, FGF-basic, EGF; downregulated by vitamin D3 or DMN; labile at pH 2
Interaction	T cells, B cells, macrophages, NK cells, endothelial cells, fibroblasts
Ligand/Receptor	IFN- γ R α (CDw119) dimerized with IFN- γ R β (AF-1)
Cell Type	Tregs
Biology Area	Cell Biology, Immunology, Neuroinflammation, Neuroscience
Molecular Family	Cytokines/Chemokines

Antigen References

1. Fitzgerald K, *et al.* Eds. 2001. The Cytokine FactsBook. Academic Press, San Diego.
2. De Maeyer E, *et al.* 1992. *Curr Opin Immunol.* 4:321-6.
3. Farrar M, *et al.* 1993. *Annu Rev Immunol.* 11:571-611.
4. Gray P, *et al.* 1987. *Lymphokines.* 13:151.

Gene ID [3458](#)


Related Protocols

- [Intracellular Flow Cytometry Staining Protocol](#)

Other Formats

PE anti-human IFN- γ , Brilliant Violet 421™ anti-human IFN- γ Antibody, Brilliant Violet 605™ anti-human IFN- γ Antibody, PE/Cyanine7 anti-human IFN- γ

Product Data

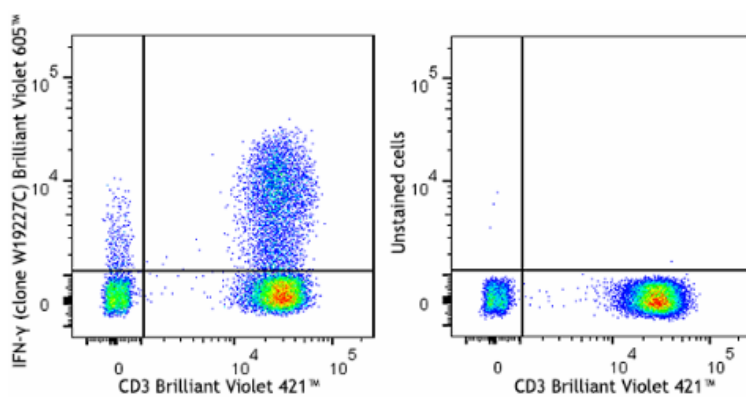

Human peripheral blood cells were stimulated with Cell Activation Cocktail (with Brefeldin A) (Cat. No. 423303) for 6 hours and surface stained with anti-human CD3 (clone UCHT1) Brilliant Violet 421™. Cells were then fixed and permeabilized with Cyto-Fast™ Fix/Perm Buffer Set (Cat. No. 426803) and intracellularly stained with anti-human IFN- γ (clone W19227C) Brilliant Violet 605™ (left) or left stained with anti-human CD3 (clone UCHT1) Brilliant Violet 421™ only (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biollegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587