

Brilliant Violet 650™ anti-mouse CD19 Antibody

Catalog# / Size	152427 / 50 µg
Clone	1D3/CD19
Regulatory Status	RUO
Other Names	B4, CVID3
Isotype	Rat IgG2a, κ
Description	CD19 is a 95 kD glycoprotein also known as B4. It is a member of the Ig superfamily, expressed on all pro-B to mature B cells (during development) and follicular dendritic cells. Plasma cells do not express CD19. CD19, in association with CD21 and CD81, forms a molecular complex integral to B cell activation.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Rat
Immunogen	Transfected k562 cells expressing the extracellular domain of murine CD19.
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 650™ under optimal conditions.
Concentration	0.2 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is ≤ 0.25 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>Brilliant Violet 650™ excites at 405 nm and emits at 645 nm. The bandpass filter 660/20 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 650™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)
Application Notes	Additional reported applications (for the relevant formats) include: blocking ^{1,2} and immunoprecipitation ¹ .

Application References

(PubMed link indicates
BioLegend citation)

1. Krop I, et al. 1996. *Eur. J. Immunol.* 26:238. (FC, Blocking, IP)
2. Abraham P, et al. 2014. *Clin. Exp. Immunol.* 175:181. (Blocking)

RRID

AB_3675002 (BioLegend Cat. No. 152427)

Antigen Details

Structure	Ig superfamily, associates with CD21 and CD81, 95 kD.
Distribution	Pro-B cells to mature B cells (during development), follicular dendritic cells.
Function	Modulates B cell activation and differentiation.
Ligand/Receptor	CD21, CD81, Leu-13.
Cell Type	B cells, Dendritic cells
Biology Area	Costimulatory Molecules, Immunology
Molecular Family	CD Molecules
Antigen References	<ol style="list-style-type: none">1. Fearon DT. 1993. <i>Curr. Opin. Immunol.</i> 5:341.2. Krop I, et al. 1996. <i>Eur. J. Immunol.</i> 26:238.3. Krop I, et al. 1996. <i>J. Immunol.</i> 157:48.4. Tedder TF, et al. 1994. <i>Immunol. Today</i> 15:437.
Gene ID	12478


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Purified anti-mouse CD19, FITC anti-mouse CD19, PerCP/Cyanine5.5 anti-mouse CD19, PE anti-mouse CD19, APC anti-mouse CD19, APC/Cyanine7 anti-mouse CD19, Alexa Fluor® 700 anti-mouse CD19, Pacific Blue™ anti-mouse CD19, PE/Cyanine7 anti-mouse CD19, Biotin anti-mouse CD19, Brilliant Violet 421™ anti-mouse CD19, Brilliant Violet 510™ anti-mouse CD19, Brilliant Violet 605™ anti-mouse CD19, Brilliant Violet 650™ anti-mouse CD19, Brilliant Violet 711™ anti-mouse CD19, Brilliant Violet 785™ anti-mouse CD19, Spark PLUS B550™ anti-mouse CD19, Spark YG™ 593 anti-mouse CD19 (Flexi-Fluor™) Antibody

Product Data

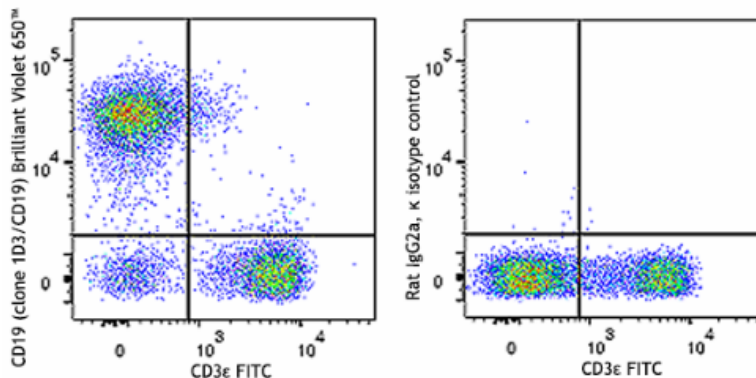

C57BL/6 mouse splenocytes were stained with anti-mouse CD3ε (clone 145-2C11) FITC and anti-mouse CD19 (clone 1D3/CD19) Brilliant Violet 650™ (left) or rat IgG2a, κ Brilliant Violet 650™ isotype control (right).

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587