

Brilliant Violet 750™ anti-human CD103 (Integrin αE) Antibody

Catalog# / Size	350263 / 25 tests 350264 / 100 tests
Clone	Ber-ACT8
Regulatory Status	RUO
Workshop	V A067
Other Names	Integrin alpha E (ITGAE)
Isotype	Mouse IgG1, κ
Description	CD103 is a type I transmembrane glycoprotein also known as αE integrin, integrin αIEL chain, and human mucosal lymphocyte antigen 1. It belongs to the integrin family and is primarily found on intestinal intraepithelial lymphocytes (IEL). CD103 is also expressed on a subpopulation of lamina propria T cells, epithelial dendritic cells, lamina propria-derived dendritic cells, and a small subset of peripheral lymphocytes. Treg cells express high level of CD103. Hairy cell leukemia has also been shown to express CD103. The expression of CD103 on lymphocytes can be induced upon activation and TGF-β stimulation. In association with integrin β7, CD103 is expressed as an αE/β7 heterodimer. Mature CD103 protein can be cleaved into 2 chains, a 150 kD (C-terminal) chain and a 25 kD (N-terminal) chain, which remain linked by disulfide bonds. CD103 binds to E-cadherin and mediates homing of lymphocytes to the intestinal epithelium.

Product Details

Verified Reactivity	Human
Reported Reactivity	Cynomolgus
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	HTLV-1 induced human T cell line MAPS16
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 750™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.

Brilliant Violet 750™ excites at 405 nm and emits at 750 nm. The bandpass filter 780/60 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 750™ is a trademark of Sirigen Group Ltd.

[Learn more about Brilliant Violet™.](#)

This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.

Excitation Laser	Violet Laser (405 nm)
Application Notes	Additional reported applications (for the relevant formats) include: Western Blotting ¹ , immunoprecipitation ¹ , and immunohistochemical staining of frozen tissue sections ¹ .
Application References (PubMed link indicates BioLegend citation)	1. Kruschwitz M, <i>et al.</i> 1991. <i>J. Clin. Pathol.</i> 44:636. (WB, IP, IHC-F) 2. Yoshino N, <i>et al.</i> 2000. <i>Exp. Anim. (Tokyo)</i> 49:97. (FC)
RRID	AB_3717156 (BioLegend Cat. No. 350263) AB_3717156 (BioLegend Cat. No. 350264)

Antigen Details

Structure	Type I transmembrane glycoprotein, integrin family; can be cleaved into 150 kD and 25 kD chains; associated with $\beta 7$ integrin
Distribution	Majority of intestinal intraepithelial lymphocytes (IEL), subpopulation of lamina propria T cells, epithelial dendritic cells, small subset of peripheral lymphocytes, Treg cells; expressed on hairy cell leukemia
Function	Retention and activation of CD103 ⁺ lymphocytes in the intestinal epithelium, regulation of tissue-specific T cell homing
Ligand/Receptor	E-Cadherin
Cell Sources	Integrin $\beta 7$
Cell Targets	Integrin $\beta 7$
Cell Type	Dendritic cells, Lymphocytes, T cells, Tregs
Biology Area	Cell Biology, Immunology, Neuroscience, Synaptic Biology
Molecular Family	Adhesion Molecules, CD Molecules
Antigen References	1. Parker CM, <i>et al.</i> 1992. <i>P. Natl. Acad. Sci. USA</i> 89:1924. 2. Kruschwitz M, <i>et al.</i> 1991. <i>J. Clin. Pathol.</i> 44:636. 3. Schon MP, <i>et al.</i> 1999. <i>J. Immunol.</i> 162:6641. 4. Shaw SK, <i>et al.</i> 1994. <i>J. Biol. Chem.</i> 269:6016.
Gene ID	3682

Related Protocols


- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

PE/Fire™ 700 anti-human CD103 (Integrin αE) Antibody, TotalSeq™-D0145 anti-human CD103 (Integrin αE), PE/Fire™ 640 anti-human CD103 (Integrin αE), FITC anti-human CD103, PE anti-human CD103, Spark YG™ 581 anti-human CD103 (Integrin αE), Biotin anti-human CD103 (Integrin αE), PerCP/Cyanine5.5 anti-human CD103 (Integrin αE), Alexa Fluor® 647 anti-human CD103 (Integrin αE), Purified anti-human CD103 (Integrin αE), FITC anti-human CD103 (Integrin αE), PE anti-human CD103 (Integrin αE), PE/Cyanine7 anti-human CD103 (Integrin αE), Alexa Fluor® 488 anti-human CD103 (Integrin αE), Brilliant Violet 421™ anti-human CD103 (Integrin αE), APC anti-human CD103 (Integrin αE), Spark Blue™ 574 anti-human CD103 (Integrin αE) (Flexi-Fluor™), PerCP/Fire™ 806 anti-human CD103 (Integrin, αE) Antibody, APC anti-human CD103, Spark NIR™ 685 anti-human CD103 (Integrin αE), Spark Blue™ 550 anti-human CD103 (Integrin αE), PE/Fire™ 810 anti-human CD103 (Integrin αE), GMP FITC anti-human CD103 (Integrin αE), GMP PE anti-human CD103 (Integrin αE), Spark Red™ 718 anti-human CD103 (Integrin αE) (Flexi-

Fluor™), Brilliant Violet 605™ anti-human CD103 (Integrin αE), Spark YG™ 593 anti-human CD103 (Integrin αE) (Flexi-Fluor™), TotalSeq™-A0145 anti-human CD103 (Integrin αE), Brilliant Violet 711™ anti-human CD103 (Integrin αE), TotalSeq™-C0145 anti-human CD103 (Integrin αE), GMP APC anti-human CD103 (Integrin αE), PE/Fire™ 744 anti-human CD103 (Integrin αE), APC/Cyanine7 anti-human CD103 (Integrin αE), PE/Dazzle™ 594 anti-human CD103 (Integrin αE), PerCP/Fire™ 780 anti-human CD103 (Integrin, αE), StarBright UltraViolet 795 anti-human CD103 (Integrin αE), Spark PLUS UV395™ anti-human CD103 (Integrin αE) Antibody, Brilliant Violet 750™ anti-human CD103 (Integrin αE), Spark PLUS V475™ anti-human CD103 (Integrin αE), Brilliant Violet 785™ anti-human CD103 (Integrin αE), TotalSeq™-B0145 anti-human CD103 (Integrin αE), APC/Fire™ 750 anti-human CD103 (Integrin αE)

Product Data

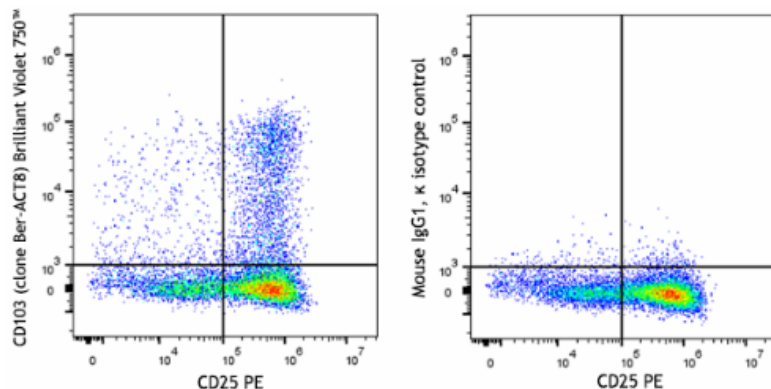

PHA-stimulated (3 days) human peripheral blood mononuclear cells were stained with anti-human CD25 (clone BC96) PE and anti-human CD103 (Integrin αE) (clone Ber-ACT8) Brilliant Violet 750™ (left) or mouse IgG1, κ Brilliant Violet 750™ isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587