

Brilliant Violet 750™ anti-human CD33 Antibody

Catalog# / Size	366649 / 25 tests
Clone	P67.6
Regulatory Status	RUO
Other Names	Siglec-3, gp67, p67
Isotype	Mouse IgG1, κ
Description	CD33, also known as Siglec-3, gp67, and p67, is a 67 kD type I transmembrane glycoprotein. It is a sialoadhesion immunoglobulin superfamily member, which is expressed on myeloid progenitors, monocytes, granulocytes, dendritic cells, and mast cells. CD33 is absent on normal platelets, lymphocytes, erythrocytes, and hematopoietic stem cells. CD33 functions as a sialic acid-dependent cell adhesion molecule with carbohydrate/lectin binding activity.

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	FMY9S5 cells expressing CD33 gene.
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 750™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. Brilliant Violet 750™ excites at 405 nm and emits at 750 nm. The bandpass filter 780/60 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 750™ is a trademark of Sirigen Group Ltd. Learn more about Brilliant Violet™. This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.
Excitation Laser	Violet Laser (405 nm)

Application References

1. Hoyer J, *et al.* 2008. *Am. J. Clin. Pathol.* 129:316.

(PubMed link indicates
BioLegend citation)

Antigen Details

Structure	Ig superfamily, sialoadhesins, type I glycoprotein, 67 kD
Distribution	Myeloid progenitor, monocytes, granulocytes, dendritic cells, and mast cells.
Function	Adhesion and lectin binding activity.
Ligand/Receptor	Sugar chains containing sialic acid.
Cell Type	Dendritic cells, Granulocytes, Hematopoietic stem and progenitors, Mast cells, Monocytes
Biology Area	Cell Biology, Immunology, Neuroinflammation, Neuroscience
Molecular Family	CD Molecules, Siglec Molecules
Antigen References	1. Favaloro E, <i>et al.</i> 1988. <i>Br. J. Haematol.</i> 69:163. 2. Freeman S, <i>et al.</i> 1995. <i>Blood</i> 85:2005.
Gene ID	945


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

TotalSeq™-D0052 anti-human CD33, Spark YG™ 581 anti-human CD33, Spark Violet™ 423 anti-human CD33, Spark Blue™ 574 anti-human CD33, Spark Red™ 718 anti-human CD33 (Flexi-Fluor™), APC anti-human CD33, PE anti-human CD33, FITC anti-human CD33, Purified anti-CD33, PE/Cyanine7 anti-human CD33, Biotin anti-human CD33, Brilliant Violet 421™ anti-human CD33, Alexa Fluor® 647 anti-human CD33, Brilliant Violet 711™ anti-human CD33, Brilliant Violet 570™ anti-human CD33, Brilliant Violet 650™ anti-human CD33, Purified anti-human CD33, Brilliant Violet 510™ anti-human CD33, Brilliant Violet 605™ anti-human CD33, APC/Cyanine7 anti-human CD33, PerCP/Cyanine5.5 anti-human CD33, APC/Fire™ 750 anti-human CD33, TotalSeq™-C0052 anti-human CD33, Spark Blue™ 550 anti-human CD33 (Flexi-Fluor™), Brilliant Violet 750™ anti-human CD33, Spark YG™ 593 anti-human CD33 (Flexi-Fluor™), Spark NIR™ 685 anti-human CD33 (Flexi-Fluor™), TotalSeq™-A0052 anti-human CD33, TotalSeq™-B0052 anti-human CD33

Product Data

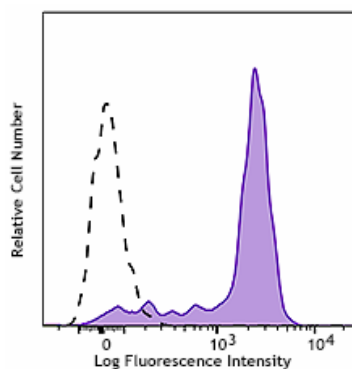

Human peripheral blood monocytes were stained with anti-human CD33 (clone P67.6) Brilliant Violet 750™ (filled histogram) or mouse IgG1, κ Brilliant Violet 750™ isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form,

used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587