

Brilliant Violet 750™ anti-human TIGIT (VSTM3) Antibody

Catalog# / Size	373751 / 25 tests 373752 / 100 tests
Clone	A15153G
Regulatory Status	RUO
Other Names	T-cell immunoreceptor with Ig and ITIM domains, VSIG9, VSTM3, WUCAM
Isotype	Mouse IgG2a, κ
Description	T cell immunoreceptor with Ig and ITIM domains (TIGIT), also known as VSTM3 or WUCAM, is a 26 kD, type I transmembrane protein and is a member of the PVR (poliovirus receptor) family of immunoglobulin-like domain containing proteins. TIGIT is expressed on activated T cells, follicular T helper, memory, and regulatory T cells as well as on NK cells. TIGIT is a negative regulator of NK and T cell activation. Expression of TIGIT is associated with decreased functionality of CD8 T cells in chronic viral infection and tumors. TIGIT also promotes the differentiation of tolerogenic phenotype in dendritic cells with an increased secretion of IL-10 and a diminished production of IL-12.

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Recombinant Human TIGIT.
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 750™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.

Brilliant Violet 750™ excites at 405 nm and emits at 750 nm. The bandpass filter 780/60 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 750™ is a trademark of Sirigen Group Ltd.

[Learn more about Brilliant Violet™.](#)

This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly

prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.

Excitation Laser	Violet Laser (405 nm)
Application Notes	This clone can suppress anti-CD3 induced T cell proliferation <i>in vitro</i> based on in-house testing. This clone has been tested in-house and determined to not be suitable for applications in immunohistochemistry of paraffin-embedded tissue sections (IHC-P). Additional reported applications (for the relevant formats) include: Blocking¹.
Application References	1. Stamm H, <i>et al.</i> 2018. <i>Oncogene</i> . Pubmed
(PubMed link indicates BioLegend citation)	
RRID	AB_3717202 (BioLegend Cat. No. 373751) AB_3717202 (BioLegend Cat. No. 373752)

Antigen Details

Structure	26kD; type I transmembrane protein, Ig-like V-type domain, ITIM motif.
Distribution	Activated T cells, Regulatory T cells (Treg), Follicular Helper T cells (TFH), NK cells.
Function	Cell signaling, negative regulation of T cells, T cell tolerance, T cell anergy.
Ligand/Receptor	CD155 (PVR), CD112 (PVRL2, NECTIN-2).
Cell Type	NK cells, T cells, Tfh, Tregs
Biology Area	Cell Adhesion, Cell Biology, Immunology, Inhibitory Molecules, Signal Transduction
Molecular Family	Adhesion Molecules, Immune Checkpoint Receptors
Antigen References	1. Stanitsky N, <i>et al.</i> 2009. <i>Proc. Natl. Acad. Sci.</i> 106:17858. 2. Yu X, <i>et al.</i> 2009. <i>Nat. Immunol.</i> 10:48. 3. Johnston R, <i>et al.</i> 2014. <i>Cancer Cell.</i> 26:923.
Gene ID	201633


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Brilliant Violet 711™ anti-human TIGIT (VSTM3), PE/Fire™ 640 anti-human TIGIT (VSTM3), PE/Cyanine5 anti-human TIGIT (VSTM3), Brilliant Violet 785™ anti-human TIGIT (VSTM3) Antibody, PerCP/Fire™ 806 anti-human TIGIT (VSTM3), TotalSeq™-D0089 anti-human TIGIT (VSTM3), PerCP/Fire™ 780 anti-human TIGIT (VSTM3), Spark Red™ 718 anti-human TIGIT (VSTM3) (Flexi-Fluor™), Brilliant Violet 421™ anti-human TIGIT (VSTM3), APC/Fire™ 750 anti-human TIGIT (VSTM3), APC anti-human TIGIT (VSTM3), PE anti-human TIGIT (VSTM3), Brilliant Violet 605™ anti-human TIGIT (VSTM3), PE/Dazzle™ 594 anti-human TIGIT (VSTM3), Alexa Fluor® 647 anti-human TIGIT (VSTM3), Purified anti-human TIGIT (VSTM3), Biotin anti-human TIGIT (VSTM3), PE/Fire™ 810 anti-human TIGIT (VSTM3), PE/Fire™ 700 anti-human TIGIT (VSTM3), Brilliant Violet 650™ anti-human TIGIT (VSTM3) Antibody, PE/Fire™ 744 anti-human TIGIT (VSTM3) Antibody, TotalSeq™-C0089 anti-human TIGIT (VSTM3), PerCP/Cyanine5.5 anti-human TIGIT (VSTM3), PE/Cyanine7 anti-human TIGIT (VSTM3), Ultra-LEAF™ Purified anti-human TIGIT (VSTM3), APC/Cyanine7 anti-human TIGIT (VSTM3), TotalSeq™-B0089 anti-human TIGIT (VSTM3), KIRAVIA Blue 520™ anti-human TIGIT (VSTM3), StarBright UltraViolet 740 anti-human TIGIT (VSTM3), Brilliant Violet 750™ anti-human TIGIT (VSTM3), Spark PLUS V475™ anti-human TIGIT (VSTM3), TotalSeq™-A0089 anti-human TIGIT (VSTM3), Brilliant Violet 510™ anti-human TIGIT (VSTM3)

Product Data

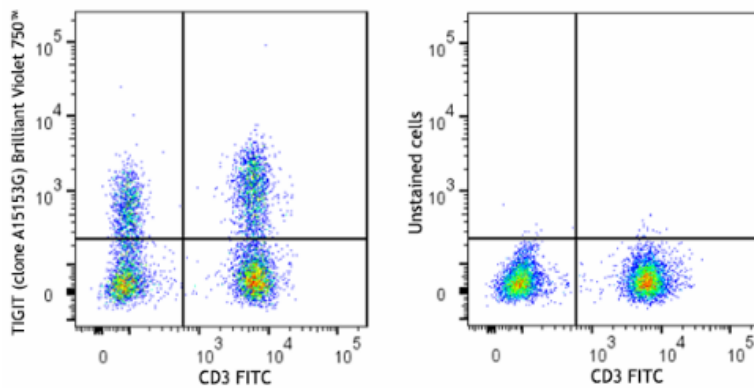

Human peripheral blood lymphocytes were stained with anti-human CD3 (clone UCHT1) FITC and anti-human TIGIT (VSTM3) (clone A15153G) Brilliant Violet 750™ (left) or stained with anti-human CD3 (clone UCHT1) FITC only (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587