

Brilliant Violet 750™ anti-mouse Ly-6A/E (Sca-1) Antibody

Catalog# / Size 108159 / 50 µg

Clone D7

Regulatory Status RUO

Other Names Sca-1

Isotype Rat IgG2a, κ

Description Ly-6A/E, also known as Sca-1, is an 18 kD member of the Ly-6 multigene family. Ly6A/E is a glycosylphosphatidylinositol (GPI)-linked protein expressed on hematopoietic stem cells. In mice expressing the Ly-6.2 haplotype (e.g., AKR, C57BL, C57BR, DBA/2, SJL, SWR, and 129), Ly-6A/E is also expressed on peripheral B lymphocytes and thymic and peripheral T lymphocytes. Strains expressing the Ly-6.1 haplotype (e.g., BALB/c, CBA, C3H/He, DBA/1, and NZB) have low Ly-6A/E expression on resting peripheral lymphocytes. The expression of Ly-6A/E on lymphocytes is upregulated upon activation from both Ly6.1 and Ly6.2 haplotype mice. Ly-6A/E is thought to be involved in the regulation of both T and B cell responses.

Product Details

Verified Reactivity Mouse

Antibody Type Monoclonal

Host Species Rat

Immunogen IL-2-dependent mouse T-cell line (CTL-L)

Formulation Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)

Preparation The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 750™ under optimal conditions.

Concentration 0.2 mg/mL

Storage & Handling The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. **Do not freeze.**

Application [FC - Quality tested](#)

Recommended Usage Each lot of this antibody is quality control tested by [immunofluorescent staining with flow cytometric analysis](#). For flow cytometric staining, the suggested use of this reagent is ≤ 0.125 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application.

Brilliant Violet 750™ excites at 405 nm and emits at 750 nm. The bandpass filter 780/60 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 750™ is a trademark of Sirigen Group Ltd.

[Learn more about Brilliant Violet™.](#)

This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.

Excitation Laser Violet Laser (405 nm)

Application Notes

The D7 antibody has been reported to induce T cell activation and inhibit TCR-induced IL-2 production. Additional reported applications (for the relevant formats) include: Western blotting^{1,2}, immunoprecipitation¹, *in vitro* lymphocyte activation³⁻⁶, induction of redirected lysis⁷, induction of T cell inhibitory signalling⁸, immunofluorescence⁹, and immunohistochemical staining of acetone-fixed frozen sections¹³ and Bouin-fixed, paraffin-embedded samples⁹.

The two Sca-1 recognizing clones D7 and [E13-161.7](#) have been shown to bind distinct epitopes due to the inability of D7 to block the binding of E13-161.7.¹⁴

Application References

(PubMed link indicates BioLegend citation)

1. Ortega G, *et al.* 1986. *J. Immunol.* 137:3240. (WB, IP)
2. Palfree RGE, *et al.* 1986. *Immunogenetics* 23:197. (WB)
3. Codias EK, *et al.* 1990. *J. Immunol.* 144:2197.
4. Malek TR, *et al.* 1986. *J. Exp. Med.* 164:709.
5. Codias EK, *et al.* 1990. *J. Immunol.* 145:1407.
6. Ivanov V, *et al.* 1994. *J. Immunol.* 153:2394.
7. Karlhofer FM, *et al.* 1991. *J. Immunol.* 146:3662.
8. Fleming T, *et al.* 1994. *J. Immunol.* 153:1955.
9. van Bragt MPA, *et al.* 2005. *Biol. Reprod.* 73:634. (IF, IHC)
10. Umland O, *et al.* 2007. *J. Immunol.* 178:4147.
11. Cridland SO, *et al.* 2009. *Blood Cell. Mol. Dis.* 45:149. (FC) [PubMed](#)
12. Pronk CJ, *et al.* 2011. *J. Exp Med.* [PubMed](#)

[See More](#)

RRID

AB_3698894 (BioLegend Cat. No. 108159)

Antigen Details

Structure	Ly-6 multigene family, 18 kD
Distribution	Hematopoietic stem cells, activated T cells and B cells, subset of resting B cells and T cells
Function	Regulates B and T cell responses
Cell Type	B cells, Hematopoietic stem and progenitors, Mesenchymal Stem Cells, T cells
Biology Area	Immunology, Stem Cells
Antigen References	<ol style="list-style-type: none">1. Rock KL, <i>et al.</i> 1989. <i>Immunol. Rev.</i> 111:195.2. Morrison SJ, <i>et al.</i> 1994. <i>Immunity</i> 1:661.3. Spangrude GJ, <i>et al.</i> 1988. <i>J. Immunol.</i> 141:3697.4. Malek T, <i>et al.</i> 1986. <i>J. Exp. Med.</i> 164:709.
Gene ID	110454


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-mouse Ly-6A/E (Sca-1), Biotin anti-mouse Ly-6A/E (Sca-1), FITC anti-mouse Ly-6A/E (Sca-1), PE anti-mouse Ly-6A/E (Sca-1), PE/Cyanine5 anti-mouse Ly-6A/E (Sca-1), Purified anti-mouse Ly-6A/E (Sca-1), PE/Cyanine7 anti-mouse Ly-6A/E (Sca-1), Alexa Fluor® 488 anti-mouse Ly-6A/E (Sca-1), Alexa Fluor® 647 anti-mouse Ly-6A/E (Sca-1), Pacific Blue™ anti-mouse Ly-6A/E (Sca-1), Brilliant Violet 421™ anti-mouse Ly-6A/E (Sca-1), PerCP anti-mouse Ly-6A/E (Sca-1), PerCP/Cyanine5.5 anti-mouse Ly-6A/E (Sca-1), APC/Cyanine7 anti-mouse Ly-6A/E (Sca-1), Brilliant Violet 510™ anti-mouse Ly-6A/E (Sca-1), Brilliant Violet 711™ anti-mouse Ly-6A/E (Sca-1), Brilliant Violet 605™ anti-mouse Ly-6A/E (Sca-1), Purified anti-mouse Ly-6A/E (Sca-1) (Maxpar® Ready), PE/Dazzle™ 594 anti-mouse Ly-6A/E (Sca-1), Brilliant Violet 785™ anti-mouse Ly-6A/E (Sca-1), Alexa Fluor® 700 anti-mouse Ly-6A/E (Sca-1), Brilliant Violet 650™ anti-mouse Ly-6A/E (Sca-1), APC/Fire™ 750 anti-mouse Ly-6A/E (Sca-1), TotalSeq™-A0130 anti-mouse Ly-6A/E (Sca-1), TotalSeq™-B0130 anti-mouse Ly-6A/E (Sca-1), TotalSeq™-C0130 anti-mouse Ly-6A/E (Sca-1), Spark NIR™ 685 anti-mouse Ly-6A/E (Sca-1), Spark PLUS UV395™ anti-mouse Ly-6A/E (Sca-1), PerCP/Fire™ 806 anti-mouse Ly-6A/E (Sca-1) Antibody, Brilliant Violet 750™ anti-mouse Ly-6A/E (Sca-1)

Product Data

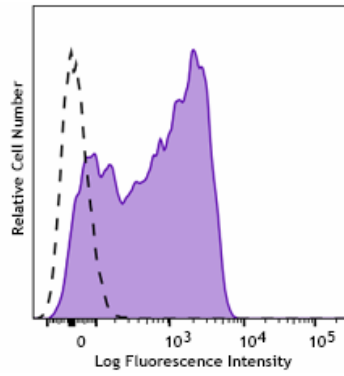

C57BL/6 mouse splenocytes were stained with anti-mouse Ly-6A/E (Sca-1) (clone D7) Brilliant Violet 750™ (filled histogram) or rat IgG2a, κ Brilliant Violet 750™ isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587