

Brilliant Violet 785™ anti-human CD66b Antibody

Catalog# / Size	392927 / 25 tests 392928 / 100 tests
Clone	6/40c
Regulatory Status	RUO
Other Names	CEACAM8, CD67, CGM6, NCA-95
Isotype	Mouse IgG1, κ
Description	CD66b is a 95-100 kD glycosylphosphatidylinositol (GPI)-linked protein also known as CD67, CGM6, and NCA-95. CD66b is a member of the immunoglobulin superfamily, carcinoembryonic antigen (CEA)-like subfamily. CD66b, expressed on granulocytes, has been reported to induce activation in neutrophils and to be involved in heterophilic adhesion with CD66c.

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Recombinant soluble human CEACAM8-Fc produced in HEK293 cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 785™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>Brilliant Violet 785™ excites at 405 nm and emits at 785 nm. The bandpass filter 780/60 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 785™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)

Application References

(PubMed link indicates BioLegend citation)

1. Singer, *et al.* 2014. *PLoS One.* 9(4):e94106. PMID: 24743304
2. Klaile, *et al.* 2013. *Respir Res.* 14:85. PMID 23941132

RRID

AB_3675124 (BioLegend Cat. No. 392927)
AB_3675124 (BioLegend Cat. No. 392928)

Antigen Details

Structure	Ig superfamily, CEA antigen group, GPI-linked glycoprotein
Distribution	Granulocytes
Function	Cell adhesion, neutrophil activation
Ligand/Receptor	CD66c
Cell Type	Granulocytes
Biology Area	Immunology
Molecular Family	Adhesion Molecules, CD Molecules

Antigen References

1. Kuijpers T, *et al.* 1993. *J Immunol.* 151:4934.
2. Kuroki M, *et al.* 1992. *J Leuk Biol.* 52:551.

Gene ID

[1088](#)


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Purified anti-human CD66b, PE anti-human CD66b, TotalSeq™-A0166 anti-human CD66b, Alexa Fluor® 594 anti-human CD66b, TotalSeq™-C0166 anti-human CD66b, Alexa Fluor® 647 anti-human CD66b, TotalSeq™-B0166 anti-human CD66b, PE/Fire™ 640 anti-human CD66b Antibody, Brilliant Violet 421™ anti-human CD66b Antibody, TotalSeq™-D0166 anti-human CD66b, TotalSeq™-Bn0166 anti-human CD66b, Brilliant Violet 605™ anti-human CD66b, Brilliant Violet 650™ anti-human CD66b, Brilliant Violet 785™ anti-human CD66b, Brilliant Violet 510™ anti-human CD66b

Product Data

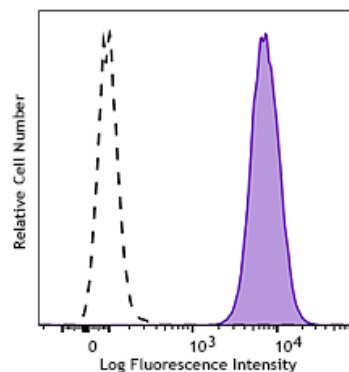

Human peripheral blood granulocytes were stained with anti-human CD66b (clone 6/40c) Brilliant Violet 785™ (filled histogram) or mouse IgG1, κ Brilliant Violet 785™ isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form,

used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587