

Pacific Blue™ anti-human CD49d Antibody

Catalog# / Size	304351 / 25 tests 304352 / 100 tests
Clone	9F10
Regulatory Status	RUO
Workshop	V S215
Other Names	VLA-4 α chain, α 4 integrin, Integrin α 4 chain, ITGA4
Isotype	Mouse IgG1, κ

Description CD49d is a 150 kD α integrin chain known as α 4 integrin or VLA-4 α chain. It forms a heterodimer with either integrin β 1 (α 4 β 1, VLA-4) or β 7 (α 4 β 7). CD49d is expressed broadly on T lymphocytes, B lymphocytes, monocytes, thymocytes, eosinophils, basophils, mast cells, NK cells, dendritic cells, and some non-hematopoietic cells, but not on normal red blood cells, platelets or neutrophils. VLA-4 binds to VCAM-1 (CD106) and fibronectin. α 4 β 7 is the receptor for VCAM-1 and MAdCAM-1. CD49d participates in mononuclear cell trafficking to endothelial sites of inflammation and has roles in cell-cell interactions and cell adhesion to extracellular matrices. CD49d is involved in lymphocyte migration, T cell activation, and hematopoietic stem cell differentiation. CD49d is a marker to isolate pure populations of Treg cells due to its absence on Foxp3⁺ cells.

Product Details

Verified Reactivity	Human, Cynomolgus, Rhesus
Reported Reactivity	African Green, Baboon, Cat, Cow, Chimpanzee, Common Marmoset, Dog, Horse, Sheep, Squirrel Monkey
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Pacific Blue™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 μ L per million cells in 100 μ L staining volume or 5 μ L per 100 μ L of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * Pacific Blue™ has a maximum emission of 455 nm when it is excited at 405 nm. Prior to using Pacific Blue™ conjugate for flow cytometric analysis, please verify your flow cytometer's capability of exciting and detecting the fluorochrome. Alexa Fluor® and Pacific Blue™ are trademarks of Life Technologies Corporation. View full statement regarding label licenses
Excitation Laser	Violet Laser (405 nm)

Application Notes	Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen tissue sections, and <i>in vitro</i> T cell costimulation ^{2,3} . The Ultra-LEAF™ Purified antibody (Endotoxin < 0.01 EU/μg, Azide-Free, 0.2 μm filtered) is recommended for functional assays (Cat. No. 304339 and 304340).
Application References	<ol style="list-style-type: none"> Schlossman S, <i>et al.</i> Eds. 1995. Leucocyte Typing V. Oxford University Press. New York. Jeong SH, <i>et al.</i> 2004. <i>J. Virol.</i> 78:6995. (Costim) Vogel TU, <i>et al.</i> 2002. <i>J. Immunol.</i> 169:4511. (Costim) Kleinewietfeld M, <i>et al.</i> 2009. <i>Blood</i> 113:827. (FC) PubMed Palacios F, <i>et al.</i> 2010. <i>Blood</i> 115:4488. PubMed Yoshino N, <i>et al.</i> 2000. <i>Exp. Anim. (Tokyo)</i> 49:97. (FC) Sestak K, <i>et al.</i> 2007. <i>Vet. Immunol. Immunopathol.</i> 119:21. Mattapallil MJ, <i>et al.</i> 2011. <i>J. Immunol.</i> 187:1977. PubMed
(PubMed link indicates BioLegend citation)	
RRID	AB_3698958 (BioLegend Cat. No. 304351) AB_3698958 (BioLegend Cat. No. 304352)

Antigen Details

Structure	Integrin, type I transmembrane glycoprotein, 150 kD.
Distribution	T cells, B cells, NK , dendritic cells, thymocytes, monocytes, eosinophils, mast cells.
Function	Lymphocyte migration, T cell activation, stem cell differentiation.
Ligand/Receptor	Fibronectin, VCAM-1, MAdCAM-1
Cell Type	B cells, Dendritic cells, Eosinophils, Mast cells, Monocytes, NK cells, T cells, Thymocytes, Tregs
Biology Area	Cell Adhesion, Cell Biology, Immunology, Innate Immunity
Molecular Family	Adhesion Molecules, CD Molecules
Antigen References	<ol style="list-style-type: none"> Elices M, Ed. 1995. <i>Springer Semin. Immunopathol.</i> 16(4). Lobb RR and Helmer ME. <i>et al.</i> 1994. <i>J. Clin. Invest.</i> 94:1722.
Gene ID	3676


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-human CD49d, PE anti-human CD49d, PE/Cyanine5 anti-human CD49d, Purified anti-human CD49d, Alexa Fluor® 594 anti-human CD49d, PerCP/Cyanine5.5 anti-human CD49d, PE/Cyanine7 anti-human CD49d, FITC anti-human CD49d, Brilliant Violet 510™ anti-human CD49d, Brilliant Violet 421™ anti-human CD49d, Purified anti-human CD49d (Maxpar® Ready), Brilliant Violet 605™ anti-human CD49d, PE/Dazzle™ 594 anti-human CD49d, APC/Cyanine7 anti-human CD49d, Brilliant Violet 711™ anti-human CD49d, Alexa Fluor® 647 anti-human CD49d, Biotin anti-human CD49d, TotalSeq™-A0576 anti-human CD49d, Brilliant Violet 785™ anti-human CD49d, APC/Fire™ 750 anti-human CD49d, Ultra-LEAF™ Purified anti-human CD49d, TotalSeq™-C0576 anti-human CD49d, TotalSeq™-B0576 anti-human CD49d Antibody, TotalSeq™-D0576 anti-human CD49d, Spark Red™ 718 anti-human CD49d (Flexi-Fluor™), Spark Blue™ 574 anti-human CD49d (Flexi-Fluor™), Spark Blue™ 550 anti-human CD49d (Flexi-Fluor™), GMP APC/Fire™ 750 anti-human CD49d, Pacific Blue™ anti-human CD49d Antibody

Product Data

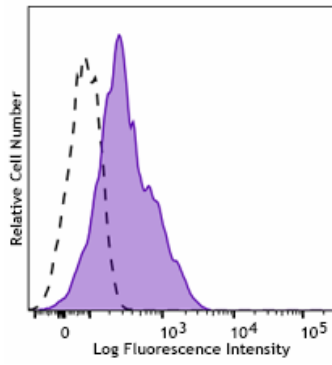

Human peripheral blood lymphocytes were stained with anti-human CD49d (clone 9F10) Pacific Blue™ (filled histogram) or mouse IgG1, κ Pacific Blue™ isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587