

PE anti-mouse Granzyme B Antibody

Catalog# / Size	168803 / 25 µg 168804 / 100 µg
Clone	S23006A
Regulatory Status	RUO
Other Names	Granzyme-2, serine protease B, CCP1, Asp-ase, CTLA-1
Isotype	Rat IgG1, λ
Description	Granzyme B, also known as granzyme-2, serine protease B, CCP1, Asp-ase, and CTLA-1, is a 32 kD serine protease. Granzyme B is abundantly stored in the granules of cytotoxic T lymphocytes and NK cells. Low level of expression has been reported in granulocytes, B cells, and activated dendritic cells. Granzyme B is crucial for rapid induction of cell death and apoptosis through interaction with mannose-6-phosphate receptor.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Rat
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Preparation	The antibody was purified by affinity chromatography and conjugated with PE under optimal conditions.
Concentration	0.2 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	ICFC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by intracellular immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is ≤ 0.25 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application.
Excitation Laser	Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm)
Application Notes	Clones S23006A and GB11 partially block each other.

Antigen Details

Structure	32 kD serine protease
Distribution	Cytotoxic T-cells and NK cells, low on granulocytes, B cells and activated dendritic cells
Function	Induction of cell death and apoptosis
Ligand/Receptor	Mannose-6-phosphate receptor
Cell Type	Neutrophils, T cells
Biology Area	Cell Biology, Immunology, Innate Immunity, Neuroscience

Molecular Family

Enzymes and Regulators, Proteases

Antigen References

1. Estebanez-Perpina E, *et al.* 2000. *Biol Chem.* 381:1203.
2. Griffiths GM. and S. Isaaz, *et al.* 1993. *J. Cell Biol.* 120:885.
3. Spaeny-Dekking EH, *et al.* 1998. *J. Immunol.* 160:3610.
4. Wagner C, *et al.* 2008. *Mol. Immunol.* 45:1761.

Gene ID

[14939](#)


Related Protocols

- [Intracellular Flow Cytometry Staining Protocol](#)

Other Formats

Purified anti-mouse Granzyme B, PE anti-mouse Granzyme B, Alexa Fluor® 647 anti-mouse Granzyme B

Product Data

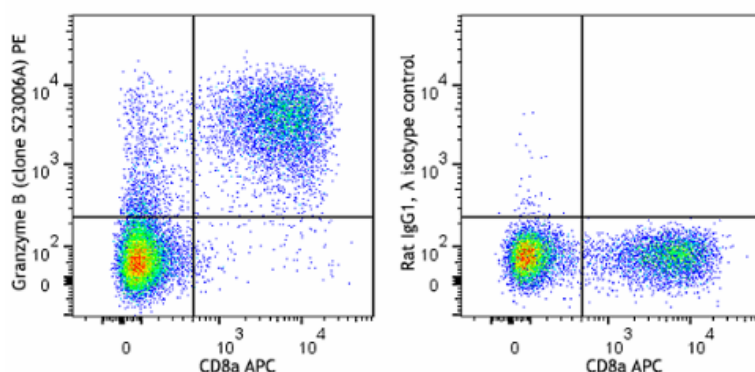

BALB/c mouse splenocytes were stimulated with plate-bound anti-mouse TCR β (Cat. No. 109254) for 2 days, and surface stained with anti-mouse CD8a (clone 53-6.7) APC. Cells were then fixed and permeabilized with Cyto-Fast™ Fix/Perm Buffer Set (Cat. No. 426803) and intracellularly stained with anti-mouse Granzyme B (clone S23006A) PE (left) or rat IgG1, λ PE isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587