

PE anti-Rb Phospho (Ser795) Antibody

Catalog# / Size	640455 / 25 tests 640456 / 100 tests
Clone	A24001E
Regulatory Status	RUO
Other Names	RB1, Retinoblastoma-associated protein, p105-Rb, p110-RB1
Isotype	Mouse IgG1, κ
Description	The Retinoblastoma protein (Rb) is a tumor suppressor that regulates cell cycle progression from G1 phase into the S phase. This regulation is crucial for preventing uncontrolled cell proliferation, which can lead to cancer. Hence, mutations or dysfunctions in the RB1 gene are associated with several cancers including retinoblastoma, lung, bladder, and breast cancers. Phosphorylation at serine 795 of the Rb protein plays a crucial role in regulating its function. The serine 795 phosphorylation leads to Rb dissociation from E2F transcription factor and allowing E2F to induce the transcription of cell cycle genes. This regulatory mechanism is vital for normal cell function and its dysregulation can contribute to cancer development.

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Synthetic peptide of human Rb phosphorylated at serine 795
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with PE under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	ICFC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by intracellular immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 μ L per million cells in 100 μ L staining volume or 5 μ L per 100 μ L of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.
Excitation Laser	Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm)
Additional Product Notes	For intracellular flow cytometry (ICFC) staining fixation and permeabilization we recommend use True-Phos™ Perm Buffer Set (Cat. No. 425401). We do not recommend fixation permeabilization using neither the Cyto-Fast™ Fix/Perm buffer Set (Cat. No. 426803) nor the True-Nuclear™ Transcription Factor Buffer Set (Cat. No. 424401).
RRID	AB_3720507 (BioLegend Cat. No. 640455) AB_3720507 (BioLegend Cat. No. 640456)

Antigen Details

Structure	The Rb protein is 928 amino acids in length and has a molecular weight of approximately 106 kDa.
Distribution	Nucleus, cytoplasmic during keratinocyte differentiation
Function	• Tumor Suppression: Rb acts as a tumor suppressor by inhibiting excessive cell growth and proliferation. • Cell Cycle Regulation: Rb controls progression from the G1 phase to the S phase of the cell cycle by binding to and inhibiting E2F transcription factors. • Transcription Regulation: Rb regulates the transcription of genes necessary for DNA synthesis and cell cycle progression. • Apoptosis: Rb can also induce apoptosis under certain conditions, contributing to its role in preventing cancer.
Interaction	• E2F Transcription Factors: Rb binds to E2F transcription factors, inhibiting their activity and preventing the transcription of genes required for cell cycle progression. • DP Proteins: Rb also interacts with DP proteins, which form heterodimers with E2F, enhancing the regulation of E2F activity. • Cyclin-Dependent Kinases (CDKs): CDKs phosphorylate Rb, leading to its inactivation and allowing the cell cycle to progress. • Histone Deacetylases (HDACs): Rb recruits HDACs to E2F target genes, leading to chromatin remodeling and transcriptional repression. • Chromatin Remodeling Enzymes: Rb interacts with various chromatin remodeling enzymes, such as methylases and acetylases, to regulate gene expression.
Biology Area	Apoptosis/Tumor Suppressors/Cell Death, Cancer Biomarkers, Cell Biology, Cell Cycle/DNA Replication, Cell Death, Cell Proliferation and Viability, DNA Repair/Replication, Signal Transduction, Transcription Factors
Molecular Family	Phospho-Proteins, Tumor Suppressors
Antigen References	1. Rezaeian A, <i>et al.</i> 2023. <i>Adv Protein Chem Struct Biol.</i> 135:179. 2. Janostiak R, <i>et al.</i> 2022. <i>Cancers (Basel).</i> 14:1265. 3. Indovina P, <i>et al.</i> 2015. <i>Oncotarget.</i> 6:17873. 4. Linn P, <i>et al.</i> 2021. <i>Cancers (Basel).</i> 13:3737.
Regulation	• Phosphorylation: Rb is regulated by phosphorylation, particularly by cyclin-dependent kinases (CDKs). Phosphorylation at specific sites, such as serine 795, inactivates Rb, allowing cell cycle progression. • Dephosphorylation: Dephosphorylation of Rb reactivates its tumor suppressor functions, inhibiting cell cycle progression. • Gene Mutations: Mutations in the RB1 gene can lead to loss of Rb function, contributing to the development of cancers such as retinoblastoma, lung, bladder, and breast cancers.
Gene ID	5925


Related Protocols

- [Intracellular Flow Cytometry Staining Protocol](#)

Other Formats

Purified anti-Rb Phospho (Ser795), Alexa Fluor® 647 anti-Rb Phospho (Ser795), PE anti-Rb Phospho (Ser795)

Product Data

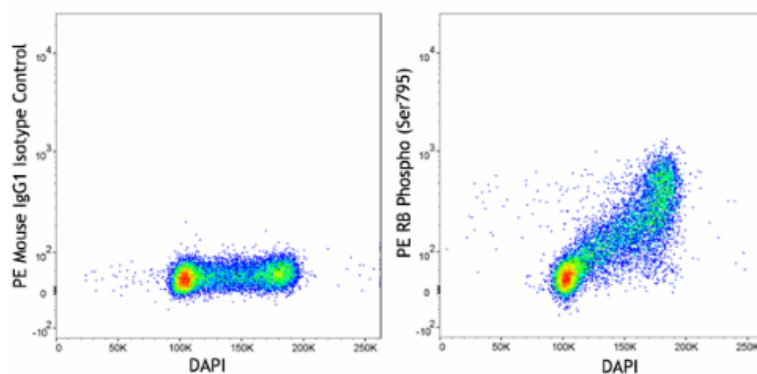

HeLa cells treated with Nocodazole to arrest cell cycle at G2/M and induce phosphorylation, were fixed and permeabilized using the True-Phos™ Perm Buffer Set (Cat. No. 425401) and intracellularly stained with PE anti-RB Phospho (Ser795) (clone A24001E) (right panel) or PE Mouse IgG1, κ Isotype Control (Cat. No. 400139) (left panel). Cell cycle analysis was performed by DAPI staining (Cat. No. 422801).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587