

PE anti-SOX17 Recombinant Antibody

Catalog# / Size	628163 / 25 tests 628164 / 100 tests
Clone	O92B6.Rec
Regulatory Status	RUO
Other Names	Transcription factor SOX17, Transcription factor SOX 17, Transcription factor SOX-17, VUR3, SRY-Box17, SRY-Related HMG-Box Transcription Factor 17
Isotype	Mouse IgG1, κ
Description	SOX17 is an important member of the SOX family of transcription factors, which are defined by possession of a homologous sequence called the HMG (high mobility group) box. SOX genes are primarily involved in developmental processes including germ layer formation, cell-type specification, and organogenesis. SOX17 activates endodermal target genes, and is required for the formation of endoderm and blood vessel endothelium. SOX17 also plays critical roles in the maintenance of fetal and neonatal hematopoietic stem cells. Additionally, SOX17 is a regulator of human primordial germ cell fate. SOX17 also acts as a tumor suppressor through suppression of Wnt signaling, and a reduction of SOX17 expression is found in gastrointestinal and colon carcinomas. Additionally, mutations in SOX17 can induce stem-like characteristics in tumor cells by causing SOX17 to dimerize with OCT4 and induce pluripotency transcriptional programs.

Product Details

Verified Reactivity	Human
Antibody Type	Recombinant
Host Species	Mouse
Immunogen	Recombinant protein expressed in E. Coli corresponding to residues 177-414 of human SOX17
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with PE under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool).
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	ICFC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by intracellular flow cytometry using our True-Phos™ Perm Buffer in Cell Suspensions Protocol. For flow cytometric staining, the suggested use of this reagent is ≤ 0.5 μ g per million cells in 100 μ L volume. It is recommended that the reagent be titrated for optimal performance for each application.
Excitation Laser	Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm)
Application Notes	This antibody (clone O92B6.Rec) does not cross-react with mouse SOX17 in western blot (WB).
RRID	AB_3720591 (BioLegend Cat. No. 628163) AB_3720591 (BioLegend Cat. No. 628164)

Antigen Details

Structure	SOX17 is a 414 amino acid protein with a predicted molecular weight of 44.1 kD.
Distribution	Nucleus
Function	Transcription factor
Interaction	Interacts with CTNNB1, LEF1, TCF4
Biology Area	Cell Biology, Signal Transduction, Transcription Factors

Antigen References

1. Wegner M. 1999. *Nucleic Acids Res.* 27:1409-20.
2. Du YC, *et al.* 2009. *Gastroenterology.* 137:1346-57.
3. Zorn AM, *et al.* 1999. *Mol Cell.* 4:487-98.
4. Sinner D, *et al.* 2004. *Development.* 131:3069-80.
5. Hudson C, *et al.* 1997. *Cell.* 91:397-405.
6. Matsui T, *et al.* 2006. *J Cell Sci.* 119:3513-26.
7. Kim I, *et al.* 2007. *Cell.* 130:470-83.
8. Zhang W, *et al.* 2008. *Cancer Res.* 68:2764-72.
9. Tan DS, *et al.* 2020. *Semin Cancer Biol.* 67:65-73.
10. Irie N, *et al.* 2015. *Cell.* 160:253-68.

Gene ID [64321](#)


Related Protocols

- [Intracellular Staining With True-Phos™ Perm Buffer in Cell Suspensions Protocol](#)
- [Intracellular Staining With True-Phos™ Perm Buffer in Whole Blood](#)
- [Intracellular Flow Cytometry Staining Protocol](#)

Other Formats

Purified anti-SOX17 Recombinant Antibody, Alexa Fluor® 647 anti-SOX17 Recombinant Antibody (For Flow Cytometry), Alexa Fluor® 488 anti-SOX17 Recombinant Antibody, PE anti-SOX17 Recombinant Antibody, Alexa Fluor® 647 anti-SOX17 Recombinant Antibody (for IHC and ICC)

Product Data

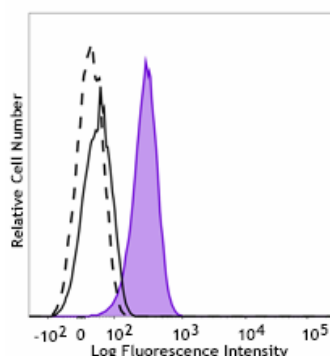

OVCAR3 cells (positive control, filled histogram) and HeLa cells (negative control, open histogram) were fixed and permeabilized using the True-Phos™ Perm Buffer (Cat. No. 425401) and intracellularly stained with PE anti-SOX17 Recombinant (clone O92B6.Rec), or PE mouse IgG1, κ isotype control (dashed histogram) (Cat No. 400112).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended

use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com

Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587