

## PE/Fire™ 700 anti-human CD86 Antibody

<b>Catalog# / Size</b>	305463 / 25 tests 305464 / 100 tests
<b>Clone</b>	IT2.2
<b>Regulatory Status</b>	RUO
<b>Workshop</b>	VI CD86.8
<b>Other Names</b>	B7-2, B70
<b>Isotype</b>	Mouse IgG2b, κ
<b>Description</b>	CD86 is an 80 kD immunoglobulin superfamily member also known as B7-2, B70, and Ly-58. CD86 is expressed on activated B and T cells, monocytes/macrophages, dendritic cells, and astrocytes. CD86, along with CD80, is the ligand of CD28 and CD152 (CTLA-4). CD86 is expressed earlier in the immune response than CD80. CD86 has also been shown to be involved in immunoglobulin class-switching and triggering of NK cell-mediated cytotoxicity. CD86 binds to CD28 to transduce costimulatory signals for T cell activation, proliferation, and cytokine production. CD86 can bind to CD152 as well, also known as CTLA-4, to deliver an inhibitory signal to T cells.

### Product Details

---

<b>Verified Reactivity</b>	Human, Cynomolgus, Rhesus
<b>Reported Reactivity</b>	African Green, Baboon, Capuchin Monkey, Common Marmoset, Cotton-topped Tamarin, Chimpanzee
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with PE/Fire™ 700 under optimal conditions.
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our <a href="#">Certificate of Analysis</a> online tool.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.  * PE/Fire™ 700 has a maximum excitation of 565 nm and a maximum emission of 695 nm.
<b>Excitation Laser</b>	Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm)
<b>Application Notes</b>	Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen tissue sections <sup>6</sup> , Western blotting <sup>3</sup> , and blocking of T cell activation <sup>2,4,5</sup> . The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for functional assays (Cat. Nos. 305449 & 305450).

## Application References

(PubMed link indicates BioLegend citation)

1. Kishimoto T, *et al.* Eds. 1997. Leucocyte Typing VI. Garland Publishing Inc. London.
2. Dieu M. 1998. *J. Exp. Med.* 188:373. (Block)
3. Esser M, *et al.* 2001. *J. Virol.* 75:6173. (WB)
4. Jeannin P, *et al.* 1999. *J. Immunol.* 162:2044. (Block)
5. Kapsogeorgou EK, *et al.* 2001. *J. Immunol.* 166:3107. (Block)
6. Geissmann F, *et al.* 2001. *Blood* 97:1241. (IHC)

RRID AB\_3662214 (BioLegend Cat. No. 305463)  
AB\_3662214 (BioLegend Cat. No. 305464)

## Antigen Details

<b>Structure</b>	Ig superfamily, single-chain transmembrane glycoprotein, 80 kD
<b>Distribution</b>	Monocytes/macrophages, activated B cells and T cells, dendritic cells
<b>Function</b>	T cell costimulation
<b>Ligand/Receptor</b>	CD28, CD152
<b>Cell Type</b>	B cells, Dendritic cells, Macrophages, Monocytes, T cells, Tregs
<b>Biology Area</b>	Cell Biology, Costimulatory Molecules, Immunology, Neuroscience, Neuroscience Cell Markers
<b>Molecular Family</b>	CD Molecules, Immune Checkpoint Receptors
<b>Antigen References</b>	1. Hathcock K, <i>et al.</i> 1996. <i>Adv. Immunol.</i> 62:131. 2. June C, <i>et al.</i> 1994. <i>Immunol. Today</i> 15:321.
<b>Gene ID</b>	<a href="#">942</a>


## Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

## Other Formats

Biotin anti-human CD86, PE anti-human CD86, PE/Cyanine5 anti-human CD86, Purified anti-human CD86, APC anti-human CD86, Alexa Fluor® 488 anti-human CD86, Alexa Fluor® 647 anti-human CD86, Pacific Blue™ anti-human CD86, PE/Cyanine7 anti-human CD86, PerCP/Cyanine5.5 anti-human CD86, Brilliant Violet 421™ anti-human CD86, Brilliant Violet 650™ anti-human CD86, Brilliant Violet 605™ anti-human CD86, Brilliant Violet 510™ anti-human CD86, PE/Dazzle™ 594 anti-human CD86, Purified anti-human CD86 (Maxpar® Ready), Brilliant Violet 711™ anti-human CD86, Brilliant Violet 785™ anti-human CD86, TotalSeq™-A0006 anti-human CD86, TotalSeq™-B0006 anti-human CD86, TotalSeq™-C0006 anti-human CD86, Ultra-LEAF™ Purified anti-human CD86, KIRAVIA Blue 520™ anti-human CD86, TotalSeq™-D0006 anti-human CD86, PE/Fire™ 810 anti-human CD86, PE/Fire™ 640 anti-human CD86, PerCP/Fire™ 806 anti-human CD86, PerCP/Fire™ 780 anti-human CD86, PE/Fire™ 700 anti-human CD86

## Product Data

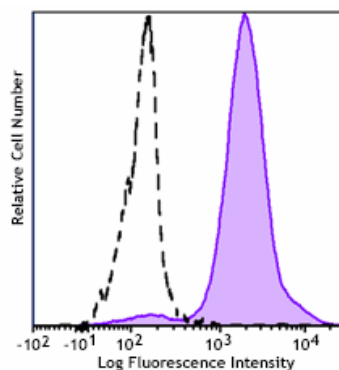

Human peripheral blood monocytes were stained with anti-human CD86 (clone IT2.2) PE/Fire™ 700 (filled histogram) or mouse IgG2b, κ PE/Fire™ 700 isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at [www.biolegend.com/terms](http://www.biolegend.com/terms) ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587