

PE/Fire™ 744 anti-human CD14 Antibody

Catalog# / Size	301877 / 25 tests 301878 / 100 tests
Clone	M5E2
Regulatory Status	RUO
Workshop	III 329
Other Names	LPS receptor
Isotype	Mouse IgG2a, κ
Description	CD14 is a 53-55 kD glycosylphosphatidylinositol (GPI)-linked membrane glycoprotein also known as LPS receptor. CD14 is expressed at high levels on monocytes and macrophages, and at lower levels on granulocytes. Some dendritic cell populations such as interfollicular dendritic cells, reticular dendritic cells, and Langerhans cells have also been reported to express CD14. As a high-affinity receptor for LPS, CD14 is involved in the clearance of gram-negative pathogens, and in the upregulation of adhesion molecules and expression of cytokines in monocytes and neutrophils.

Product Details

Verified Reactivity	Human, Cynomolgus, Rhesus
Reported Reactivity	African Green, Capuchin Monkey, Cow, Chimpanzee, Common Marmoset, Cotton-topped Tamarin, Dog, Pigtailed Macaque, Squirrel Monkey
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Full-length human CD14 protein
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with PE/Fire™ 744 under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. True-Stain Monocyte Blocker™ (Cat No. 426102) is recommended to minimize non-specific staining of PE/Fire™ 744 on monocytes and macrophages. * PE/Fire™ 744 has a maximum excitation of 565 nm and a maximum emission of 744 nm.
Excitation Laser	Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm)
Application Notes	The M5E2 antibody inhibits monocyte activation and cytokine production induced by LPS. Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen sections, blocking of LPS stimulation ⁴ , and immunofluorescence

microscopy⁵. Clone M5E2 is not recommended for immunohistochemical staining of formalin-fixed paraffin-embedded sections. The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/μg, Azide-Free, 0.2 μm filtered) is recommended for functional assays (Cat. No. 301861 and 301862).

Application References

(PubMed link indicates
BioLegend citation)

1. McMichael A, *et al.* 1987. Leucocyte Typing III. Oxford University Press. New York.
2. Knapp W, *et al.* Eds. 1989. Leucocyte Typing IV. Oxford University Press. New York. (IHC-F)
3. Schlossman S, *et al.* Eds. 1995. Leucocyte Typing V. Oxford University Press. New York.
4. Power CP, *et al.* 2004. *J. Immunol.* 173:5229. (Block)
5. Williams KC, *et al.* 2001. *J. Exp. Med.* 193:905.
6. Iwamoto S, *et al.* 2007. *J. Immunol.* 179:1449. (FC) [PubMed](#)
7. Santer DM, *et al.* 2010. *J. Immunol.* 485:4739. [PubMed](#)
8. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
9. Zizzo G, *et al.* 2012. *J. Immunol.* 189:3508. [PubMed](#)
10. Stoeckius M, *et al.* 2017. *Nat. Methods.* 14:865. (PG)
11. Peterson VM, *et al.* 2017. *Nat. Biotechnol.* 35:936. (PG)

RRID

AB_3720499 (BioLegend Cat. No. 301877)

AB_3720499 (BioLegend Cat. No. 301878)

Antigen Details

Structure	GPI-linked membrane glycoprotein, 53-55 kD
Distribution	Monocytes, macrophages, granulocytes (low)
Function	LPS receptor, clearance of Gram-negative pathogens
Ligand/Receptor	LPS
Cell Type	Granulocytes, Macrophages, Monocytes, Neutrophils
Biology Area	Cell Biology, Immunology, Innate Immunity, Neuroinflammation, Neuroscience
Molecular Family	CD Molecules
Antigen References	1. Stocks S, <i>et al.</i> 1990. <i>Biochem. J.</i> 268:275. 2. Wright S, <i>et al.</i> 1990. <i>Science</i> 249:1434.
Gene ID	929


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

TotalSeq™-D0081 anti-human CD14, APC/Fire™ 750 anti-human CD14, GMP FITC anti-human CD14, GMP APC anti-human CD14, GMP Pacific Blue™ anti-human CD14, GMP PE anti-human CD14, GMP APC/Fire™ 750 anti-human CD14, PerCP/Cyanine5.5 anti-human CD14, GMP PE/Dazzle™ 594 anti-human CD14, APC/Fire™ 810 anti-human CD14, PE/Fire™ 700 anti-human CD14, PerCP anti-human CD14, PE/Cyanine7 anti-human CD14, FITC anti-human CD14, PE/Dazzle™ 594 anti-human CD14, APC anti-human CD14, Brilliant Violet 605™ anti-human CD14, PE anti-human CD14, Purified anti-human CD14, Brilliant Violet 510™ anti-human CD14, Purified anti-human CD14 (Maxpar® Ready), Ultra-LEAF™ Purified anti-human CD14, Spark Violet™ 538 anti-human CD14, Spark Violet™ 500 anti-human CD14, Alexa Fluor® 488 anti-human CD14, Pacific Blue™ anti-human CD14, APC/Cyanine7 anti-human CD14, Biotin anti-human CD14, TotalSeq™-A0081 anti-human CD14, PE/Cyanine5 anti-human CD14, Brilliant Violet 711™ anti-human CD14, Brilliant Violet 785™ anti-human CD14, PE/Fire™ 744 anti-human CD14, Alexa Fluor® 647 anti-human CD14, Brilliant Violet 421™ anti-human CD14, Brilliant Violet 570™ anti-human CD14, Brilliant Violet 650™ anti-human CD14, GMP PE/Cyanine7 anti-human CD14, Spark Red™ 718 anti-human CD14, GMP PerCP/Cyanine5.5 anti-human CD14, Alexa Fluor® 700 anti-human CD14, TotalSeq™-B0081 anti-human CD14, TotalSeq™-C0081 anti-human CD14

Product Data

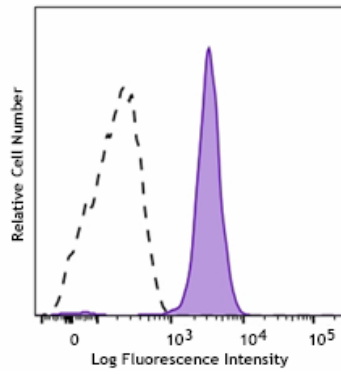

Human peripheral blood monocytes were stained with anti-human CD14 (clone M5E2) PE/Fire™ 744 (filled histogram) or mouse IgG2a, κ PE/Fire™ 744 isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587