

PE/Fire™ 744 anti-human CD45 Antibody

Catalog# / Size	304040 / 25 tests 304046 / 100 tests
Clone	HI30
Regulatory Status	RUO
Workshop	IV N816
Other Names	LCA, T200
Isotype	Mouse IgG1, κ
Description	CD45 is a 180-240 kD single chain type I membrane glycoprotein also known as leukocyte common antigen (LCA) and T200. It is a tyrosine phosphatase expressed on the plasma membrane of all hematopoietic cells, except erythrocytes and platelets. CD45 is a signaling molecule that regulates a variety of cellular processes including cell growth, differentiation, cell cycle, and oncogenic transformation. CD45 plays a critical role in T and B cell antigen receptor-mediated activation by dephosphorylating substrates including p56Lck, p59Fyn, and other Src family kinases. CD45 non-covalently associates with lymphocyte phosphatase-associated phosphoprotein (LPAP) on T and B lymphocytes. CD45 has been reported to bind galectin-1 and to be associated with several other cell surface antigens including CD1, CD2, CD3, and CD4.

Product Details

Verified Reactivity	Human
Reported Reactivity	Chimpanzee
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with PE/Fire™ 744 under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. True-Stain Monocyte Blocker™ (Cat No. 426102) is recommended to minimize non-specific staining of PE/Fire™ 744 on monocytes and macrophages. * PE/Fire™ 744 has a maximum excitation of 565 nm and a maximum emission of 744 nm.
Excitation Laser	Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm)
Application Notes	Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen tissue sections and formalin-fixed paraffin-embedded tissue sections ⁹ , inhibition of CD45 functions ⁴ , immunofluorescence ¹¹ , Western blotting ³ , and spatial biology

(IBEX)^{16,17}.

It was found that the HI30 clone and the 2D1 clone can cross block each other's binding.

Application References

(PubMed link indicates
BioLegend citation)

1. Knapp W, *et al.* 1989. Leucocyte Typing IV. Oxford University Press. New York.
2. Kishihara K, *et al.* 1993. *Cell* 74:143.
3. Esser M, *et al.* 2001. *J. Virol.* 75:6173. (WB)
4. Yamada T, *et al.* 2002. *J. Biol. Chem.* 277:28830.
5. Nagano M, *et al.* 2007. *Blood* 110:151.
6. Jiang Q, *et al.* 2008. *Blood* 112:2858. [PubMed](#)
7. Morozov A, *et al.* 2010. *Clin Cancer Res.* 16:5630. [PubMed](#)
8. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
9. Friedman T, *et al.* 1999. *J. Immunol.* 162:5256. (IHC)
10. Oeztuerk-Winder F, *et al.* 2012. *EMBO J.* 31:3431. (FC) [PubMed](#)
11. Rees LE, *et al.* 2003. *Clin. Exp. Immunol.* 134:497. (IF)
12. Lee J, *et al.* 2015. *J Exp Med.* 212:385. [PubMed](#)

[See More](#)

RRID

AB_3720497 (BioLegend Cat. No. 304040)

AB_3720497 (BioLegend Cat. No. 304046)

Antigen Details

Structure	Tyrosine phosphatases, type I transmembrane protein, 180-240 kD (multiple isoforms)
Distribution	Hematopoietic cells, not expressed in circulating erythrocytes or platelets
Function	TCR and BCR mediated activation
Ligand/Receptor	Galectin-1, CD2, CD3, CD4
Cell Type	Hematopoietic stem and progenitors, Mesenchymal Stem Cells
Biology Area	Cell Biology, Immunology, Inhibitory Molecules, Innate Immunity, Neuroscience, Neuroscience Cell Markers, Stem Cells
Molecular Family	CD Molecules
Antigen References	1. Thomas M. 1989. <i>Annu. Rev. Immunol.</i> 7:339. 2. Trowbridge I, <i>et al.</i> 1994. <i>Annu. Rev. Immunol.</i> 12:85.
Gene ID	5788

Related Protocols


- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Alexa Fluor® 660 anti-human CD45 Antibody, Spark Violet™ 538 anti-human CD45, GMP APC/Fire™ 750 anti-human CD45, GMP PerCP anti-human CD45, Spark UV™ 387 anti-human CD45, GMP FITC anti-human CD45, GMP Pacific Blue™ anti-human CD45, GMP PE/Dazzle™ 594 anti-human CD45, GMP PE/Cyanine7 anti-human CD45, TotalSeq™-D0391 anti-human CD45, Spark YG™ 593 anti-human CD45, Spark Violet™ 500 anti-human CD45, PE/Fire™ 810 anti-human CD45, Spark NIR™ 685 anti-human CD45, GMP PE anti-human CD45, GMP APC anti-human CD45, Alexa Fluor® 594 anti-human CD45, Pacific Blue™ anti-human CD45, PE/Cyanine7 anti-human CD45, PE/Dazzle™ 594 anti-human CD45, PE/Cyanine5 anti-human CD45, APC anti-human CD45, Biotin anti-human CD45, FITC anti-human CD45, PE anti-human CD45, Brilliant Violet 421™ anti-human CD45, Brilliant Violet 605™ anti-human CD45, Purified anti-human CD45, Brilliant Violet 650™ anti-human CD45, GMP PerCP/Cyanine5.5 anti-human CD45, Brilliant Violet 785™ anti-human CD45, Spark Blue™ 515 anti-human CD45, PE/Fire™ 744 anti-human CD45, Spark PLUS UV395™ anti-human CD45, Alexa Fluor® 488 anti-human CD45, Spark PLUS B550™ anti-human CD45 Antibody, Spark PLUS V475™ anti-human CD45 Antibody, Alexa Fluor® 700 anti-human CD45, PerCP anti-human CD45, TotalSeq™-C0391 anti-human CD45, APC/Fire™ 750 anti-human CD45, PerCP/Cyanine5.5 anti-human CD45, PE/Fire™ 640 anti-human CD45, APC/Fire™ 810 anti-human CD45, Brilliant Violet 510™ anti-human CD45, Purified anti-human CD45 (Maxpar® Ready), Brilliant Violet 711™ anti-human CD45, GMP Spark Violet™ 538 anti-human CD45, Alexa Fluor® 647 anti-human CD45, Ultra-LEAF™ Purified anti-human CD45 (Flexi-Pure™), APC/Cyanine7 anti-human CD45, Spark YG™ 570 anti-human CD45, PE/Fire™ 700 anti-human CD45,

Brilliant Violet 570™ anti-human CD45, Spark Violet™ 423 anti-human CD45, Spark YG™ 581 anti-human CD45 (Flexi-Fluor™), TotalSeq™-A0391 anti-human CD45, TotalSeq™-B0391 anti-human CD45

Product Data

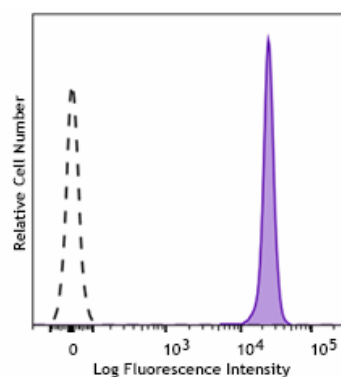

Human peripheral blood lymphocytes were stained with anti-human CD45 (clone HI30) PE/Fire™ 744 (filled histogram) or mouse IgG1, κ PE/Fire™ 744 isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587