

PerCP/Cyanine5.5 anti-mouse CX3CR1 Recombinant Antibody

Catalog# / Size	153709 / 25 µg 153710 / 100 µg
Clone	QA16A03
Regulatory Status	RUO
Workshop	HLDA11
Other Names	Chemokine (C-X3-C motif) receptor 1, Fractalkine receptor, GPR13, CCRL1, CMKBRL1, CMKDR1, V28, CD373
Isotype	Mouse IgG1, κ
Description	CX3CR1 is a 40 kD, G-protein coupled receptor, with seven transmembrane regions. CX3CR1 is expressed by resident and alternatively activated macrophages (M2), a subset of monocytes, dendritic cells (DCs), NK cells, a subset of memory T cells, and mast cells. CX3CR1 is involved in cell recruitment during inflammation and participates in cell adhesion and extravasation from blood vessels. Its ligand is CX3CL1, also known as fractalkine or neurotactin. CX3CR1 is also a coreceptor for HIV1 and variations in this gene leads to increased susceptibility to HIV. In the brain, it is expressed by glial cells, which interact with CX3CL1 expressed by neurons.

Product Details

Verified Reactivity	Mouse
Antibody Type	Recombinant
Host Species	Mouse
Immunogen	Mouse CX3CR1-transfected cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Preparation	The antibody was purified by affinity chromatography and conjugated with PerCP/Cyanine5.5 under optimal conditions.
Concentration	0.2 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is ≤ 0.25 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application. * PerCP/Cyanine5.5 has a maximum absorption of 482 nm and a maximum emission of 690 nm.
Excitation Laser	Blue Laser (488 nm)
Application Notes	This recombinant antibody clone is derived from clone SA011F11.

Antigen Details

Structure	G-protein coupled receptor, 7 transmembrane regions, 40 kD.
------------------	---

Distribution	Resident and alternatively activated macrophages (M2), subset of monocytes, dendritic cells (DCs), natural killer cells (NK), subset of memory T cells, mast cells, glial cells, and astrocytes.
Function	Involved in cell recruitment during inflammation, participates in cell adhesion and extravasation from blood vessels, coreceptor for HIV1.
Ligand/Receptor	CX3CL1 (fractalkine, neurotactin).
Cell Type	Astrocytes, Dendritic cells, Macrophages, Mast cells, Microglia, Monocytes, NK cells, T cells
Biology Area	Cell Biology, Immunology, Neuroinflammation, Neuroscience, Neuroscience Cell Markers
Molecular Family	Cytokine/Chemokine Receptors, GPCR
Antigen References	1. Hamon P, <i>et al.</i> 2017. <i>Blood</i>.129:1296 2. Afik R, <i>et al.</i> 2016. <i>J Exp Med.</i> 213:2315 3. Mass E, <i>et al.</i> 2016. <i>Science.</i> 353 pii: aaf4238 4. Safaiyan S, <i>et al.</i> 2016. <i>Nat Neurosci.</i> 19:995 5. Peng J, <i>et al.</i> 2016. <i>Nat Commun.</i> 7:12029
Gene ID	13051


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

PE anti-mouse CX3CR1 Recombinant Antibody, PE/Cyanine7 anti-mouse CX3CR1 Recombinant Antibody, PerCP/Cyanine5.5 anti-mouse CX3CR1 Recombinant Antibody, APC anti-mouse CX3CR1 Recombinant Antibody

Product Data

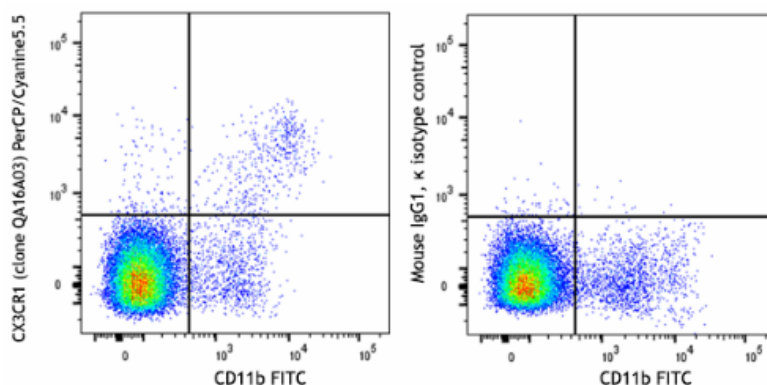

C57BL/6 mouse splenocytes was stained with anti-mouse/human CD11b (clone M1/70) FITC and anti-mouse CX3CR1 recombinant (clone QA16A03) PerCP/Cyanine5.5 (left) or mouse IgG1, κ PerCP/Cyanine5.5 isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587