

## PerCP/Fire™ 780 anti-human CD197 (CCR7) Antibody

<b>Catalog# / Size</b>	353283 / 25 tests 353284 / 100 tests
<b>Clone</b>	G043H7
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	BLR2, CDw197, EBI1, CMKBR7
<b>Isotype</b>	Mouse IgG2a, κ
<b>Description</b>	CCR7, also known as CD197, is a chemokine receptor that binds CCL19 and CCL21. CCR7 and its ligands link innate and adaptive immunity by affecting interactions between T cells and dendritic cells and their downstream effect. Naïve T cells enter the lymph node through high endothelial venules, which express CCL21. Dendritic cells and macrophages enter the lymph node through afferent lymphatics. The encounter of T cells and dendritic cells in the T cell zone is CCR7-dependent. In addition, during immunological surveillance, B cells recirculate between B-cell-rich compartments (follicles or B cell zones) in secondary lymphoid organs, surveying for antigen. After antigen binding, B cells move to the boundary of B and T zones to interact with T-helper cells; this B cell migration is directed by CCR7 and its ligands. CCR7-positive cancer cell expression has been associated with lymph node metastasis.

### Product Details

---

<b>Verified Reactivity</b>	Human
<b>Reported Reactivity</b>	African Green, Baboon, Cynomolgus, Rhesus
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Immunogen</b>	CCR7-transfected cells
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with PerCP/Fire™ 780 under optimal conditions.
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our <a href="#">Certificate of Analysis</a> online tool.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.  * PerCP/Fire™ 780 has a maximum excitation of 478 nm and a maximum emission of 780 nm.
<b>Excitation Laser</b>	Blue Laser (488 nm)
<b>RRID</b>	AB_3675074 (BioLegend Cat. No. 353283) AB_3675074 (BioLegend Cat. No. 353284)

### Antigen Details

---

<b>Structure</b>	Chemokine receptor, G protein-coupled receptors (GPCR), seven transmembrane receptor.
<b>Distribution</b>	T cells, B cells, NK, dendritic cells.
<b>Function</b>	The chemokine receptor CCR7 plays a pivotal role in the homing of naïve T cells and regulatory T cells to secondary lymphoid organs, and the migration of dendritic cells into afferent lymphatic vessels.
<b>Ligand/Receptor</b>	CCL19 and CCL21.
<b>Cell Type</b>	B cells, Dendritic cells, NK cells, T cells
<b>Biology Area</b>	Immunology
<b>Molecular Family</b>	CD Molecules, Cytokine/Chemokine Receptors, GPCR
<b>Antigen References</b>	<ol style="list-style-type: none"> <li>1. Yanagihara S, <i>et al.</i> 1998. <i>J. Immunol.</i> 161:3096.</li> <li>2. Charo IF, <i>et al.</i> 2006. <i>N. Engl. J. Med.</i> 354:610.</li> <li>3. Reif K, <i>et al.</i> 2002. <i>Nature</i> 416:94.</li> <li>4. Nakata B, <i>et al.</i> 2008. <i>Oncology</i> 74:69.</li> <li>5. Brodie T. <i>et al.</i> 2013. <i>Cytometry A.</i> 6: 530-2. <a href="#">PubMed</a></li> <li>6. Graves A.J. <i>et al.</i> 2014. <i>Cytometry A.</i> 7: 576–9 <a href="#">PubMed</a></li> <li>7. Moncunill G. <i>et al.</i> 2014. <i>Cytometry A.</i> 12: 995-8 <a href="#">PubMed</a></li> </ol>
<b>Gene ID</b>	<a href="#">1236</a>


## Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

## Other Formats

Purified anti-human CD197 (CCR7), Alexa Fluor® 488 anti-human CD197 (CCR7), Brilliant Violet 421™ anti-human CD197 (CCR7), PE anti-human CD197 (CCR7), APC/Cyanine7 anti-human CD197 (CCR7), Pacific Blue™ anti-human CD197 (CCR7), APC anti-human CD197 (CCR7), FITC anti-human CD197 (CCR7), Alexa Fluor® 647 anti-human CD197 (CCR7), PerCP/Cyanine5.5 anti-human CD197 (CCR7), Brilliant Violet 605™ anti-human CD197 (CCR7), PE/Cyanine7 anti-human CD197 (CCR7), Brilliant Violet 711™ anti-human CD197 (CCR7), Brilliant Violet 785™ anti-human CD197 (CCR7), Brilliant Violet 510™ anti-human CD197 (CCR7), Brilliant Violet 650™ anti-human CD197 (CCR7), PE/Dazzle™ 594 anti-human CD197 (CCR7), Biotin anti-human CD197 (CCR7), Purified anti-human CD197 (CCR7) (Maxpar® Ready), PerCP anti-human CD197 (CCR7), Alexa Fluor® 700 anti-human CD197 (CCR7), APC/Fire™ 750 anti-human CD197 (CCR7), TotalSeq™-A0148 anti-human CD197 (CCR7), TotalSeq™-B0148 anti-human CD197 (CCR7), TotalSeq™-C0148 anti-human CD197 (CCR7), Brilliant Violet 750™ anti-human CD197 (CCR7), Ultra-LEAF™ Purified anti-human CD197 (CCR7), Spark NIR™ 685 anti-human CD197 (CCR7), KIRAVIA Blue 520™ anti-human CD197 (CCR7), PE/Fire™ 640 anti-human CD197 (CCR7), Spark YG™ 581 anti-human CD197 (CCR7), APC/Fire™ 810 anti-human CD197 (CCR7) Antibody, TotalSeq™-D0148 anti-human CD197 (CCR7), PE/Fire™ 810 anti-human CD197 (CCR7) Antibody, PE/Cyanine5 anti-human CD197 (CCR7), Spark YG™ 593 anti-human CD197 (CCR7), PE/Fire™ 744 anti-human CD197 (CCR7), Spark Red™ 718 anti-human CD197 (CCR7), PE anti-human CD197, Spark PLUS UV395™ anti-human CD197 (CCR7), PerCP/Fire™ 780 anti-human CD197 (CCR7), Brilliant Violet 570™ anti-human CD197 (CCR7) Antibody

## Product Data

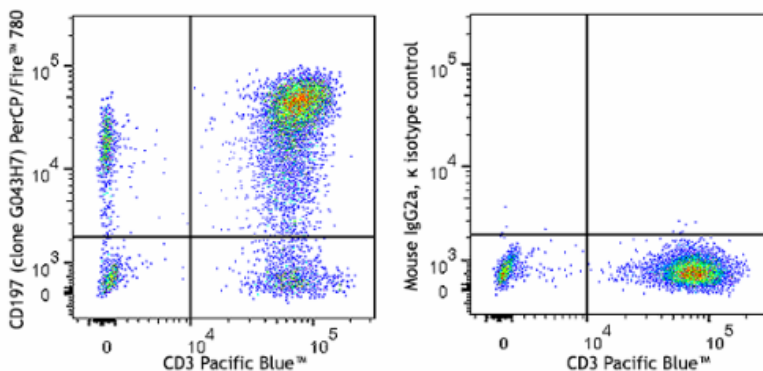

Human peripheral blood lymphocytes were stained with anti-human CD3 (clone UCHT1) Pacific Blue™ and anti-human CD197 (CCR7) (clone G043H7) PerCP/Fire™ 780 (left) or mouse IgG2a, κ PerCP/Fire™ 780 isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at [www.biolegend.com/terms](http://www.biolegend.com/terms) ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587