

PerCP/Fire™ 780 anti-human CD64 Antibody

Catalog# / Size	399511 / 25 tests 399512 / 100 tests
Clone	S18012C
Regulatory Status	RUO
Other Names	FcγRI, FcR I
Isotype	Mouse IgG1, κ
Description	CD64 is a 72 kD single chain type I glycoprotein also known as FcγRI and FcR I. CD64 is a member of the immunoglobulin superfamily and is expressed on monocytes/macrophages, dendritic cells, and activated granulocytes. The expression can be upregulated by IFN-γ stimulation. CD64 binds IgG immune complex. It plays a role in antigen capture, phagocytosis of IgG/antigen complexes, and antibody-dependent cellular cytotoxicity (ADCC).

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Recombinant human CD64-Fc (1257-FC-050)
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with PerCP/Fire™ 780 under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * PerCP/Fire 780 has a maximum excitation of 478 nm and a maximum emission of 780 nm.
Excitation Laser	Blue Laser (488 nm)
Application Notes	In-house testing suggests that clones S18012C and 10.1 do not cross-block each other and recognize different epitopes. Clone 10.1 recognizes the EC3 epitope of CD64. While both CD64 subunits contain EC3 domains, clone 10.1 preferentially binds to CD64A (FcγRIA) but not CD64B (FcγRIB).
RRID	AB_3662297 (BioLegend Cat. No. 399511) AB_3662297 (BioLegend Cat. No. 399512)

Antigen Details

Structure	Ig superfamily, type I glycoprotein, 72 kD
Distribution	Monocytes, macrophages, dendritic cells, activated granulocytes
Function	Phagocytosis, ADCC

Ligand/Receptor	IgG receptor
Cell Type	Dendritic cells, Granulocytes, Macrophages, Monocytes
Biology Area	Immunology, Inhibitory Molecules
Molecular Family	CD Molecules, Fc Receptors

Antigen References


1. HHulett MD & Hogarth PM. 1994. *Adv Immunol.* 57:1-127.
2. van de Winkel JG & Capel PJ. 1993. *Immunol Today.* 14:215-21.

Gene ID [2209](#)


Other Formats

FITC anti-human CD64, PE anti-human CD64, Purified anti-human CD64, APC anti-human CD64, PE/Cyanine7 anti-human CD64, PerCP/Fire™ 780 anti-human CD64, Spark Blue™ 550 anti-human CD64 (Flexi-Fluor™) Antibody , Spark Blue™ 574 anti-human CD64 (Flexi-Fluor™) Antibody, Spark Red™ 718 anti-human CD64 (Flexi-Fluor™) Antibody

Product Data

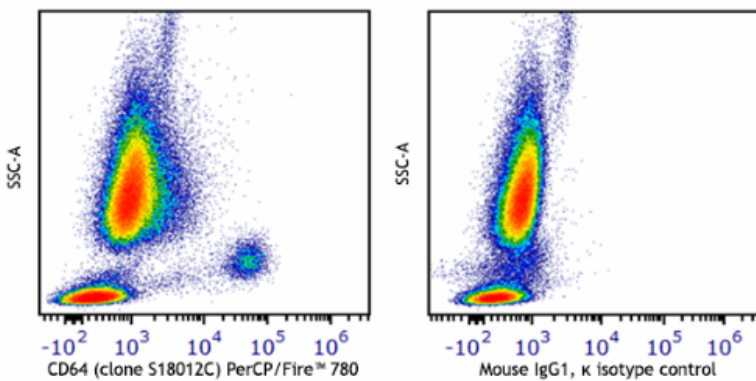

Human peripheral blood cells were stained with anti-human CD64 (clone S18012C) PerCP/Fire™ 780 (left) or mouse IgG1, κ PerCP/Fire™ 780 isotype control (right).

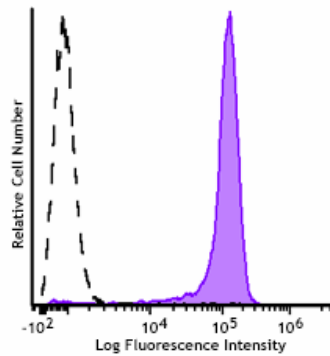

Human peripheral blood monocytes were stained with anti-human CD64 (clone S18012C) PerCP/Fire™ 780 (filled histogram) or mouse IgG1, κ PerCP/Fire™ 780 isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587