

PerCP/Fire™ 806 anti-human CD20 Antibody

Catalog# / Size	375537 / 25 tests 375538 / 100 tests
Clone	S18015E
Regulatory Status	RUO
Other Names	B1, Bp35
Isotype	Mouse IgG2a, κ
Description	CD20 (also known as B1 and Bp35) is a four transmembrane spanning protein measuring 33-37 kD. Expressed on pre-B cells, resting and activated B cells (not plasma cells), some follicular dendritic cells, and at low levels on a T cell subset. Involved in B-cell activation and is associated with a variety of Src family kinases.

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Human CD20 transfectant
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with PerCP/Fire™ 806 under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool).
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * PerCP/Fire™ 806 has a maximum excitation of 478 nm and a maximum emission of 806 nm.
Excitation Laser	Blue Laser (488 nm)
RRID	AB_3717196 (BioLegend Cat. No. 375537) AB_3717196 (BioLegend Cat. No. 375538)

Antigen Details

Structure	Four transmembrane protein (TM4SF), heavily phosphorylated after activation, 33-37 kD
Distribution	B cell, T cell subsets
Function	B cell activation

Interaction	The CD20 molecule is involved in B-cell activation and is associated with various Src family kinases (Lyn, Lck, Fyn).
Ligand/Receptor	Src family tyrosine kinases, MHC class I, II, CD53, CD81, CD82
Cell Type	B cells, T cells
Biology Area	Costimulatory Molecules, Immunology
Molecular Family	CD Molecules
Antigen References	1. Hultin LE, <i>et al.</i> 1993. <i>Cytometry</i> 14:196. 2. Tedder T and Engel P. 1994. <i>Immunol. Today</i> 15:450.
Gene ID	931


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

PerCP/Cyanine5.5 anti-human CD20 Antibody, FITC anti-human CD20 Antibody, PE/Dazzle™ 594 anti-human CD20, Brilliant Violet 421™ anti-human CD20, Pacific Blue™ anti-human CD20, PE/Cyanine7 anti-human CD20, Spark UV™ 387 anti-human CD20, APC/Fire™ 810 anti-human CD20, APC/Fire™ 750 anti-human CD20, Alexa Fluor® 647 anti-human CD20, Spark Violet™ 500 anti-human CD20, PerCP anti-human CD20, Alexa Fluor® 700 anti-human CD20, PerCP/Fire™ 780 anti-human CD20 Antibody, Spark PLUS UV395™ anti-human CD20, PerCP/Fire™ 806 anti-human CD20 Antibody, PE anti-human CD20 Antibody, APC anti-human CD20 Antibody, PE/Fire™ 810 anti-human CD20, Purified anti-human CD20

Product Data

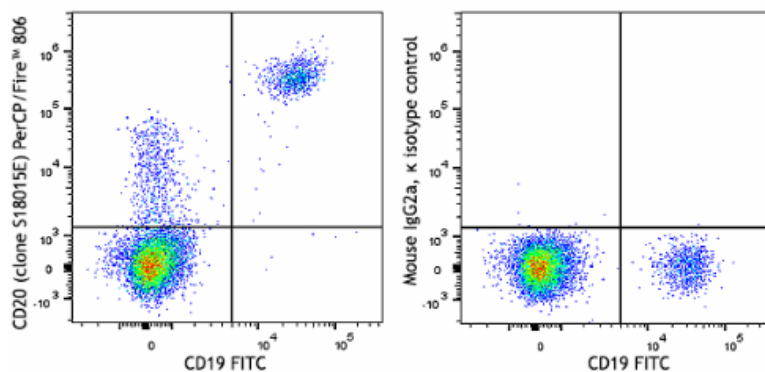

Human peripheral blood lymphocytes were stained with anti-human CD19 (clone HIB19) FITC and anti-human CD20 (clone S18015E) PerCP/Fire™ 806 (left) or mouse IgG2a, κ PerCP/Fire™ 806 isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587