

## PerCP/Fire™ 806 anti-human CD27 Antibody

<b>Catalog# / Size</b>	302873 / 25 tests 302874 / 100 tests
<b>Clone</b>	O323
<b>Regulatory Status</b>	RUO
<b>Workshop</b>	IV T-186
<b>Other Names</b>	S152, T14, TNFRSF7
<b>Isotype</b>	Mouse IgG1, κ
<b>Description</b>	CD27 is a 50-55 kD type I membrane protein also known as S152 and T14. It is a lymphocyte-specific member of the TNF-receptor superfamily. CD27 is expressed on medullary thymocytes, virtually all mature T cells, some B cells, and NK cells. CD27 binds to CD70 and plays an important role in costimulation of T cell activation, and regulation of B cell differentiation and proliferation. The cytoplasmic domains of CD27 have also been shown to interact with TRAF2 and TRAF5 to elicit NF-κB and SAPK/JNK activation.

### Product Details

---

<b>Verified Reactivity</b>	Human, Cynomolgus, Rhesus
<b>Reported Reactivity</b>	African Green, Baboon, Squirrel Monkey
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with PerCP/Fire™ 806 under optimal conditions.
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our <a href="#">Certificate of Analysis</a> online tool.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.  * PerCP/Fire™ 806 has a maximum excitation of 478 nm and a maximum emission of 806 nm.
<b>Excitation Laser</b>	Blue Laser (488 nm)
<b>Application References</b>	1. Knapp W, <i>et al.</i> Eds. 1989. Leucocyte Typing IV. Oxford University Press. New York. 2. Correia DV, <i>et al.</i> 2011. <i>Blood</i> 118:992. (FC) <a href="#">PubMed</a>
<b>(PubMed link indicates BioLegend citation)</b>	
<b>RRID</b>	AB_3698954 (BioLegend Cat. No. 302873) AB_3698954 (BioLegend Cat. No. 302874)

### Antigen Details

---

<b>Structure</b>	TNF-R superfamily, type I transmembrane glycoprotein, 50-55 kD
<b>Distribution</b>	Medullary thymocytes, T and B cell subsets, NK cells
<b>Function</b>	T cell costimulation
<b>Ligand/Receptor</b>	CD70
<b>Cell Type</b>	B cells, NK cells, T cells, Thymocytes, Tregs
<b>Biology Area</b>	Costimulatory Molecules, Immunology
<b>Molecular Family</b>	CD Molecules
<b>Antigen References</b>	1. Hintzen R, <i>et al.</i> 1994. <i>Immunol. Today</i> 15:307. 2. Agematsu K, <i>et al.</i> 1995. <i>J. Immunol.</i> 154:3627.
<b>Gene ID</b>	<a href="#">939</a>


## Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

## Other Formats

APC anti-human CD27, Biotin anti-human CD27, FITC anti-human CD27, PE anti-human CD27, Purified anti-human CD27, Alexa Fluor® 647 anti-human CD27, Alexa Fluor® 700 anti-human CD27, APC/Cyanine7 anti-human CD27, PerCP anti-human CD27, PerCP/Cyanine5.5 anti-human CD27, Pacific Blue™ anti-human CD27, Brilliant Violet 421™ anti-human CD27, Brilliant Violet 570™ anti-human CD27, Brilliant Violet 650™ anti-human CD27, Brilliant Violet 605™ anti-human CD27, Brilliant Violet 711™ anti-human CD27, Brilliant Violet 785™ anti-human CD27, Brilliant Violet 510™ anti-human CD27, PE/Cyanine7 anti-human CD27, Purified anti-human CD27 (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD27, APC/Fire™ 750 anti-human CD27, TotalSeq™-A0154 anti-human CD27, Brilliant Violet 750™ anti-human CD27, TotalSeq™-B0154 anti-human CD27, TotalSeq™-C0154 anti-human CD27, Spark NIR™ 685 anti-human CD27, PE/Fire™ 810 anti-human CD27, TotalSeq™-D0154 anti-human CD27, APC/Fire™ 810 anti-human CD27, PE/Cyanine5 anti-human CD27 Antibody, Spark UV™ 387 anti-human CD27, Spark Blue™ 515 anti-human CD27, Spark Red™ 718 anti-human CD27, Alexa Fluor® 488 anti-human CD27, PerCP/Fire™ 806 anti-human CD27 Antibody, Spark YG™ 581 anti-human CD27 (Flexi-Fluor™), Spark YG™ 593 anti-human CD27 (Flexi-Fluor™) Antibody

## Product Data

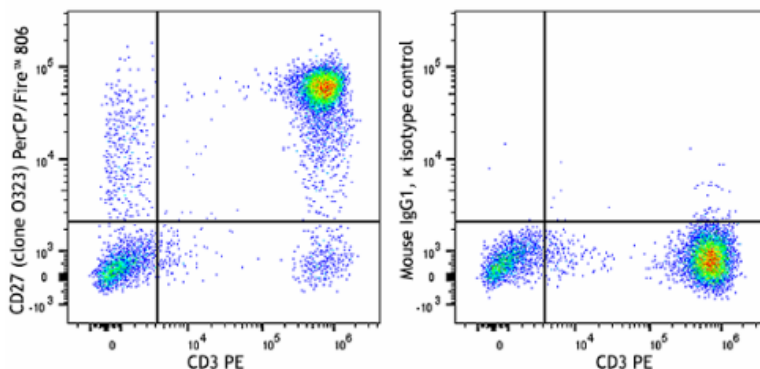

Human peripheral blood lymphocytes were stained with anti-human CD3 (clone UCHT1) PE and anti-human CD27 (clone O323) PerCP/Fire™ 806 (left) or mouse IgG1, κ PerCP/Fire™ 806 isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at [www.biolegend.com/terms](http://www.biolegend.com/terms) ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587