

PerCP/Fire™ 806 anti-human CD38 Antibody

Catalog# / Size	397127 / 25 tests 397128 / 100 tests
Clone	S17015A
Regulatory Status	RUO
Other Names	HB7, T10, ADP-ribosyl cyclase
Isotype	Mouse IgG1, κ
Description	CD38 is a 45 kD type II transmembrane glycoprotein also known as T10. It is an ADP-ribosyl hydrolase expressed at variable levels on hematopoietic cells and in some non-hematopoietic tissues (such as brain, muscles, and kidney). In humans, it is expressed at high levels on plasma cells and activated T and B cells. By functioning as both a cyclase and a hydrolase, CD38 mediates lymphocyte activation, adhesion, and the metabolism of cADPR and NAADP. CD31 is the ligand of CD38.

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Human CD38 transfectants
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with PerCP/Fire™ 806 under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * PerCP/Fire™ 806 has a maximum excitation of 478 nm and a maximum emission of 806 nm.
Excitation Laser	Blue Laser (488 nm)
Application Notes	S17015A does not cross-block binding of clones HIT2, HB-7 or S17015F also raised against human CD38 based on in-house testing.
RRID	AB_3675128 (BioLegend Cat. No. 397127) AB_3675128 (BioLegend Cat. No. 397128)

Antigen Details

Structure	ADP-ribosyl cyclase, ectoenzyme, type II glycoprotein, 45 kD
Distribution	T cells, B cells, NK, myeloid, plasma, and dendritic cells

Function	Ecto-ADP-ribosyl cyclase, calcium signaling, cell activation
Ligand/Receptor	CD31, hyaluronic acid
Cell Type	B cells, Dendritic cells, NK cells, Plasma cells, T cells
Biology Area	Cell Adhesion, Cell Biology, Immunology
Antigen References	<ol style="list-style-type: none"> 1. Ferrero E, <i>et al.</i> 1999. <i>J Leukoc Biol.</i> 65:151. 2. Lund F, <i>et al.</i> 1995. <i>Immunol Today.</i> 16:469.

Gene ID [952](#)


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

PE anti-human CD38, Purified anti-human CD38, APC anti-human CD38, APC/Fire™ 750 anti-human CD38, Alexa Fluor® 700 anti-human CD38, PerCP/Cyanine5.5 anti-human CD38, FITC anti-human CD38, PE/Cyanine7 anti-human CD38, APC/Cyanine7 anti-human CD38, PE/Fire™ 700 anti-human CD38, APC/Fire™ 810 anti-human CD38 Antibody, Brilliant Violet 421™ anti-human CD38 Antibody, Spark PLUS UV395™ anti-human CD38, PerCP/Fire™ 806 anti-human CD38

Product Data

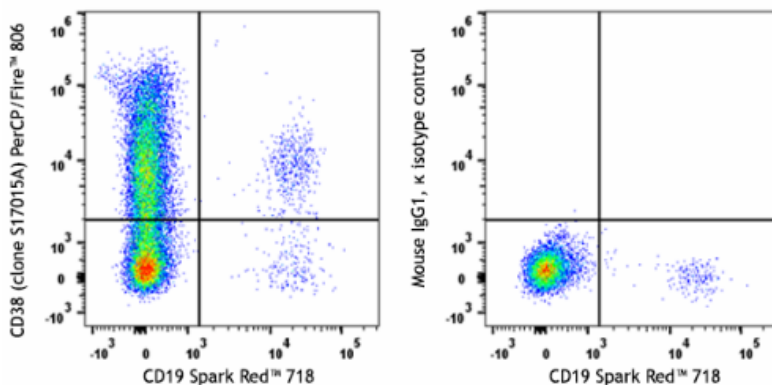

Human peripheral blood lymphocytes were stained with anti-human CD19 (clone HIB19) Spark Red™ 718 and anti-human CD38 (clone S17015A) PerCP/Fire™ 806 (left) or mouse IgG1, κ PerCP/Fire™ 806 isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587