

PerCP/Fire™ 806 anti-human CD45RB Recombinant Antibody

Catalog# / Size	377711 / 25 tests 377712 / 100 tests
Clone	QA20A68
Regulatory Status	RUO
Isotype	Mouse IgG1, κ
Description	CD45RB is a 180-240 kD single chain type I membrane glycoprotein. It is an exon 5 splice variant of the tyrosine phosphatase CD45. The CD45RB isoform is expressed on thymocytes, T cell subsets, B cells, monocytes/macrophages, granulocytes, and dendritic cells. CD45RB enhances both T cell receptor and B cell receptor signaling mediated activation. CD45RB negatively regulates protein kinase activity, and negatively regulates cytokine-mediated signaling pathway.

Product Details

Verified Reactivity	Human
Antibody Type	Recombinant
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with PerCP/Fire™ 806 under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * PerCP/Fire™ 806 has a maximum excitation of 478 nm and a maximum emission of 806 nm.
Excitation Laser	Blue Laser (488 nm)
Application Notes	Clone QA20A68 can completely block the binding of MEM-55 to target cells.
RRID	AB_3683412 (BioLegend Cat. No. 377711) AB_3683412 (BioLegend Cat. No. 377712)

Antigen Details

Structure	Tyrosine phosphatases, type I transmembrane, 180-240 kD (exon 5 splicing of CD45 gene)
Distribution	Highly expressed in bone marrow, colon, lung, and most normal epithelial cells. Also highly expressed on carcinomas of gastrointestinal origin. Expressed during early erythropoiesis.
Function	Enhances TCR and BCR signaling. Upon T-cell activation, recruits and dephosphorylates SKAP1 and FYN.

Interaction	Human
Ligand/Receptor	Galectin-1, CD2, CD3, CD4
Cell Type	B cells, Dendritic cells, Granulocytes, Macrophages, Monocytes, T cells, Thymocytes, Tregs
Biology Area	Cell Biology, Immunology, Inhibitory Molecules, Neuroscience, Neuroscience Cell Markers
Molecular Family	CD Molecules

Antigen References

1. Thomas ML. 1989. *Annu Rev Immunol.* 7:339-69.
2. Trowbridge IS, et al. 1994. *Annu Rev Immunol.* 12:85-116.

Gene ID [5788](#)


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Purified anti-human CD45RB Recombinant Antibody, PE anti-human CD45RB Recombinant Antibody, FITC anti-human CD45RB Recombinant Antibody, Alexa Fluor® 594 anti-human CD45RB Recombinant Antibody, APC anti-human CD45RB Recombinant Antibody, Spark PLUS UV395™ anti-human CD45RB Recombinant Antibody, Spark Red™ 718 anti-human CD45RB Recombinant Antibody, PerCP/Fire™ 806 anti-human CD45RB Recombinant Antibody

Product Data

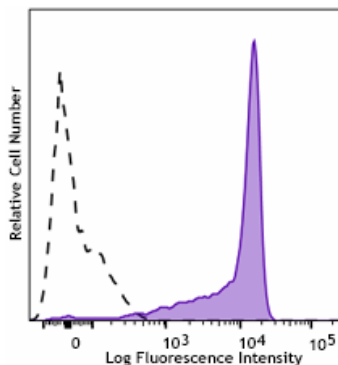

Human peripheral blood lymphocytes were stained with anti-human CD45RB recombinant (clone QA20A68) PerCP/Fire™ 806 (filled histogram) or mouse IgG1, κ PerCP/Fire™ 806 isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587