

PerCP/Fire™ 806 anti-human CD45RO Antibody

| | |
|--------------------------|---|
| Catalog# / Size | 304275 / 25 tests 304276 / 100 tests |
| Clone | UCHL1 |
| Regulatory Status | RUO |
| Workshop | IV N31 |
| Other Names | CD45RO |
| Isotype | Mouse IgG2a, κ |
| Description | CD45RO is a 180 kD single chain membrane glycoprotein. It is a splice variant of tyrosine phosphatase CD45, lacking the A, B, and C determinants. The CD45RO isoform is expressed on activated and memory T cells, some B cell subsets, activated monocytes/macrophages, and granulocytes. CD45RO enhances both T cell receptor and B cell receptor signaling mediated activation. CD45 and its isoforms non-covalently associate with lymphocyte phosphatase-associated phosphoprotein (LPAP) on T and B lymphocytes. CD45 has been reported to be associated with several other cell surface antigens including CD1, CD2, CD3, and CD4. CD45 has also been reported to bind galectin-1 and CD22. CD45 isoform expression can change in response to cytokines. |

Product Details

| | |
|-------------------------------|---|
| Verified Reactivity | Human |
| Reported Reactivity | Chimpanzee, Cynomolgus, Common Marmoset |
| Antibody Type | Monoclonal |
| Host Species | Mouse |
| Immunogen | IL-2 dependent T cell line, CA1 |
| Formulation | Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA) |
| Preparation | The antibody was purified by affinity chromatography and conjugated with PerCP/Fire™ 806 under optimal conditions. |
| Concentration | Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.) |
| Storage & Handling | The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze. |
| Application | FC - Quality tested |
| Recommended Usage | Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * PerCP/Fire™ 806 has a maximum excitation of 478 nm and a maximum emission of 806 nm. |
| Excitation Laser | Blue Laser (488 nm) |
| Application Notes | The UCHL1 antibody is commonly used in combination with antibodies against CD45RA to discern memory and naïve T cells. Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen tissue sections ⁵ and formalin-fixed paraffin-embedded tissue sections ⁴ , Western blotting ² , and immunoprecipitation ³ . |

Application References

(PubMed link indicates
BioLegend citation)

1. Knapp W, *et al.* Eds. 1989. Leucocyte Typing IV. Oxford University Press. New York. (FC)
2. Ishii T, *et al.* 2001. *P. Natl. Acad. Sci. USA* 98:12138. (WB)
3. Ponsford M, *et al.* 2001. *Clin. Exp. Immunol.* 124:315. (IP)
4. Yamada M, *et al.* 1996. *Stroke* 27:1155. (IHC)
5. Sakkas LI, *et al.* 1998. *Clin. Diagn. Lab. Immunol.* 5:430. (IHC)
6. Baba N, *et al.* 2010. *Int. Immunol.* 22:237. [PubMed](#)
7. Thakral D, *et al.* 2008. *J. Immunol.* 180:7431. (FC) [PubMed](#)
8. Weiss L, *et al.* 2010. *P. Natl. Acad. Sci. USA* 107:10632. [PubMed](#)
9. Wu YY, *et al.* 2007. *Infect. Immun.* 75:4357. [PubMed](#)
10. Mozaffarian N, *et al.* 2008. *Rheumatology* 47:1335. [PubMed](#)
11. Roque S, *et al.* 2007. *J. Immunol.* 178:8028. [PubMed](#)
12. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)

[See More](#)

RRID

AB_3675030 (BioLegend Cat. No. 304275)
AB_3675030 (BioLegend Cat. No. 304276)

Antigen Details

| | |
|---------------------------|--|
| Structure | Tyrosine phosphatases, type I transmembrane, 180 kD (isoform of CD45 containing none of the A, B, or C determinants) |
| Distribution | Activated and memory T cells, B cell subsets, monocytes, macrophages, granulocytes |
| Function | Enhances TCR and BCR signaling |
| Ligand/Receptor | CD22 |
| Cell Type | B cells, Granulocytes, Macrophages, Mesenchymal Stem Cells, Monocytes, T cells, Tregs |
| Biology Area | Cell Biology, Immunology, Inhibitory Molecules, Neuroscience, Neuroscience Cell Markers, Stem Cells |
| Molecular Family | CD Molecules |
| Antigen References | 1. Thomas M. 1989. <i>Annu. Rev. Immunol.</i> 7:339. 2. Trowbridge I, <i>et al.</i> 1994. <i>Annu. Rev. Immunol.</i> 12:85. |
| Gene ID | 5788 |


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-human CD45RO, FITC anti-human CD45RO, PE anti-human CD45RO, PE/Cyanine5 anti-human CD45RO, Purified anti-human CD45RO, Alexa Fluor® 488 anti-human CD45RO, Pacific Blue™ anti-human CD45RO, Alexa Fluor® 700 anti-human CD45RO, Biotin anti-human CD45RO, Brilliant Violet 421™ anti-human CD45RO, PerCP/Cyanine5.5 anti-human CD45RO, Brilliant Violet 570™ anti-human CD45RO, Brilliant Violet 605™ anti-human CD45RO, APC/Cyanine7 anti-human CD45RO, PE/Cyanine7 anti-human CD45RO, Brilliant Violet 650™ anti-human CD45RO, Brilliant Violet 711™ anti-human CD45RO, Brilliant Violet 785™ anti-human CD45RO, Alexa Fluor® 594 anti-human CD45RO, Purified anti-human CD45RO (Maxpar® Ready), Brilliant Violet 510™ anti-human CD45RO, PE/Dazzle™ 594 anti-human CD45RO, APC/Fire™ 750 anti-human CD45RO, PerCP anti-human CD45RO, TotalSeq™-A0087 anti-human CD45RO, TotalSeq™-B0087 anti-human CD45RO, TotalSeq™-C0087 anti-human CD45RO, Brilliant Violet 750™ anti-human CD45RO, PE/Fire™ 640 anti-human CD45RO, TotalSeq™-D0087 anti-human CD45RO, Spark Violet™ 423 anti-human CD45RO, GMP APC anti-human CD45RO, TotalSeq™-Bn1369 anti-human CD45RO, GMP PE anti-human CD45RO, Spark UV™ 387 anti-human CD45RO, APC/Fire™ 810 anti-human CD45RO, PerCP/Fire™ 806 anti-human CD45RO, PerCP/Fire™ 780 anti-human CD45RO Antibody, Spark YG™ 581 anti-human CD45RO (Flexi-Fluor™), Spark YG™ 593 anti-human CD45RO (Flexi-Fluor™) Antibody

Product Data

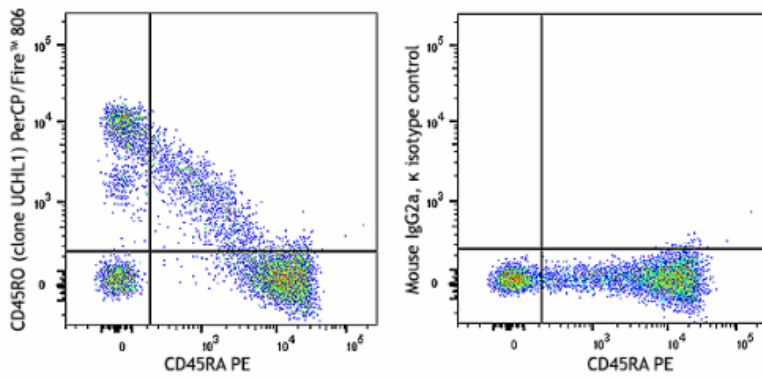

Human peripheral blood lymphocytes were stained with anti-human CD45RA (clone HI100) PE and anti-human CD45RO (clone UCHL1) PerCP/Fire™ 806 (left) or mouse IgG2a, κ PerCP/Fire™ 806 isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587