

PerCP/Fire™ 806 anti-mouse CD25 Antibody

Catalog# / Size	102085 / 25 µg 102086 / 100 µg
Clone	PC61
Regulatory Status	RUO
Other Names	IL-2R α , Ly-43, p55, Tac
Isotype	Rat IgG1, λ
Description	CD25 is a 55 kD glycoprotein also known as the low affinity IL-2R α , Ly-43, p55, or Tac. It is expressed on activated T and B cells, thymocyte subsets, pre-B cells, and T regulatory cells. In association with CD122 (IL-2R β) and CD132 (common γ chain), CD25 forms the high affinity signaling IL-2 receptor.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Rat
Immunogen	IL-2-dependent cytolytic mouse T-cell clone B6.1
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Preparation	The antibody was purified by affinity chromatography and conjugated with PerCP/Fire™ 806 under optimal conditions.
Concentration	0.2 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is \leq 0.125 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application. * PerCP/Fire™ 806 has a maximum excitation of 478 nm and a maximum emission of 806 nm.
Excitation Laser	Blue Laser (488 nm)
Application Notes	Additional reported applications (for the relevant formats) include: immunoprecipitation ^{1,2} , <i>in vitro</i> blocking of IL-2 binding to low- and high-affinity receptors ¹⁻⁴ , growth inhibition of IL-2-dependent T-cell lines ¹⁻⁴ , <i>in vivo</i> depletion of CD25 ⁺ CD4 ⁺ Treg cells ^{5-8,10} , and immunohistochemical staining of acetone-fixed frozen sections ² . PC61 antibody recognizes a different epitope than 3C7 antibody (Cat. No. 101902). For <i>in vivo</i> studies or highly sensitive assays, we recommend Ultra-LEAF™ purified antibody (Cat. No. 102040) with endotoxin < 0.01 EU/µg, Azide-Free, 0.2 µm filtered.

Application References

(PubMed link indicates BioLegend citation)

1. Lowenthal JW, *et al.* 1985. *Nature* 315:669. (IP, Block)
2. Ceredig R, *et al.* 1985. *Nature* 314:98. (IP, IHC, Block)
3. Lowenthal JW, *et al.* 1985. *J. Immunol.* 135:3988. (Block)
4. Moreau JL, *et al.* 1987. *Eur. J. Immunol.* 17:929. (Block)
5. Takahashi T, *et al.* 2000. *J. Exp. Med.* 192:303. (Deplete)
6. Onizuka S, *et al.* 1999. *Cancer Res.* 59:3128. (Deplete)
7. Lei TC, *et al.* 2005. *Blood* 105:4865. (Deplete)
8. Pasare C, *et al.* 2004. *Immunity* 21:733. (Deplete)
9. León-Ponte M, *et al.* 2007. *Blood* 109:3139.
10. Cao OW, *et al.* 2007. *Blood* doi:10.1182/blood-2007-02-073304. (Deplete)
11. Benson MJ, *et al.* 2007. *J. Exp. Med.* doi:10.1084/jem.20070719.
12. Liu F, *et al.* 2011. *Arch Toxicol.* 85:1383. [PubMed](#)

[See More](#)

RRID

AB_3674963 (BioLegend Cat. No. 102085)
AB_3674963 (BioLegend Cat. No. 102086)

Antigen Details

Structure	Forms high affinity IL-2R with IL-2R β (CD122) and IL-2R γ (CD132), 55 kD
Distribution	Activated T cells and B cells, thymocyte subset, pre-B cells, T regulatory cells
Function	IL-2 receptor
Ligand/Receptor	IL-2
Cell Type	B cells, T cells, Thymocytes, Tregs
Biology Area	Immunology
Molecular Family	CD Molecules, Cytokine/Chemokine Receptors

Antigen References

1. Taniguchi T, *et al.* 1993. *Cell* 73:5-8.
2. Waldmann TA. 1991. *J Biol Chem.* 266:2681-4.
3. Read S, *et al.* 2000. *J Exp Med.* 192:295-302.
4. Lowenthal JW, *et al.* 1985. *J Immunol.* 135:3988-94.

Gene ID

[16184](#)


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-mouse CD25, Biotin anti-mouse CD25, FITC anti-mouse CD25, PE anti-mouse CD25, PE/Cyanine5 anti-mouse CD25, Purified anti-mouse CD25, PE/Cyanine7 anti-mouse CD25, Alexa Fluor® 488 anti-mouse CD25, Alexa Fluor® 647 anti-mouse CD25, Pacific Blue™ anti-mouse CD25, Alexa Fluor® 700 anti-mouse CD25, APC/Cyanine7 anti-mouse CD25, PerCP/Cyanine5.5 anti-mouse CD25, PerCP anti-mouse CD25, Brilliant Violet 421™ anti-mouse CD25, Brilliant Violet 605™ anti-mouse CD25, Brilliant Violet 650™ anti-mouse CD25, Ultra-LEAF™ Purified anti-mouse CD25, Brilliant Violet 510™ anti-mouse CD25, PE/Dazzle™ 594 anti-mouse CD25, Brilliant Violet 711™ anti-mouse CD25, Brilliant Violet 785™ anti-mouse CD25, Alexa Fluor® 594 anti-mouse CD25, APC/Fire™ 750 anti-mouse CD25, TotalSeq™-A0097 anti-mouse CD25, KIRAVIA Blue 520™ anti-mouse CD25, TotalSeq™-B0097 anti-mouse CD25, TotalSeq™-C0097 anti-mouse CD25, Spark NIR™ 685 anti-mouse CD25 Antibody, PE/Fire™ 640 anti-mouse CD25, Spark YG™ 581 anti-mouse CD25, APC/Fire™ 810 anti-mouse CD25, Brilliant Violet 750™ anti-mouse CD25, PerCP/Fire™ 780 anti-mouse CD25, PE/Fire™ 700 anti-mouse CD25, Spark PLUS UV395™ anti-mouse CD25, Spark Blue™ 574 anti-mouse CD25 (Flexi-Fluor™), PerCP/Fire™ 806 anti-mouse CD25 Antibody, Spark PLUS B550™ anti-mouse CD25, Spark YG™ 593 anti-mouse CD25 (Flexi-Fluor™) Antibody, Spark UV™ 387 anti-mouse CD25 (Flexi-Fluor™), Spark PLUS V475™ anti-mouse CD25 Antibody

Product Data

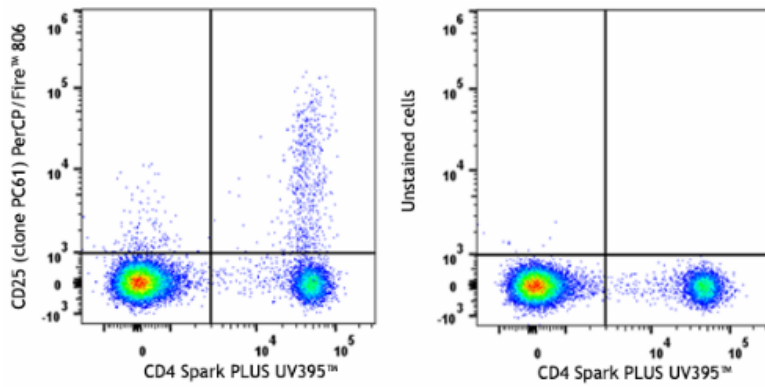

C57BL/6 splenocytes were stained with anti-mouse CD4 (clone GK1.5) Spark PLUS UV395™ and anti-mouse CD25 (clone PC61) PerCP/Fire™ 806 (left) or stained with anti-mouse CD4 (clone GK1.5) Spark PLUS UV395™ only (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587