

## PerCP/Fire™ 806 anti-mouse CD4 Antibody

<b>Catalog# / Size</b>	100497 / 25 µg 100498 / 100 µg
<b>Clone</b>	GK1.5
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	L3T4, T4
<b>Isotype</b>	Rat IgG2b, κ
<b>Description</b>	CD4 is a 55 kD protein also known as L3T4 or T4. It is a member of the Ig superfamily, primarily expressed on most thymocytes, a subset of T cells, and weakly on macrophages and dendritic cells. It acts as a coreceptor with the TCR during T cell activation and thymic differentiation by binding MHC class II and associating with the protein tyrosin kinase, Lck.

### Product Details

<b>Verified Reactivity</b>	Mouse
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Rat
<b>Immunogen</b>	Mouse CTL clone V4
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with PerCP/Fire™ 806 under optimal conditions.
<b>Concentration</b>	0.2 mg/mL
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is ≤ 0.25 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application.  * PerCP/Fire™ 806 has a maximum excitation of 478 nm and a maximum emission of 806 nm.
<b>Excitation Laser</b>	Blue Laser (488 nm)
<b>Application Notes</b>	Additional reported applications (for the relevant formats) include: blocking of CD4 <sup>+</sup> T cell activation <sup>1,4,11</sup> , thymocyte costimulation <sup>3</sup> , <i>in vitro</i> and <i>in vivo</i> depletion <sup>2,5-8</sup> , blocking of egg-sperm cell adhesion <sup>1,4</sup> , immunohistochemical staining of acetone-fixed frozen sections <sup>9,10</sup> , immunoprecipitation <sup>1,2</sup> , and spatial biology (IBEX) <sup>12,13</sup> . The GK1.5 antibody is able to block CD4 mediated cell adhesion and T cell activation. Binding of GK1.5 antibody to CD4 T cells can be blocked by RM4-5 antibody, but not RM4-4 antibody. For <i>in vivo</i> studies or highly sensitive assays, we recommend Ultra-LEAF™ purified antibody (Cat. No. 100442) with a lower endotoxin limit than standard LEAF™ purified antibodies (Endotoxin < 0.01 EU/µg).

## Application References

(PubMed link indicates  
BioLegend citation)

1. Dialynas DP, *et al.* 1983. *J. Immunol.* 131:2445. (Block, IP)
2. Dialynas DP, *et al.* 1983. *Immunol. Rev.* 74:29. (IP, Deplete)
3. Wu L, *et al.* 1991. *J. Exp. Med.* 174:1617. (Costim)
4. Godfrey DI, *et al.* 1994. *J. Immunol.* 152:4783. (Block)
5. Gavett SH, *et al.* 1994. *Am. J. Respir. Cell. Mol. Biol.* 10:587. (Deplete)
6. Schuyler M, *et al.* 1994. *Am. J. Respir. Crit. Care Med.* 149:1286. (Deplete)
7. Ghobrial RR, *et al.* 1989. *Clin. Immunol. Immunopathol.* 52:486. (Deplete)
8. Israelski DM, *et al.* 1989. *J. Immunol.* 142:954. (Deplete)
9. Zheng B, *et al.* 1996. *J. Exp. Med.* 184:1083. (IHC)
10. Frei K, *et al.* 1997. *J. Exp. Med.* 185:2177. (IHC)
11. Felix NJ, *et al.* 2007. *Nat. Immunol.* 8:388. (Block)
12. Radtke AJ, *et al.* 2020. *Proc Natl Acad Sci U S A.* 117:33455-65. (SB) [PubMed](#)

[See More](#)

## RRID

AB\_3698888 (BioLegend Cat. No. 100497)  
AB\_3698888 (BioLegend Cat. No. 100498)

## Antigen Details

---

<b>Structure</b>	Ig superfamily, 55 kD
<b>Distribution</b>	Majority of thymocytes, T cell subset
<b>Function</b>	TCR co-receptor, T cell activation
<b>Ligand/Receptor</b>	MHC class II molecule
<b>Cell Type</b>	Dendritic cells, T cells, Thymocytes, Tregs
<b>Biology Area</b>	Immunology
<b>Molecular Family</b>	CD Molecules
<b>Antigen References</b>	<ol style="list-style-type: none"><li>1. Barclay A, <i>et al.</i> 1997. <i>The Leukocyte Antigen FactsBook</i> Academic Press.</li><li>2. Bierer BE, <i>et al.</i> 1989. <i>Annu. Rev. Immunol.</i> 7:579.</li><li>3. Janeway CA. 1992. <i>Annu. Rev. Immunol.</i> 10:645.</li></ol>
<b>Gene ID</b>	<a href="#">12504</a>

## Related Protocols

---

- [Cell Surface Flow Cytometry Staining Protocol](#)


## Other Formats

---

APC anti-mouse CD4, Biotin anti-mouse CD4, FITC anti-mouse CD4, PE anti-mouse CD4, PE/Cyanine5 anti-mouse CD4, Purified anti-mouse CD4, PE/Cyanine7 anti-mouse CD4, APC/Cyanine7 anti-mouse CD4, Alexa Fluor® 647 anti-mouse CD4, Alexa Fluor® 488 anti-mouse CD4, Pacific Blue™ anti-mouse CD4, Alexa Fluor® 700 anti-mouse CD4, PerCP anti-mouse CD4, PerCP/Cyanine5.5 anti-mouse CD4, Brilliant Violet 421™ anti-mouse CD4, Ultra-LEAF™ Purified anti-mouse CD4, Alexa Fluor® 594 anti-mouse CD4, Brilliant Violet 711™ anti-mouse CD4, Brilliant Violet 510™ anti-mouse CD4, Brilliant Violet 605™ anti-mouse CD4, Brilliant Violet 785™ anti-mouse CD4, PE/Dazzle™ 594 anti-mouse CD4, APC/Fire™ 750 anti-mouse CD4, GoInVivo™ Purified anti-mouse CD4, Brilliant Violet 750™ anti-mouse CD4, Brilliant Violet 650™ anti-mouse CD4, Spark Blue™ 550 anti-mouse CD4, Spark NIR™ 685 anti-mouse CD4, KIRAVIA Blue 520™ anti-mouse CD4, PE/Fire™ 640 anti-mouse CD4, APC/Fire™ 810 anti-mouse CD4, PE/Fire™ 700 anti-mouse CD4, Spark Violet™ 538 anti-mouse CD4, Spark YG™ 593 anti-mouse CD4, Spark Blue™ 574 anti-mouse CD4 Antibody, Spark UV™ 387 anti-mouse CD4, Spark Blue™ 515 anti-mouse CD4, Spark PLUS UV395™ anti-mouse CD4, Spark Red™ 718 anti-mouse CD4 (Flexi-Fluor™), PerCP/Fire™ 806 anti-mouse CD4 Antibody, PerCP/Fire™ 780 anti-mouse CD4 Antibody, Spark YG™ 581 anti-mouse CD4 (Flexi-Fluor™), StarBright UltraViolet 795 anti-mouse CD4, Spark PLUS V475™ anti-mouse CD4

## Product Data

---

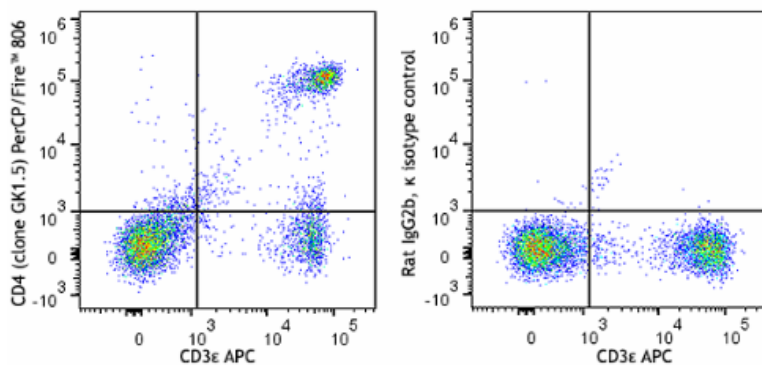

C57BL/6 mouse splenocytes were stained with anti-mouse CD3ε (clone 145-2C11) APC and anti-mouse CD4 (clone GK1.5) PerCP/Fire™ 806 (left) or rat IgG2b, κ PerCP/Fire™ 806 isotype control (right)

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at [www.biolegend.com/terms](http://www.biolegend.com/terms) ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587