

## PerCP/Fire™ 806 anti-mouse Ly-6G Antibody

<b>Catalog# / Size</b>	127683 / 25 µg 127684 / 100 µg
<b>Clone</b>	1A8
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	Lymphocyte antigen 6 complex, locus G
<b>Isotype</b>	Rat IgG2a, κ
<b>Description</b>	Lymphocyte antigen 6 complex, locus G (Ly-6G), a 21-25 kD GPI-anchored protein, is expressed on the majority of myeloid cells in bone marrow and peripheral granulocytes.

### Product Details

---

<b>Verified Reactivity</b>	Mouse
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Rat
<b>Immunogen</b>	Ly-6G transfected EL-4J cell line.
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with PerCP/Fire™ 806 under optimal conditions.
<b>Concentration</b>	0.2 mg/mL
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is ≤ 0.5 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application.  * PerCP/Fire™ 806 has a maximum excitation of 478 nm and a maximum emission of 806 nm.
<b>Excitation Laser</b>	Blue Laser (488 nm)
<b>Application Notes</b>	While 1A8 recognizes only Ly-6G, clone RB6-8C5 recognizes both Ly-6G and Ly-6C. Clone RB6-8C5 binds with high affinity to mouse Ly-6G molecules and to a lower extent to Ly-6C <sup>15</sup> . Clone RB6-8C5 impairs the binding of anti-mouse Ly-6G clone 1A8 <sup>15</sup> . However, clone RB6-8C5 is able to stain in the presence of anti-mouse Ly-6C clone HK1.4 <sup>16</sup> .  Additional reported applications (for the relevant formats) include: immunohistochemistry <sup>9</sup> of frozen sections <sup>10</sup> and paraffin-embedded sections <sup>11</sup> , depletion <sup>4, 12-14</sup> , and spatial biology (IBEX) <sup>20,21</sup> . The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for <i>in vivo</i> studies or highly sensitive assays (Cat. No. 127632, 127649, 127650, 127661 and 127662).

## Application References

(PubMed link indicates  
BioLegend citation)

1. Fleming TJ, *et al.* 1993. *J. Immunol.* 151:2399. (FC)
2. Daley JM, *et al.* 2008. *J. Leukocyte Biol.* 83:1. (FC)
3. Dietlin TA, *et al.* 2007. *J. Leukocyte Biol.* 81:1205. (FC)
4. Daley J, *et al.* 2007. *J. Leukocyte Biol.* doi:10.1189. (Deplete) [PubMed](#)
5. Tadagavadi RK, *et al.* 2010. *J. Immunol.* 185:4904. [PubMed](#)
6. Sumagin R, *et al.* 2010. *J. Immunol.* 185:7057. [PubMed](#)
7. Guiducci C, *et al.* 2010. *J. Exp Med.* 207:2931. [PubMed](#)
8. Fujita M, *et al.* 2011. *Cancer Res.* 71:2664. [PubMed](#)
9. Van Leeuwen, *et al.* 2008. *Arterioscler. Thromb. Vasc. Biol.* 28:84. (IHC)
10. Kowanez M, *et al.* 2010. *P. Natl. Acad. Sci. USA* 107:21248. [supplementary data] (IHC)
11. Esbona K, *et al.* 2016. *Breast Cancer Res.* 18:35. (IHC)
12. Wojtasiak M, *et al.* 2010. *J. Gen. Virol.* 91:2158. (FC, Deplete)

[See More](#)

## RRID

AB\_3674983 (BioLegend Cat. No. 127683)  
AB\_3674983 (BioLegend Cat. No. 127684)

## Antigen Details

---

<b>Structure</b>	A 21-35 kD GPI-anchored membrane protein
<b>Distribution</b>	Expressed on the majority of myeloid cells in bone marrow and peripheral granulocytes. The monoclonal antibody RB6-8C5 recognizes both Ly-6G and Ly-6C.
<b>Cell Type</b>	Granulocytes, Macrophages, Monocytes
<b>Biology Area</b>	Immunology, Innate Immunity
<b>Antigen References</b>	Fleming TJ, <i>et al.</i> 1993. <i>J. Immunol.</i> 151:2399.
<b>Gene ID</b>	<a href="#">546644</a>

## Related Protocols

---

- [Cell Surface Flow Cytometry Staining Protocol](#)


## Other Formats

---

Alexa Fluor® 594 anti-mouse Ly-6G, Purified anti-mouse Ly-6G, Biotin anti-mouse Ly-6G, FITC anti-mouse Ly-6G, PE anti-mouse Ly-6G, Alexa Fluor® 647 anti-mouse Ly-6G, Pacific Blue™ anti-mouse Ly-6G, APC anti-mouse Ly-6G, PerCP/Cyanine5.5 anti-mouse Ly-6G, PE/Cyanine7 anti-mouse Ly-6G, Alexa Fluor® 700 anti-mouse Ly-6G, APC/Cyanine7 anti-mouse Ly-6G, Alexa Fluor® 488 anti-mouse Ly-6G, Brilliant Violet 421™ anti-mouse Ly-6G, Brilliant Violet 570™ anti-mouse Ly-6G, Ultra-LEAF™ Purified anti-mouse Ly-6G, Brilliant Violet 510™ anti-mouse Ly-6G, Purified anti-mouse Ly-6G (Maxpar® Ready), Brilliant Violet 650™ anti-mouse Ly-6G, Brilliant Violet 711™ anti-mouse Ly-6G, Brilliant Violet 605™ anti-mouse Ly-6G, Brilliant Violet 785™ anti-mouse Ly-6G, PE/Dazzle™ 594 anti-mouse Ly-6G, APC/Fire™ 750 anti-mouse Ly-6G, PerCP anti-mouse Ly-6G, TotalSeq™-A0015 anti-mouse Ly-6G, TotalSeq™-C0015 anti-mouse Ly-6G, TotalSeq™-B0015 anti-mouse Ly-6G, Spark Blue™ 550 anti-mouse Ly-6G, Spark NIR™ 685 anti-mouse Ly-6G, Spark YG™ 593 anti-mouse Ly-6G, APC/Fire™ 810 anti-mouse Ly-6G Antibody, PE/Cyanine5 anti-mouse Ly-6G, PE/Fire™ 810 anti-mouse Ly-6G Antibody, Spark UV™ 387 anti-mouse Ly-6G, PE/Fire™ 640 anti-mouse Ly-6G, Spark YG™ 570 anti-mouse Ly-6G, Spark Red™ 718 anti-mouse Ly-6G (Flexi-Fluor™), Spark Blue™ 574 anti-mouse Ly-6G (Flexi-Fluor™), PerCP/Fire™ 806 anti-mouse Ly-6G, StarBright UltraViolet 575 anti-mouse Ly-6G

## Product Data

---

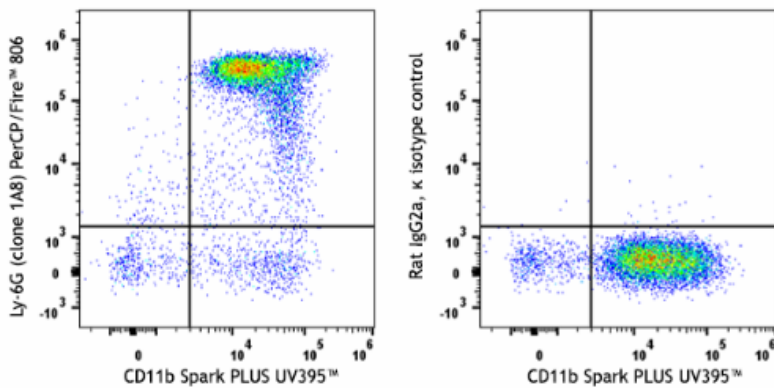

C57BL/6 mouse bone marrow cells were stained with anti-mouse CD11b (clone M1/70) Spark PLUS UV395™ and anti-mouse Ly-6G (clone 1A8) PerCP/Fire™ 806 (left) or rat IgG2a, κ PerCP/Fire™ 806 isotype control (right). Data shown were gated on the myeloid population.

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at [www.biolegend.com/terms](http://www.biolegend.com/terms) ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587