

Spark Blue™ 574 anti-mouse CD90.2 (Thy1.2) Antibody

Catalog# / Size	105359 / 25 µg 105360 / 100 µg
Clone	30-H12
Regulatory Status	RUO
Other Names	Thy-1.2
Isotype	Rat IgG2b, κ
Description	CD90.2 is a 25-35 kD immunoglobulin superfamily member also known as Thy1.2. It is expressed on hematopoietic stem cells and neurons, all thymocytes, and peripheral T cells in Thy1.2 bearing mouse strains (Balb/c, CBA/J, C3H/He, C57BL/6, DBA, NZB/BL). CD90.2 is a glycosylphosphatidylinositol (GPI)-anchored membrane glycoprotein involved in signal transduction. CD90.2 is involved in costimulation of lymphocyte proliferation and induction of hematopoietic stem cells differentiation. CD90.2 has been shown to interact with CD45. The 30-H12 antibody has been reported to induce Ca ²⁺ flux in thymocytes and, in combination with antibody against the CD3/TCR complex, promote thymocyte apoptosis and inhibit CD3-mediated proliferative responses of mature T lymphocytes.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Rat
Immunogen	Mouse thymus or spleen
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Preparation	The antibody was purified by affinity chromatography and conjugated with Spark Blue™ 574 under optimal conditions.
Concentration	0.5 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is ≤ 0.5 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application. * Spark Blue™ 574 has a maximum excitation of 506 nm and a maximum emission of 574 nm.
Excitation Laser	Blue Laser (488 nm)
Application Notes	Additional reported applications (for the relevant formats) include: <i>in vivo</i> and <i>in vitro</i> depletion ^{1,2,7} , costimulation of CD3/TCR-mediated signal transduction ^{3,4} , and immunohistochemical staining ⁵ of acetone-fixed frozen sections. The 30-H12 antibody does not react with Thy-1.1 alloantigen of the AKR/J and PL strains. To reduce non-specific binding to cells bearing Fc-receptors, pre-incubation of cells with anti-mouse CD16/CD32, clone 93 (Cat. No. 101301 & 101302) is recommended prior to immunofluorescent staining. The Ultra-LEAF™ purified antibody (Endotoxin <0.01 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for functional assays (Cat. Nos. 105351 & 105352).

Application References

(PubMed link indicates
BioLegend citation)

1. Hathcock KS. 1991. Current Protocols in Immunology. 3.4.1. (Deplete)
2. Seaman WE. 1983. *J. Immunol.* 130:1713. (Deplete)
3. Nakashima I, *et al.* 1991. *J. Immunol.* 147:1153. (Costim)
4. Nakashima I, *et al.* 1993. *J. Immunol.* 151:3511. (Costim)
5. Ledbetter JA, *et al.* 1980. *J. Exp. Med.* 152:280. (IHC)
6. Hardy B, *et al.* 2005. *Int. Immunol.* 17:615.
7. Drobyski W, *et al.* 1996. *Blood* 87:5355. (Deplete)
8. Dyer KD, *et al.* 2007. *J. Immunol.* 179:1693. (FC) [PubMed](#)
9. Sungur CM, *et al.* 2013. *PNAS.* 110:7401. [PubMed](#)

Antigen Details

Structure	Ig superfamily, 25-35 kD
Distribution	Hematopoietic stem cells and neurons, all thymocytes, peripheral T cells of the Thy-1.2 bearing mice
Function	Lymphocyte costimulation, proliferation and differentiation of hematopoietic stem cells
Ligand/Receptor	CD45
Cell Type	Hematopoietic stem and progenitors, Neurons, T cells, Thymocytes
Biology Area	Immunology
Molecular Family	CD Molecules
Antigen References	1. Barclay A, <i>et al.</i> 1997. The Leukocyte Antigen FactsBook Academic Press. 2. Craig W, <i>et al.</i> 1993. <i>J. Exp. Med.</i> 177:1331. 3. Reif AE and Schlesinger M. 1989. Cell Surface Antigen Thy-1. 4. Mayani H, <i>et al.</i> 1994. <i>Blood</i> 83:2410.
Gene ID	21838


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-mouse CD90.2 (Thy1.2), PE/Cyanine5 anti-mouse CD90.2 (Thy1.2), PE/Fire™ 810 anti-mouse CD90.2 (Thy1.2), Purified anti-mouse CD90.2 (Thy1.2) (Maxpar® Ready), Biotin anti-mouse CD90.2 (Thy1.2), Spark Red™ 718 anti-mouse CD90.2 (Thy1.2) (Flexi-Fluor™), FITC anti-mouse CD90.2 (Thy1.2), Brilliant Violet 510™ anti-mouse CD90.2 (Thy1.2), Brilliant Violet 421™ anti-mouse CD90.2 (Thy1.2), PerCP/Cyanine5.5 anti-mouse CD90.2 (Thy1.2), PE/Cyanine7 anti-mouse CD90.2 (Thy1.2), APC/Cyanine7 anti-mouse CD90.2 (Thy1.2), Brilliant Violet 570™ anti-mouse CD90.2 (Thy1.2), Spark Blue™ 574 anti-mouse CD90.2 (Thy1.2), PE anti-mouse CD90.2 (Thy1.2), Purified anti-mouse CD90.2 (Thy1.2), Spark PLUS B550™ anti-mouse CD90.2 (Thy1.2) (Flexi-Fluor™), Spark PLUS UV395™ anti-mouse CD90.2 (Thy1.2) (Flexi-Fluor™), Spark NIR™ 685 anti-mouse CD90.2 (Thy1.2) (Flexi-Fluor™), Alexa Fluor® 488 anti-mouse CD90.2 (Thy1.2), TotalSeq™-A0075 anti-mouse CD90.2 (Thy1.2), Brilliant Violet 605™ anti-mouse CD90.2 (Thy1.2), PE/Dazzle™ 594 anti-mouse CD90.2 (Thy1.2), Ultra-LEAF™ Purified anti-mouse CD90.2 (Thy1.2), TotalSeq™-C0075 anti-mouse CD90.2 (Thy1.2), TotalSeq™-B0075 anti-mouse CD90.2 (Thy1.2), Spark YG™ 581 anti-mouse CD90.2 (Thy1.2) (Flexi-Fluor™), PerCP anti-mouse CD90.2 (Thy1.2), APC/Fire™ 750 anti-mouse CD90.2 (Thy1.2), Pacific Blue™ anti-mouse CD90.2 (Thy1.2), Brilliant Violet 785™ anti-mouse CD90.2 (Thy1.2), Spark PLUS V475™ anti-mouse CD90.2 (Thy1.2) (Flexi-Fluor™), Spark YG™ 593 anti-mouse CD90.2 (Thy1.2) (Flexi-Fluor™), Alexa Fluor® 647 anti-mouse CD90.2 (Thy1.2), Alexa Fluor® 700 anti-mouse CD90.2 (Thy1.2), Brilliant Violet 711™ anti-mouse CD90.2 (Thy1.2)

Product Data

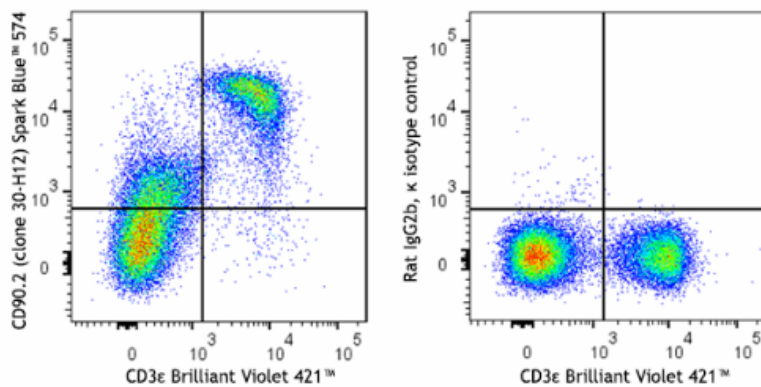

C57BL/6 mouse splenocytes were stained with anti-mouse CD3ε (clone 145-2C11) Brilliant Violet 421™ and anti-mouse CD90.2 (Thy1.2) (clone 30-H12) Spark Blue™ 574 (left) or rat IgG2b, κ Spark Blue™ 574 isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587