

Spark PLUS B550™ anti-human CD223 (LAG-3) Antibody

Catalog# / Size	369355 / 25 tests 369356 / 100 tests
Clone	11C3C65
Regulatory Status	RUO
Other Names	CD223, LAG-3, LAG3, lymphocyte-activation gene-3
Isotype	Mouse IgG1, κ
Description	CD223, also known as LAG-3, is a 70 kD type I transmembrane glycoprotein that is involved in T-cell signaling. Similar to CD4, CD223 binds MHC class II, but with a higher affinity. CD223 negatively regulates T-cell activation. It is expressed by activated T-cells and natural killer cells (NKs), as well as regulatory T-cells. It is transiently expressed on the surface of activated T-cells in acute conditions but high expression is maintained under tolerizing conditions. CD223 deficiency results in reduced tumor growth. CD223 and PD-1 can act in synergy and reverse exhausted phenotypes, improve tumor rejection, and control viral load.

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Human LAG-3 transfected cells.
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Spark PLUS B550™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * Spark PLUS B550™ has a maximum excitation of 516 nm and a maximum emission of 540 nm.
Excitation Laser	Blue Laser (488 nm)
Application Notes	The staining of clone 11C3C65 cannot be blocked by clone 7H2C65, which is another anti-human CD223 (LAG-3) antibody.
RRID	AB_3675087 (BioLegend Cat. No. 369355) AB_3675087 (BioLegend Cat. No. 369356)

Antigen Details

Structure	70 kD transmembrane glycoprotein, Ig superfamily, highly homologous to CD4.
Distribution	Activated T-cells and natural killer cells (NKs) and regulatory T cells.
Function	Negatively regulates T-cell activation.

Ligand/Receptor	Binds MHC class II molecules.
Cell Type	Dendritic cells, NK cells, T cells, Tregs
Biology Area	Immunology, Inhibitory Molecules
Molecular Family	CD Molecules, Immune Checkpoint Receptors
Antigen References	<ol style="list-style-type: none"> 1. Castelli C, <i>et al.</i> 2014. <i>Oncoimmunology</i>. 3(11):e967146. 2. Poirier N, <i>et al.</i> 2011. <i>Clin. Exp. Immunol.</i> 164:265. 3. Juno JA, <i>et al.</i> 2015. <i>Retrovirology</i>. 12:17. 4. Casati C, <i>et al.</i> 2006. <i>Cancer Res.</i> 66:4450.
Gene ID	3902


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)


Other Formats

Purified anti-human CD223 (LAG-3), Alexa Fluor® 647 anti-human CD223 (LAG-3), PE anti-human CD223 (LAG-3), FITC anti-human CD223 (LAG-3), PE/Cyanine7 anti-human CD223 (LAG-3), PerCP/Cyanine5.5 anti-human CD223 (LAG-3), Brilliant Violet 421™ anti-human CD223 (LAG-3), Brilliant Violet 650™ anti-human CD223 (LAG-3), Brilliant Violet 510™ anti-human CD223 (LAG-3), Brilliant Violet 785™ anti-human CD223 (LAG-3), Brilliant Violet 711™ anti-human CD223 (LAG-3), Brilliant Violet 605™ anti-human CD223 (LAG-3), Alexa Fluor® 488 anti-human CD223 (LAG-3), Biotin anti-human CD223 (LAG-3), PE/Dazzle™ 594 anti-human CD223 (LAG-3), APC/Fire™ 750 anti-human CD223 (LAG-3), TotalSeq™-A0152 anti-human CD223 (LAG-3), TotalSeq™-C0152 anti-human CD223 (LAG-3), TotalSeq™-B0152 anti-human CD223 (LAG-3), TotalSeq™-D0152 anti-human CD223 (LAG-3), Alexa Fluor® 700 anti-human CD223 (LAG-3), Pacific Blue™ anti-human CD223 (LAG-3), PE/Cyanine5 anti-human CD223 (LAG-3), APC/Cyanine7 anti-human CD223 (LAG-3), APC/Fire™ 810 anti-human CD223 (LAG-3), Brilliant Violet 750™ anti-human CD223 (LAG-3), Spark Red™ 718 anti-human CD223 (LAG-3) (Flexi-Fluor™), Spark PLUS UV395™ anti-human CD223 (LAG-3), Spark Blue™ 574 anti-human CD223 (LAG-3) (Flexi-Fluor™), Spark PLUS B550™ anti-human CD223 (LAG-3) Antibody

Product Data

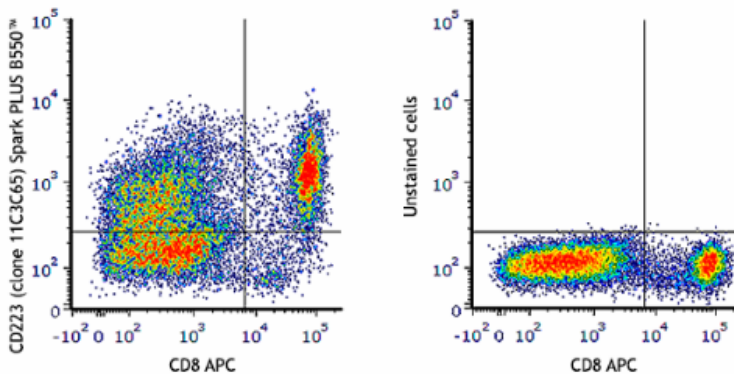

CD3/CD28/IL-2 stimulated (3 days) human peripheral blood mononuclear cells were stained with anti-human CD8 (clone SK1) APC and anti-human CD223 (LAG-3) (clone 11C3C65) Spark PLUS B550™ (left) or stained with anti-human CD8 (clone SK1) APC only (right).

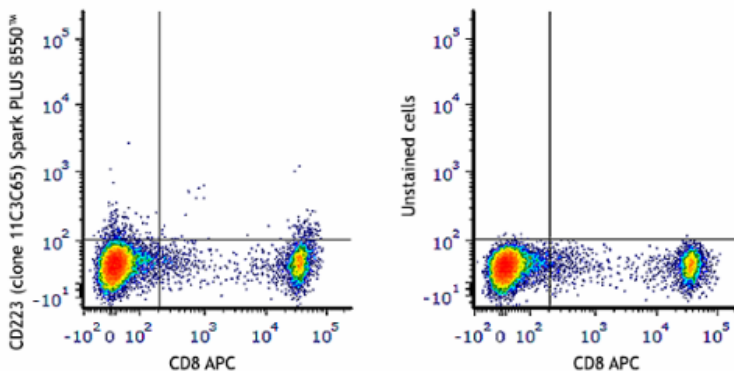

Unstimulated human peripheral blood mononuclear cells were stained with anti-human CD8 (clone SK1) APC and anti-human CD223 (LAG-3) (clone 11C3C65) Spark PLUS B550™ (left) or stained with anti-human CD8 (clone SK1) APC only (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587