

## Spark PLUS B550™ anti-mouse CD11c Antibody

<b>Catalog# / Size</b>	117381 / 25 µg 117382 / 100 µg
<b>Clone</b>	N418
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	α <sub>X</sub> integrin, integrin α <sub>X</sub> chain, CR4, p150, ITGAX
<b>Isotype</b>	Armenian Hamster IgG
<b>Description</b>	CD11c is a 150 kD glycoprotein also known as α <sub>X</sub> integrin, CR4, and p150. CD11c forms a α <sub>X</sub> β <sub>2</sub> heterodimer with β <sub>2</sub> integrin (CD18). It is primarily expressed on dendritic cells, NK cells, a subset of intestinal intraepithelial lymphocytes (IEL), and some activated T cells. The α <sub>X</sub> β <sub>2</sub> integrin plays an important role in cell-cell contact by binding its ligands: iC3b, fibrinogen, and CD54.

### Product Details

---

<b>Verified Reactivity</b>	Mouse
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Armenian Hamster
<b>Immunogen</b>	Mouse spleen dendritic cells
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with Spark PLUS B550™ under optimal conditions.
<b>Concentration</b>	0.2 mg/mL
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is ≤ 1.0 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application.  * Spark PLUS B550™ has a maximum excitation of 516 nm and a maximum emission of 540 nm.
<b>Excitation Laser</b>	Blue Laser (488 nm)
<b>Application Notes</b>	Additional reported applications (for the relevant formats) include: immunoprecipitation <sup>3</sup> , immunohistochemical staining of acetone-fixed frozen sections <sup>3</sup> , immunofluorescence microscopy <sup>5,9</sup> (Alexa Fluor® 488 conjugated N418 was used for IHC in frozen sections <sup>10</sup> ), and spatial biology (IBEX) <sup>22,23</sup> .

## Application References

(PubMed link indicates  
BioLegend citation)

1. Granucci F, *et al.* 1997. *J. Immunol.* 159:1794.
2. Stokes RW, *et al.* 1998. *J. Immunol.* 160:5514.
3. Metlay JP, *et al.* 1990. *J. Exp. Med.* 171:1753. (IHC, IP)
4. Ma XT, *et al.* 2006. *Cancer Research* 66:1169.
5. Chin RK, *et al.* 2006. *J. Immunol.* 177:290. (IF)
6. Cervantes-Barragan L, *et al.* 2007. *Blood* 109:1131. (FC) [PubMed](#)
7. Turnquist HR, *et al.* 2007. *J. Immunol.* 178:7018. (FC) [PubMed](#)
8. Benson MJ, *et al.* 2007. *J. Exp. Med.* doi:10.1084/jem.20070719. (FC) [PubMed](#)
9. You Y, *et al.* 2009. *J. Immunol.* 182:7343. (IF) [PubMed](#)
10. Roland CL, *et al.* 2009. *Mol. Cancer Res.* 8:1761. (IHC, FC) [PubMed](#)
11. Wikstrom M, *et al.* 2006. *J. Immunol.* 177:913. [PubMed](#)
12. Pericolini E, *et al.* 2008. *J. Leukocyte Biol.* 83:1286. [PubMed](#)

[See More](#)

## RRID

AB\_3698899 (BioLegend Cat. No. 117381)  
AB\_3698899 (BioLegend Cat. No. 117382)

## Antigen Details

---

<b>Structure</b>	Integrin $\alpha$ -chain, associates with integrin $\beta_2$ (CD18), 150 kD
<b>Distribution</b>	Dendritic cells, NK cells, intestinal intraepithelial lymphocytes (IEL), some activated T cells
<b>Function</b>	Cellular adhesion
<b>Ligand/Receptor</b>	iC3b, fibrinogen
<b>Cell Type</b>	Dendritic cells, Epithelial cells, NK cells, T cells, Tregs
<b>Biology Area</b>	Cell Adhesion, Cell Biology, Costimulatory Molecules, Immunology, Innate Immunity, Neuroscience, Neuroscience Cell Markers
<b>Molecular Family</b>	Adhesion Molecules, CD Molecules
<b>Antigen References</b>	<ol style="list-style-type: none"><li>1. Barclay A, <i>et al.</i> 1997. <i>The Leukocyte Antigen Facts Book</i> Academic Press.</li><li>2. Springer TA. 1994. <i>Cell</i> 76:301.</li><li>3. Lopez-Rodriguez C, <i>et al.</i> 1996. <i>J. Immunol.</i> 156:3780.</li></ol>
<b>Gene ID</b>	<a href="#">16411</a>

## Related Protocols

---

- [Cell Surface Flow Cytometry Staining Protocol](#)


## Other Formats

---

APC anti-mouse CD11c, Biotin anti-mouse CD11c, FITC anti-mouse CD11c, PE anti-mouse CD11c, Purified anti-mouse CD11c, Alexa Fluor® 488 anti-mouse CD11c, Alexa Fluor® 647 anti-mouse CD11c, PE/Cyanine5 anti-mouse CD11c, PE/Cyanine7 anti-mouse CD11c, Brilliant Violet 605™ anti-mouse CD11c, Alexa Fluor® 700 anti-mouse CD11c, Pacific Blue™ anti-mouse CD11c, APC/Cyanine7 anti-mouse CD11c, PerCP/Cyanine5.5 anti-mouse CD11c, PerCP anti-mouse CD11c, Brilliant Violet 421™ anti-mouse CD11c, Brilliant Violet 570™ anti-mouse CD11c, Brilliant Violet 785™ anti-mouse CD11c, Brilliant Violet 510™ anti-mouse CD11c, Brilliant Violet 650™ anti-mouse CD11c, Purified anti-mouse CD11c (Maxpar® Ready), Alexa Fluor® 594 anti-mouse CD11c, PE/Dazzle™ 594 anti-mouse CD11c, Brilliant Violet 711™ anti-mouse CD11c, APC/Fire™ 750 anti-mouse CD11c, TotalSeq™-A0106 anti-mouse CD11c, Brilliant Violet 750™ anti-mouse CD11c, TotalSeq™-B0106 anti-mouse CD11c, TotalSeq™-C0106 anti-mouse CD11c, KIRAVIA Blue 520™ anti-mouse CD11c, Spark Blue™ 550 anti-mouse CD11c, Spark NIR™ 685 anti-mouse CD11c, Spark UV™ 387 anti-mouse CD11c, Spark Red™ 718 anti-mouse CD11c, Spark Blue™ 515 anti-mouse CD11c, PerCP/Fire™ 806 anti-mouse CD11c, Spark PLUS UV395™ anti-mouse CD11c, Spark Blue™ 574 anti-mouse CD11c (Flexi-Fluor™), Spark PLUS B550™ anti-mouse CD11c, APC/Fire™ 810 anti-mouse CD11c Antibody, PE/Fire™ 640 anti-mouse CD11c Antibody, Spark YG™ 581 anti-mouse CD11c (Flexi-Fluor™), Spark YG™ 593 anti-mouse CD11c (Flexi-Fluor™) Antibody, Spark PLUS V475™ anti-mouse CD11c Antibody

## Product Data

---

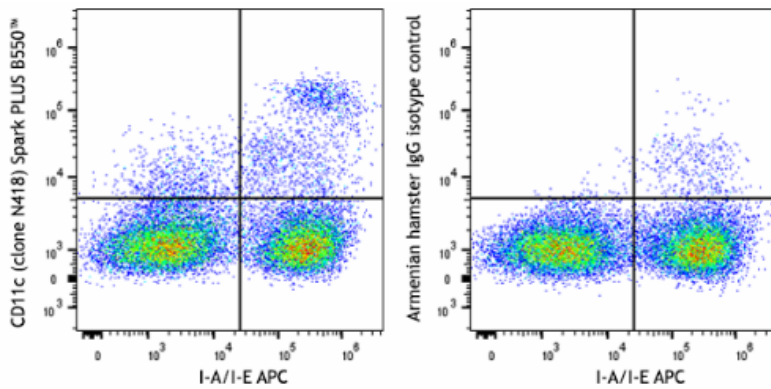

C57BL/6 mouse splenocytes were stained with anti-mouse I-A/I-E (clone M5/114.15.2) APC and anti-mouse CD11c (clone N418) Spark PLUS B550™ (left) or with Armenian hamster IgG Spark PLUS B550™ isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at [www.biolegend.com/terms](http://www.biolegend.com/terms) ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587