

## Spark PLUS B550™ anti-mouse CD69 Antibody

<b>Catalog# / Size</b>	104561 / 25 µg 104562 / 100 µg
<b>Clone</b>	H1.2F3
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	Very Early Activation Antigen (VEA), AIM, EA1, MLR3, gp34/28
<b>Isotype</b>	Armenian Hamster IgG
<b>Description</b>	CD69 is a 60 kD type II membrane protein composed of a 27/33 kD disulfide-linked homodimer, also known as Very Early Activation Antigen (VEA), AIM, EA1, MLR3, and gp34/28. It is expressed on a subset of thymocytes and platelets. CD69 is rapidly induced on activated T and B cells, neutrophils, and NK cells. It is a C-type lectin, closely related to the NKR-P1 and Ly-49 NK cell activation molecules. CD69 is involved in the early events of cell activation and thymocyte positive selection.

### Product Details

---

<b>Verified Reactivity</b>	Mouse
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Armenian Hamster
<b>Immunogen</b>	Mouse dendritic epidermal T cell line Y245
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with Spark PLUS B550™ under optimal conditions.
<b>Concentration</b>	0.2 mg/mL
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is ≤ 1.0 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application.  * Spark PLUS B550™ has a maximum excitation of 516 nm and a maximum emission of 540 nm.
<b>Excitation Laser</b>	Blue Laser (488 nm)
<b>Application Notes</b>	The H1.2F3 antibody has been reported to augment T cell activation. Additional reported applications (for the relevant formats) include: <i>in vitro</i> T cell and NK cell activation <sup>1-3</sup> , immunohistochemistry <sup>4,5</sup> , and immunoprecipitation <sup>1</sup> .  This antibody has been characterized in the literature as containing a lambda light chain.

## Application References

(PubMed link indicates  
BioLegend citation)

1. Yokoyama WM, *et al.* 1988. *J. Immunol.* 141:369. (IP)
2. Sobel ES, *et al.* 1993. *J. Immunol.* 150:673.
3. Karlhofer FM, *et al.* 1991. *J. Immunol.* 146:3662.
4. Zhou X, *et al.* 2005. *J. Biol. Chem.* 280:31240. (IHC)
5. Podd BS, *et al.* 2006. *J. Immunol.* 176:6532. (IHC)
6. Lawson BR, *et al.* 2007. *J. Immunol.* 178:5366.
7. Lee JW, *et al.* 2006. *Nature Immunol.* 8:181.
8. Epardaud M, *et al.* 2008. *Cancer Res.* 15:2972. [PubMed](#)
9. Jordan JM, *et al.* 2008. 76:3717. [PubMed](#)
10. Kenna TJ, *et al.* 2008. *Blood* 111:2091. [PubMed](#)
11. Ishikawa C, *et al.* 2013. *Biochim Biophys Acta.* 167:99. [PubMed](#)

## RRID

AB\_3674965 (BioLegend Cat. No. 104561)  
AB\_3674965 (BioLegend Cat. No. 104562)

## Antigen Details

---

<b>Structure</b>	C-type lectin, 27/33 kD
<b>Distribution</b>	Activated T cells and B cells, NK cells, granulocytes, thymocytes, platelets
<b>Function</b>	Lymphocyte activation
<b>Cell Type</b>	B cells, Granulocytes, NK cells, Platelets, T cells, Thymocytes, Tregs
<b>Biology Area</b>	Costimulatory Molecules, Immunology, Innate Immunity
<b>Molecular Family</b>	CD Molecules
<b>Antigen References</b>	<ol style="list-style-type: none"><li>1. Barclay AN, <i>et al.</i> 1997. <i>The Leukocyte Antigen FactsBook</i> Academic Press.</li><li>2. Testi R, <i>et al.</i> 1994. <i>Immunol. Today</i> 15:479.</li><li>3. Moretta A, <i>et al.</i> 1991. <i>J. Exp. Med.</i> 174:1393.</li><li>4. Yokoyama WM, <i>et al.</i> 1988. <i>J. Immunol.</i> 141:369.</li></ol>
<b>Gene ID</b>	<a href="#">12515</a>

## Related Protocols

---

- [Cell Surface Flow Cytometry Staining Protocol](#)


## Other Formats

---

Biotin anti-mouse CD69, FITC anti-mouse CD69, PE anti-mouse CD69, PE/Cyanine5 anti-mouse CD69, Purified anti-mouse CD69, PE/Cyanine7 anti-mouse CD69, APC anti-mouse CD69, Alexa Fluor® 488 anti-mouse CD69, Alexa Fluor® 647 anti-mouse CD69, PerCP anti-mouse CD69, PerCP/Cyanine5.5 anti-mouse CD69, Pacific Blue™ anti-mouse CD69, Brilliant Violet 421™ anti-mouse CD69, APC/Cyanine7 anti-mouse CD69, Brilliant Violet 605™ anti-mouse CD69, Brilliant Violet 510™ anti-mouse CD69, Purified anti-mouse CD69 (Maxpar® Ready), PE/Dazzle™ 594 anti-mouse CD69, Brilliant Violet 711™ anti-mouse CD69, Alexa Fluor® 700 anti-mouse CD69, Brilliant Violet 650™ anti-mouse CD69, Brilliant Violet 785™ anti-mouse CD69, TotalSeq™-A0197 anti-mouse CD69, APC/Fire™ 750 anti-mouse CD69, TotalSeq™-C0197 anti-mouse CD69, TotalSeq™-B0197 anti-mouse CD69, KIRAVIA Blue 520™ anti-mouse CD69, Spark NIR™ 685 anti-mouse CD69, Spark Red™ 718 anti-mouse CD69, Spark Blue™ 574 anti-mouse CD69 (Flexi-Fluor™), Spark PLUS B550™ anti-mouse CD69, Spark PLUS UV395™ anti-mouse CD69 Antibody, Spark YG™ 581 anti-mouse CD69 (Flexi-Fluor™), Spark YG™ 593 anti-mouse CD69 (Flexi-Fluor™) Antibody, Spark UV™ 387 anti-mouse CD69 (Flexi-Fluor™)

## Product Data

---

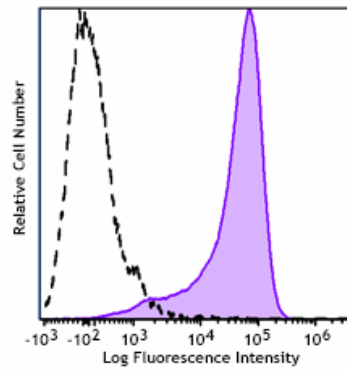

PMA+ionomycin-stimulated C57BL/6 mouse splenocytes (6 hours) were stained with anti-mouse CD69 (clone H1.2F3) Spark PLUS B550™ (filled histogram) or left unstained (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at [www.biolegend.com/terms](http://www.biolegend.com/terms) ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587