

Spark PLUS B550™ anti-mouse TIGIT (Vstm3) Antibody

Catalog# / Size	622209 / 25 µg 622210 / 100 µg
Clone	A17200C
Regulatory Status	RUO
Other Names	T-Cell Immunoreceptor with Ig and ITIM domains, T cell Ig, Immunoglobulin Domain-Containing Protein 9, V-Set, Transmembrane Domain-Containing Protein 3, WUCAM
Isotype	Rat IgG2a, κ
Description	<p>TIGIT (T cell immunoreceptor with Ig and ITIM domains precursor) contains an immunoglobulin variable domain, a transmembrane domain and an immunoreceptor tyrosine-based inhibitory motif. TIGIT is an inhibitory receptor expressed on the surface of natural killer (NK) cells. TIGIT recognizes nectin and nectin-like adhesion molecules and thus plays a critical role in the innate immune response to malignant transformation. TIGIT in conjunction with CD226 forms a pathway that has prominent similarities to the costimulatory CD28/CTLA-4 pathway. CD226 and TIGIT bind the same set of ligands (CD155 and CD112) and CD226 is a positive regulator of T cell responses, while TIGIT inhibits them. Knockdown of TIGIT expression in human CD4⁺ T cells results in an increase of both T-bet and IFN-γ mRNA and protein expression. These effects are overcome by blocking CD226 signaling indicating that TIGIT exerts immunosuppressive effects by competing with CD226 for the same ligand, CD155. The binding of CD155 to TIGIT on human dendritic cells enhances the production of interleukin 10 and diminishes the production of interleukin 12p40. TIGIT shows an inhibitory effect on cell-signaling pathway in murine models of experimental autoimmune encephalomyelitis (EAE). Loss of TIGIT expression in mice results in hyperproliferative T cell responses and increased susceptibility to EAE.</p>

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Rat
Immunogen	Mouse TIGIT Recombinant Protein
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Preparation	The antibody was purified by affinity chromatography and conjugated with Spark PLUS B550™ under optimal conditions.
Concentration	0.2 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is ≤ 0.5 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>* Spark PLUS B550™ has a maximum excitation of 516 nm and a maximum emission of 540 nm.</p>
Excitation Laser	Blue Laser (488 nm)
RRID	AB_3683473 (BioLegend Cat. No. 622209)

Antigen Details

Structure	Homodimer
Distribution	Human follicular B helper T cells (TFH), regulatory, memory and activated T cells, NK cells.
Function	Immunosuppressive effects by binding CD155 and modulating cytokine production by DC. Inhibits NK cell cytotoxicity. Inhibits T cell functions by competing with CD226. Tregs TIGIT ⁺ cells suppress proinflammatory Th1 and Th17 cells, but not Th2 cell responses.
Interaction	Dendritic cells, monocytes
Ligand/Receptor	Poliovirus receptor (PVR/CD155) and nectin-2 (PRR-2/CD112)
Cell Type	NK cells, T cells
Biology Area	Cell Adhesion, Cell Biology, Immunology, Signal Transduction
Molecular Family	Adhesion Molecules, Immune Checkpoint Receptors, Soluble Receptors
Antigen References	

1. Yu X, *et al.* 2009. *Nat. Immunol.* 10:48.
2. Stanietsky N, *et al.* 2009. *Proc. Natl. Acad. Sci. U S A.* 106:17858.
3. Joller N, *et al.* 2011. *J. Immunol.* 186:1338.
4. Tahara-Hanaoka S, *et al.* 2004. *Int. Immunol.* 16:533.
5. Lozano E, *et al.* 2012. *J. Immunol.* 188:3869.
6. Joller N, *et al.* 2014. *Immunity.* 40:569.
7. Deuss FA, *et al.* 2017. *J. Biol. Chem.* 292:11413.

Gene ID [100043314](#)


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Ultra-LEAF™ Purified anti-mouse TIGIT (Vstm3), PE anti-mouse TIGIT (Vstm3), PE/Cyanine5 anti-mouse TIGIT (Vstm3), Spark PLUS B550™ anti-mouse TIGIT (Vstm3) Antibody

Product Data

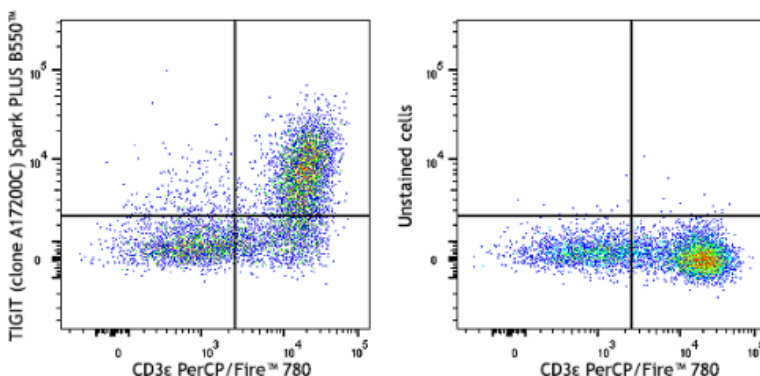

Con A and IL-2 stimulated BALB/c splenocytes (3 days) were stained with anti-mouse CD3ε (clone 145-2C11) PerCP/Fire™ 780 and anti-mouse TIGIT (Vstm3) (clone A17200C) Spark PLUS B550™ (left) or stained with anti-mouse CD3ε (clone 145-2C11) PerCP/Fire™ 780 only (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended

use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587