

Spark PLUS UV395™ anti-human CD69 Antibody

Catalog# / Size	310973 / 25 tests 310974 / 100 tests
Clone	FN50
Regulatory Status	RUO
Workshop	IV A91
Other Names	Very Early Activation Antigen (VEA), Activation inducer molecule (AIM)
Isotype	Mouse IgG1, κ
Description	CD69 is a 27-33 kD type II transmembrane protein also known as activation inducer molecule (AIM), very early activation antigen (VEA), and MLR3. It is a member of the C-type lectin family, expressed as a disulfide-linked homodimer. Other members of this receptor family include NKG2, NKR-P1 CD94, and Ly49. CD69 is transiently expressed on activated leukocytes including T cells, thymocytes, B cells, NK cells, neutrophils, and eosinophils. CD69 is constitutively expressed by a subset of medullary mature thymocytes, platelets, mantle B cells, and certain CD4 ⁺ T cells in germinal centers of normal lymph nodes. CD69 is involved in early events of lymphocyte, monocyte, and platelet activation, and has a functional role in redirected lysis mediated by activated NK cells.

Product Details

Verified Reactivity	Human
Reported Reactivity	African Green, Baboon, Chimpanzee, Cynomolgus, Pigtailed Macaque, Rhesus
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Spark PLUS UV395™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * Spark PLUS UV395™ has a maximum excitation of 355 nm and a maximum emission of 385 nm.
Excitation Laser	Ultraviolet Laser (355 nm)
Application Notes	Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen tissue sections ² , immunofluorescence microscopy ³ , and spatial biology (IBEX) ^{8,9} .

Application References

(PubMed link indicates
BioLegend citation)

1. Knapp WB, *et al.* 1989. Leucocyte Typing IV. Oxford University Press. New York.
2. Sakkas LI, *et al.* 1998. *Clin. and Diag. Lab. Immunol.* 5:430. (IHC)
3. Kim JR, *et al.* 2005. *BMC Immunol.* 6:3. (IF)
4. Verjans GM, *et al.* 2007. *P. Natl. Acad. Sci. USA* 104:3496.
5. Lu H, *et al.* 2009. *Toxicol Sci.* 112:363. (FC) [PubMed](#)
6. Thakral D, *et al.* 2008. *J. Immunol.* 180:7431. (FC) [PubMed](#)
7. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
8. Radtke AJ, *et al.* 2020. *Proc Natl Acad Sci USA.* 117:33455-33465. (SB) [PubMed](#)
9. Radtke AJ, *et al.* 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

RRID

AB_3106238 (BioLegend Cat. No. 310973)
AB_3106238 (BioLegend Cat. No. 310974)

Antigen Details

Structure	C-type lectin, type II glycoprotein, 28/32 kD
Distribution	Activated T cells, B cells, NK cells, granulocytes, thymocytes, platelets, Langerhans cells
Function	Lymphocyte, monocyte, and platelet activation, NK cell killing
Cell Type	B cells, Granulocytes, Langerhans cells, NK cells, Platelets, T cells, Thymocytes, Tregs
Biology Area	Costimulatory Molecules, Immunology
Molecular Family	CD Molecules
Antigen References	1. Schlossman S, <i>et al.</i> Eds. 1995. Leucocyte Typing V. Oxford University Press. New York. 2. Testi R, <i>et al.</i> 1994. <i>Immunol. Today</i> 15:479.
Gene ID	969


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Purified anti-human CD69, FITC anti-human CD69, PE anti-human CD69, PE/Cyanine5 anti-human CD69, APC anti-human CD69, APC/Cyanine7 anti-human CD69, PE/Cyanine7 anti-human CD69, Alexa Fluor® 488 anti-human CD69, Alexa Fluor® 647 anti-human CD69, Pacific Blue™ anti-human CD69, Alexa Fluor® 700 anti-human CD69, Biotin anti-human CD69, PerCP/Cyanine5.5 anti-human CD69, PerCP anti-human CD69, Brilliant Violet 421™ anti-human CD69, Brilliant Violet 785™ anti-human CD69, Brilliant Violet 650™ anti-human CD69, Brilliant Violet 510™ anti-human CD69, Brilliant Violet 605™ anti-human CD69, Purified anti-human CD69 (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD69, Brilliant Violet 711™ anti-human CD69, APC/Fire™ 750 anti-human CD69, TotalSeq™-A0146 anti-human CD69, TotalSeq™-B0146 anti-human CD69, TotalSeq™-C0146 anti-human CD69, Brilliant Violet 750™ anti-human CD69, KIRAVIA Blue 520™ anti-human CD69, Spark NIR™ 685 anti-human CD69 Antibody, PE/Fire™ 640 anti-human CD69, Spark YG™ 581 anti-human CD69, TotalSeq™-D0146 anti-human CD69, Spark Blue™ 550 anti-human CD69, Spark Red™ 718 anti-human CD69, GMP PE anti-human CD69, PE/Fire™ 810 anti-human CD69, PE/Fire™ 744 anti-human CD69, Spark PLUS UV395™ anti-human CD69, Spark Blue™ 574 anti-human CD69 (Flexi-Fluor™), GMP APC anti-human CD69, Spark YG™ 593 anti-human CD69 (Flexi-Fluor™) Antibody, Spark UV™ 387 anti-human CD69 (Flexi-Fluor™), PerCP/Fire™ 806 anti-human CD69 Antibody, Brilliant Violet 570™ anti-human CD69

Product Data

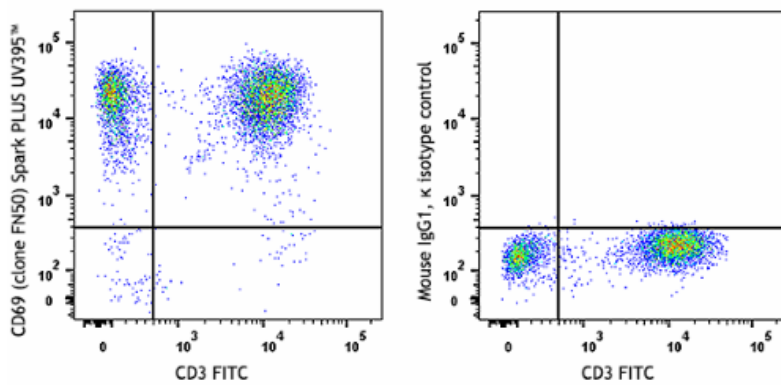

Human peripheral blood lymphocytes were stimulated with PMA/ionomycin for 6 hours and then stained with anti-human CD3 (clone UCHT1) FITC and anti-human CD69 (clone FN50) Spark PLUS UV395™ (left) or with mouse IgG1, κ Spark PLUS UV395™ isotype control.

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587