

Spark PLUS UV395™ anti-human CD90 (Thy1) Antibody

Catalog# / Size	328151 / 25 tests 328152 / 100 tests
Clone	5E10
Regulatory Status	RUO
Workshop	HCDM listed
Other Names	Thy-1, Thy1
Isotype	Mouse IgG1, κ
Description	CD90 is a 25-35 kD GPI-anchored protein, also known as Thy-1. It belongs to the Ig superfamily. Human CD90 is expressed on neuronal cells, a subset of CD34 ⁺ cells, a subset of fetal liver cells and fetal thymocytes, fibroblasts, activated endothelial cells, and some leukemia cell lines. CD34 ⁺ CD90 ⁺ cells are primitive hematopoietic stem cells. It has been reported that Thy-1 binds with β 2 and β 3 integrins and plays bimodal roles in the regulation of cell adhesion and neurite outgrowth, and inhibits hematopoietic stem cells proliferation and differentiation.

Product Details

Verified Reactivity	Human
Reported Reactivity	African Green, Baboon, Cynomolgus, Pigtailed Macaque, Rhesus, Pig
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	HEL cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Spark PLUS UV395™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 μ L per million cells in 100 μ L staining volume or 5 μ L per 100 μ L of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * Spark PLUS UV395™ has a maximum excitation of 355 nm and a maximum emission of 385 nm.
Excitation Laser	Ultraviolet Laser (355 nm)
Application Notes	Clone 5E10 recognizes an epitope on Thy-1 independent of its glycosylation, but is abolished under reducing conditions. ⁴ Additional reported (for the relevant formats) applications include: immunohistochemical staining of acetone-fixed frozen sections, immunoprecipitation ¹ , and immunofluorescence ³ .

Application References

(PubMed link indicates
BioLegend citation)

1. Craig W, *et al.* 1993. *J. Exp. Med.* 177:1331. (IP)
2. Gundlach CW 4th, *et al.* 2011. *Bioconjug. Chem.* 22:1706. (Pig Reactivity)
3. Touboul C, *et al.* 2013. *J. Transl. Med.* 11:28. (IF)
4. Bradley JE, *et al.* 2013. *Lab Invest.* 93:365. (Epitope)
5. Donnenberg VS, *et al.* 2010. *Cytometry B. Clin. Cytom.* 5:287. (IHC)

RRID

AB_3716926 (BioLegend Cat. No. 328151)
AB_3716926 (BioLegend Cat. No. 328152)

Antigen Details

Structure	25-35 kD glycoprotein, Ig superfamily
Distribution	Subset of CD34 ⁺ hematopoietic stem cells, subset of fetal thymocytes, subset of fetal liver cells, fibroblast, activated endothelial cells, neurons and some leukemia cell lines
Function	Regulate hematopoiesis and neural cell growth, cell adhesion
Ligand/Receptor	$\beta 3$ integrin, $\beta 2$ integrin
Cell Type	Endothelial cells, Fibroblasts, Hematopoietic stem and progenitors, Leukemia, Mesenchymal Stem Cells, Neurons, Thymocytes
Biology Area	Immunology, Stem Cells
Molecular Family	CD Molecules
Antigen References	1. McKenzie JL, <i>et al.</i> 1981. <i>J. Immunol.</i> 126:843. 2. Avalos AM, <i>et al.</i> 2002. <i>Biol. Res.</i> 35:231. 3. Wetzel A, <i>et al.</i> 2004. <i>J. Immunol.</i> 172:3850.
Gene ID	7070


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

TotalSeq™-D0060 anti-human CD90 (Thy1), Brilliant Violet 421™ anti-human CD90 (Thy1), Brilliant Violet 605™ anti-human CD90 (Thy1), Brilliant Violet 510™ anti-human CD90 (Thy1), Purified anti-human CD90 (Thy1) (Maxpar® Ready), Spark Blue™ 550 anti-human CD90 (Thy1) (Flexi-Fluor™), Spark PLUS UV395™ anti-human CD90 (Thy1), Spark YG™ 581 anti-human CD90 (Thy1) (Flexi-Fluor™), Spark YG™ 593 anti-human CD90 (Thy1) (Flexi-Fluor™), Spark NIR™ 685 anti-human CD90 (Thy1) (Flexi-Fluor™), APC/Cyanine7 anti-human CD90 (Thy1), PE/Dazzle™ 594 anti-human CD90 (Thy1), PE/Cyanine5 anti-human CD90 (Thy1), Alexa Fluor® 647 anti-human CD90 (Thy1), Purified anti-human CD90 (Thy1), Biotin anti-human CD90 (Thy1), FITC anti-human CD90 (Thy1), PerCP/Cyanine5.5 anti-human CD90 (Thy1), Alexa Fluor® 700 anti-human CD90 (Thy1), TotalSeq™-A0060 anti-human CD90 (Thy1), Brilliant Violet 785™ anti-human CD90 (Thy1), APC/Fire™ 750 anti-human CD90 (Thy1), Brilliant Violet 711™ anti-human CD90 (Thy1), TotalSeq™-B0060 anti-human CD90 (Thy1), PE/Cyanine7 anti-human CD90 (Thy1), Spark Blue™ 574 anti-human CD90 (Thy1) (Flexi-Fluor™), PE anti-human CD90 (Thy1), APC anti-human CD90 (Thy1), Brilliant Violet 650™ anti-human CD90 (Thy1), Spark Red™ 718 anti-human CD90 (Thy1) (Flexi-Fluor™), TotalSeq™-C0060 anti-human CD90 (Thy1)

Product Data

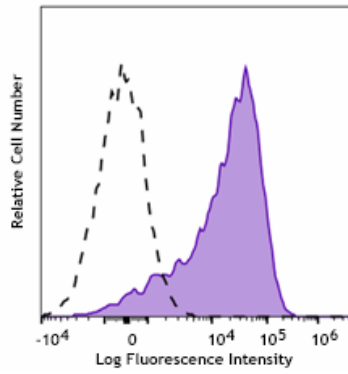

Human erythroleukemia cell line HEL was stained with anti-human CD90 (Thy1) (clone 5E10) Spark PLUS UV395™ (filled histogram) or mouse IgG1, κ Spark PLUS UV395™ isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587