

Spark PLUS UV395™ anti-human/mouse CD49f Antibody

Catalog# / Size	313647 / 25 tests 313648 / 100 tests
Clone	GoH3
Regulatory Status	RUO
Workshop	HCDM listed
Other Names	VLA-6 α chain, α_6 integrin, integrin α_6 , ITGA6
Isotype	Rat IgG2a, κ
Description	CD49f is a 120 kD integrin family member also known as VLA-6 α chain and α_6 integrin subunit. CD49f associates with either integrin β_1 (CD29) or integrin β_4 (CD104) to form receptors (VLA-6 or $\alpha_6\beta_4$ complex) for laminin and kalinin. CD49f is expressed on platelets, monocytes, T cells, placental trophoblasts, and epithelial and endothelial cells. CD49f is involved in adhesion and can act as a co-stimulatory molecule for T cell activation and proliferation.

Product Details

Verified Reactivity	Human, Mouse, Cynomolgus, Rhesus
Reported Reactivity	African Green, Baboon, Capuchin Monkey, Cat, Cow, Chimpanzee, Cynomolgus, Dog, Horse, Rabbit, Sheep, Pig
Antibody Type	Monoclonal
Host Species	Rat
Immunogen	Mouse mammary tumor cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Spark PLUS UV395™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 μ L per million cells in 100 μ L staining volume or 5 μ L per 100 μ L of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * Spark PLUS UV395™ has a maximum excitation of 355 nm and a maximum emission of 385 nm.
Excitation Laser	Ultraviolet Laser (355 nm)
Application Notes	Additional reported applications (for the relevant formats) include: immunoprecipitation ^{1,5} , <i>in vitro</i> and <i>in vivo</i> blocking of cell binding to laminin and blocking the function of integrin α_6 ^{1,4} , and immunohistochemistry of acetone-fixed frozen sections ^{2,3,5} . The GoH3 antibody has been reported to block laminin binding <i>in vitro</i> and to block integrin α_6 function <i>in vivo</i> .

Application References

(PubMed link indicates
BioLegend citation)

1. Georas SN, *et al.* 1993. *Blood* 82:2872. (IP, Block)
2. Honda T, *et al.* 1995. *J. Clin. Endocrinol. Metab.* 80:2899. (IHC)
3. Sonnenberg A, *et al.* 1986. *J. Histochem. Cytochem.* 34:1037. (IHC)
4. Nakamura K, *et al.* 1997 *Biochem. Biophys. Res. Commun.* 235:524. (Block)
5. Sonnenberg A, *et al.* 1987 *J. Biol. Chem.* 262:10376. (IP, IHC)
6. Deregibus MC, *et al.* 2007. *Blood* doi:10.1182/blood-2007-03-078709.
7. Horwitz KB, *et al.* 2008. *Proc Natl Acad Sci USA.* 105:5774. [PubMed](#)
8. Nardella C, *et al.* 2009. *Sci Signal.* 2:55. [PubMed](#)
9. Xu T, *et al.* 2010. *Mol Cancer Ther.* 9:438. [PubMed](#)
10. Stepp MA, *et al.* 2007. *J Cell Sci.* 120:2851. [PubMed](#)
11. Jo M, *et al.* 2010. *Cancer Res.* 70:8948. [PubMed](#)
12. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)

[See More](#)

RRID

AB_3716925 (BioLegend Cat. No. 313647)
AB_3716925 (BioLegend Cat. No. 313648)

Antigen Details

Structure	Integrin family, associates with $\beta 1$ or $\beta 4$, 120 kD
Distribution	Platelets, monocytes, T cells, placental trophoblasts, epithelial and endothelial cells
Function	Adhesion, receptor for laminin and kalinin; laminin binding to VLA-6 induces T cell co-stimulation for proliferation and activation
Ligand/Receptor	With integrin $\beta 1$ (CD29) forms VLA-6, with integrin $\beta 4$ (CD104) forms $\alpha \beta 4$ integrin; laminin and kalinin are ligands for these receptors
Cell Type	Embryonic Stem Cells, Endothelial cells, Epithelial cells, Monocytes, Platelets, T cells
Biology Area	Cell Adhesion, Cell Biology, Immunology, Innate Immunity, Stem Cells
Molecular Family	Adhesion Molecules, CD Molecules
Antigen References	1. Sonnenberg A, <i>et al.</i> 1990. <i>J. Cell Biol.</i> 110:2145. 2. Sonnenberg A, <i>et al.</i> 1990. <i>J. Cell. Sci.</i> 96:207. 3. Aumailley M, <i>et al.</i> 1990. <i>Exp. Cell Res.</i> 188:55. 4. Niessen CM, <i>et al.</i> 1994. <i>Exp. Cell Res.</i> 211:360.
Gene ID	16403 3655


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

TotalSeq™-D0070 anti-human/mouse CD49f, APC/Fire™ 750 anti-human/mouse CD49f, PerCP/Cyanine5.5 anti-human/mouse CD49f, Pacific Blue™ anti-human/mouse CD49f, Brilliant Violet 421™ anti-human/mouse CD49f, Brilliant Violet 510™ anti-human/mouse CD49f, Spark Blue™ 550 anti-human/mouse CD49f (Flexi-Fluor™), Spark Red™ 718 anti-human/mouse CD49f (Flexi-Fluor™), Spark PLUS UV395™ anti-human/mouse CD49f, Alexa Fluor® 488 anti-human/mouse CD49f, Alexa Fluor® 647 anti-human/mouse CD49f, PE anti-human/mouse CD49f, PE/Dazzle™ 594 anti-human/mouse CD49f, APC/Cyanine7 anti-human/mouse CD49f, TotalSeq™-C0070 anti-human/mouse CD49f, Ultra-LEAF™ Purified anti-human/mouse CD49f, TotalSeq™-B0070 anti-human/mouse CD49f, PE/Cyanine7 anti-human/mouse CD49f, Spark Blue™ 574 anti-human/mouse CD49f (Flexi-Fluor™) Antibody, TotalSeq™-A0070 anti-human/mouse CD49f, APC anti-human/mouse CD49f, Purified anti-human/mouse CD49f, Biotin anti-human/mouse CD49f, FITC anti-human/mouse CD49f

Product Data

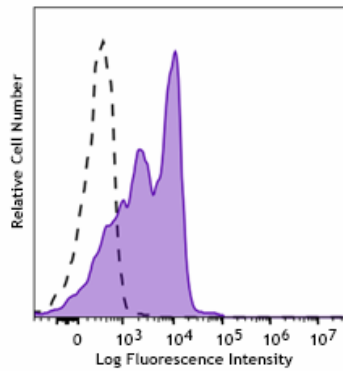

Human peripheral blood lymphocytes were stained with anti-human/mouse CD49f (clone GoH3) Spark PLUS UV395™ (filled histogram) or rat IgG2a, κ Spark PLUS UV395™ isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587