

Spark PLUS UV395™ anti-mouse CD14 Antibody

Catalog# / Size	123343 / 25 µg 123344 / 100 µg
Clone	Sa14-2
Regulatory Status	RUO
Other Names	LPS receptor
Isotype	Rat IgG2a, κ
Description	CD14 is a 53-55 kD glycosylphosphatidylinositol (GPI)-linked membrane glycoprotein also known as LPS receptor. CD14 is expressed on macrophages, dendritic cells, Kupffer cells, hepatocytes, and granulocytes. As a high-affinity receptor for LPS-LBP (LPS-binding protein) complex, CD14, in association with Toll-like Receptor 4 (TLR4) or 2 (TLR2), is involved in the clearance of gram-negative pathogens.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Rat
Immunogen	Mouse thymus or spleen
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Preparation	The antibody was purified by affinity chromatography and conjugated with Spark PLUS UV395™ under optimal conditions.
Concentration	0.2 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is ≤ 0.5 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application. * Spark PLUS UV395™ has a maximum excitation of 355 nm and a maximum emission of 385 nm.
Excitation Laser	Ultraviolet Laser (355 nm)

Antigen Details

Structure	GPI-linked membrane glycoprotein, 53-55 kD
Distribution	macrophages, dendritic cells, kupffer cells, hepatocytes
Function	LPS receptor, clearance of Gram-negative pathogens
Ligand/Receptor	LPS
Cell Type	Dendritic cells, Macrophages

Biology Area	Cell Biology, Immunology, Innate Immunity, Neuroinflammation, Neuroscience
Molecular Family	Adhesion Molecules, CD Molecules
Antigen References	1. Stocks S, <i>et al.</i> 1990. <i>Biochem. J.</i> 268:275. 2. Akashi S, <i>et al.</i> 2003. <i>J. Exp. Med.</i> 198:1035. 3. Matsuura K, <i>et al.</i> 1994. <i>J. Exp. Med.</i> 179:1671. 4. Liu S, <i>et al.</i> 1998. <i>Infect. Immun.</i> 66:5089.
Gene ID	12475


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Purified anti-mouse CD14 (Maxpar® Ready), Brilliant Violet 510™ anti-mouse CD14, Brilliant Violet 421™ anti-mouse CD14, PE/Dazzle™ 594 anti-mouse CD14, Alexa Fluor® 647 anti-mouse CD14, APC/Cyanine7 anti-mouse CD14, Purified anti-mouse CD14, FITC anti-mouse CD14, TotalSeq™-C0424 anti-mouse CD14, TotalSeq™-B0424 anti-mouse CD14, Brilliant Violet 785™ anti-mouse CD14 Antibody, Biotin anti-mouse CD14, PE anti-mouse CD14, PerCP/Cyanine5.5 anti-mouse CD14, TotalSeq™-A0424 anti-mouse CD14, PE/Cyanine7 anti-mouse CD14, Spark PLUS UV395™ anti-mouse CD14, APC anti-mouse CD14, APC/Fire™ 750 anti-mouse CD14, Brilliant Violet 605™ anti-mouse CD14

Product Data

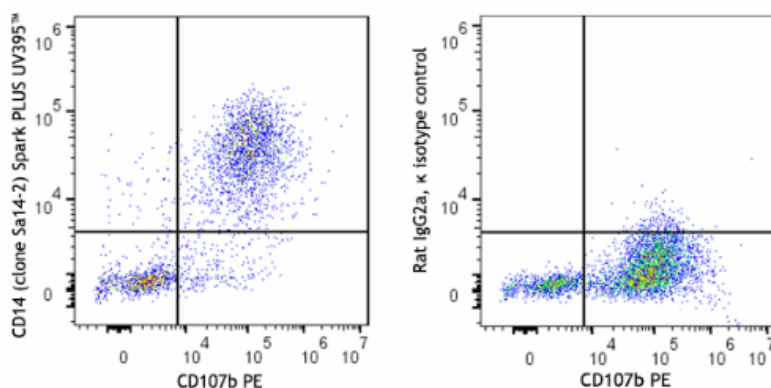

Thioglycollate-elicited BALB/c mouse peritoneal macrophages were stained with anti-mouse CD107b (clone M3/84) PE and anti-mouse CD14 (clone Sa14-2) Spark PLUS UV395™ (left) or rat IgG2a, κ Spark PLUS UV395™ isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587