

Spark PLUS UV395™ anti-mouse IFN-γ Antibody

Catalog# / Size	505871 / 25 µg 505872 / 100 µg
Clone	XMG1.2
Regulatory Status	RUO
Other Names	Interferon-γ, Immune interferon, Type II interferon, T cell interferon, Macrophage-activating factor (MAF)
Isotype	Rat IgG1, κ
Description	IFN-γ is a potent multifunctional cytokine which is secreted primarily by activated NK cells and T cells. Originally characterized based on anti-viral activities, IFN-γ also exerts anti-proliferative, immunoregulatory, and proinflammatory activities. IFN-γ can upregulate MHC class I and II antigen expression by antigen-presenting cells.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Rat
Immunogen	<i>E. coli</i> -expressed, recombinant mouse IFN-γ
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Preparation	The antibody was purified by affinity chromatography and conjugated with Spark PLUS UV395™ under optimal conditions.
Concentration	0.2 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	ICFC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by intracellular immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is ≤ 1.0 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application. * Spark PLUS UV395™ has a maximum excitation of 355 nm and a maximum emission of 385 nm.
Excitation Laser	Ultraviolet Laser (355 nm)
Application Notes	ELISA^{1-4,11,14} or ELISPOT⁵ Detection: The biotinylated XMG1.2 antibody is useful as a detection antibody for a sandwich ELISA or ELISPOT assay, when used in conjunction with purified R4-6A2 antibody (Cat. No. 505702/505706) as the capture antibody and recombinant mouse IFN-γ (Cat. No. 575309) as the standard. ELISA or ELISPOT Capture: The purified XMG1.2 antibody is useful as a capture antibody for a sandwich ELISA or ELISPOT assay, when used in conjunction with biotinylated R4-6A2 antibody (Cat. No. 505704) as the detection antibody and recombinant mouse IFN-γ (Cat. No. 575309) as the standard. The LEAF™ purified antibody is suggested for ELISPOT capture (Cat. No. 505812). Flow Cytometry^{7,8,12,13,16}: The fluorochrome-labeled XMG1.2 antibody is useful for intracellular immunofluorescent staining and flow cytometric analysis to identify IFN-γ-producing cells within mixed cell populations. Neutralization^{1-3,9,10}: The XMG1.2 antibody can neutralize the bioactivity of natural or recombinant IFN-γ. The LEAF™ purified antibody (Endotoxin <0.1 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for neutralization of mouse IFN-γ bioactivity <i>in vivo</i> and <i>in vitro</i> (Cat. No.

505812). For *in vivo* studies or highly sensitive assays, we recommend Ultra-LEAF™ purified antibody (Cat. No. 505834) with a lower endotoxin limit than standard LEAF™ purified antibodies (Endotoxin <0.01 EU/μg).

Additional reported applications (for the relevant formats) include: Western blotting, immunohistochemical staining of frozen tissue sections^{6,22,23}, and immunocytochemistry.

Note: For testing mouse IFN-γ in serum, plasma or supernatant, BioLegend's ELISA Max™ Sets (Cat. No. 430801 to 430806) are specially developed and recommended.

Application References

(PubMed link indicates BioLegend citation)

1. Abrams J, *et al.* 1992. *Immunol. Rev.* 127:5. (ELISA, Neut)
2. Sander B, *et al.* 1993. *J. Immunol. Meth.* 166:201. (ELISA, Neut)
3. Abrams J, *et al.* 1995. *Curr. Prot. Immunol.* John Wiley and Sons, New York. Unit 6.20. (ELISA, Neut)
4. Yang X, *et al.* 1993. *J. Immunoassay* 14:129. (ELISA)
5. Klinman D, *et al.* 1994. *Curr. Prot. Immunol.* John Wiley and Sons, New York. Unit 6.19. (ELISPOT)
6. Sander B, *et al.* 1991. *Immunol. Rev.* 119:65. (IHC)
7. Ferrick D, *et al.* 1995. *Nature* 373:255. (FC)
8. Ko SY, *et al.* 2005. *J. Immunol.* 175:3309. (FC) [PubMed](#)
9. Peterson KE, *et al.* 2000. *J. Virol.* 74:5363. (Neut)
10. DeKrey GK, *et al.* 1998. *Infect. Immun.* 66:827. (Neut)
11. Dzhagalov I, *et al.* 2007. *J. Immunol.* 178:2113. (ELISA)
12. Lawson BR, *et al.* 2007. *J. Immunol.* 178:5366. (FC)

[See More](#)

Antigen Details

Structure	Cytokine; dimer; 40-80 kD (Mammalian)
Bioactivity	Antiviral/antiparasitic activities; inhibits proliferation; enhances MHC class I and II expression on APCs
Cell Sources	CD8 ⁺ and CD4 ⁺ T cells, NK cells
Cell Targets	T cells, B cells, macrophages, NK cells, endothelial cells, fibroblasts
Receptors	IFN-γRα (CDw119) dimerized with IFN-γRβ (AF-1)
Cell Type	Tregs
Biology Area	Cell Biology, Immunology, Neuroinflammation, Neuroscience
Molecular Family	Cytokines/Chemokines
Antigen References	1. Fitzgerald K, <i>et al.</i> Eds. 2001. The Cytokine FactsBook. Academic Press, San Diego. 2. De Maeyer E, <i>et al.</i> 1992. <i>Curr. Opin. Immunol.</i> 4:321. 3. Farrar M, <i>et al.</i> 1993. <i>Annu. Rev. Immunol.</i> 11:571. 4. Gray P, <i>et al.</i> 1987. <i>Lymphokines</i> 13:151.
Regulation	Upregulated by IL-2, FGF-basic, EGF; downregulated by 1-α-25-Dihydroxy vitamin D3, dexamethasone
Gene ID	15978

Related Protocols


- [Intracellular Flow Cytometry Staining Protocol](#)

Other Formats

Spark NIR™ 685 anti-mouse IFN-γ, Spark UV™ 387 anti-mouse IFN-γ, Alexa Fluor® 700 anti-mouse IFN-γ, APC/Cyanine7 anti-mouse IFN-γ, Brilliant Violet 421™ anti-mouse IFN-γ, Brilliant Violet 650™ anti-mouse IFN-γ, Ultra-LEAF™ Purified anti-mouse IFN-γ, Brilliant Violet 785™ anti-mouse IFN-γ, Purified anti-mouse IFN-γ, APC anti-mouse IFN-γ, Biotin anti-mouse IFN-γ, FITC anti-mouse IFN-γ, PE/Dazzle™ 594 anti-mouse IFN-γ, Brilliant Violet 750™ anti-mouse IFN-γ, TotalSeq™-C1526 anti-mouse IFN-γ Antibody, Alexa Fluor® 488 anti-mouse IFN-γ, Alexa Fluor® 647 anti-mouse IFN-γ, Spark PLUS UV395™ anti-mouse IFN-γ, Purified anti-mouse IFN-γ (Maxpar® Ready), PE anti-mouse IFN-γ, GoInVivo™ Purified anti-mouse IFN-γ, PE/Cyanine7 anti-mouse

IFN- γ , Brilliant Violet 711™ anti-mouse IFN- γ , Brilliant Violet 605™ anti-mouse IFN- γ , Brilliant Violet 510™ anti-mouse IFN- γ , Pacific Blue™ anti-mouse IFN- γ , PerCP/Cyanine5.5 anti-mouse IFN- γ , APC/Fire™ 750 anti-mouse IFN- γ

Product Data

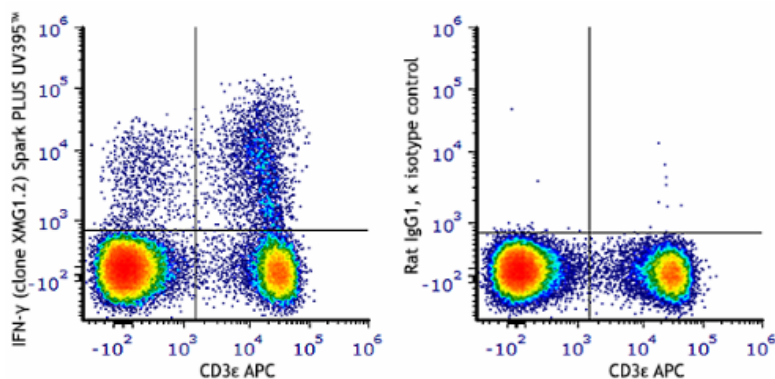

PMA+ionomycin stimulated C57BL/6 mouse splenocytes (6 hours, in the presence of monensin) were surface stained with anti-mouse CD3 ϵ (clone 145-2C11) APC. Cells were then fixed, permeabilized, and intracellularly stained with anti-mouse IFN- γ (clone XMG1.2) Spark PLUS UV395™ (left) or rat IgG1, κ Spark PLUS UV395™ isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587