

Spark PLUS V475™ anti-human CD127 (IL-7Rα) Antibody

Catalog# / Size	351387 / 25 tests
Clone	A019D5
Regulatory Status	RUO
Other Names	IL-7 receptor α chain, IL-7Rα
Isotype	Mouse IgG1, κ
Description	CD127 is a 60-90 kD type I transmembrane glycoprotein also known as IL-7 receptor α chain or IL-7Rα. It forms a heterodimer with the common γ chain (γc or CD132) which is shared with the receptors for IL-2, IL-4, IL-9, IL-13, IL-15, and IL-21. CD127 is expressed on immature B cells through early pre-B stage cells, thymocytes (except CD4/CD8 double positive thymocytes), peripheral T cells, and bone marrow stromal cells. CD127 has been reported to be a useful marker for identifying memory and effector T cells. Studies have shown that CD127 expression is down-modulated on Treg cells. It can be used as a marker for differentiation of Treg and conventional T cells. The ligation of IL-7 with its receptor is important for stimulation of mature and immature T cells as well as immature B cell proliferation and development.

Product Details

Verified Reactivity	Human
Reported Reactivity	African Green, Baboon, Cynomolgus, Rhesus
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Recombinant human CD127
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Spark PLUS V475™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our Certificate of Analysis online tool.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * Spark PLUS V475™ has a maximum excitation of 404 nm and a maximum emission of 475 nm.
Excitation Laser	Violet Laser (405 nm)
Application Notes	Additional reported (for the relevant formats) application: proteogenomics ¹ .
Application References	1. Peterson VM, <i>et al.</i> 2017. <i>Nat. Biotechnol.</i> 35:936. (PG)
(PubMed link indicates BioLegend citation)	

Antigen Details

Structure	Type I transmembrane glycoprotein, associates with CD132, 60-90 kD
Distribution	Immature B cells through early pre-B stage, thymocytes (except CD4/CD8 double positive thymocytes), peripheral T cells, bone marrow stromal cells
Function	T cell and immature B cell proliferation and development
Ligand/Receptor	IL-7
Cell Type	B cells, T cells, Thymocytes, Tregs
Biology Area	Immunology
Molecular Family	CD Molecules, Cytokine/Chemokine Receptors
Antigen References	1. Sudo T, <i>et al.</i> 1993. <i>P. Natl. Acad. Sci. USA</i> 90:9125. 2. He YW and Malek TR. 1998. <i>Crit. Rev. Immunol.</i> 18:503. 3. Huster KM, <i>et al.</i> 2004. <i>P. Natl. Acad. Sci. USA</i> 101:5610. 4. Pillai M, <i>et al.</i> 2004. <i>Leukemia Lymphoma</i> 45:2403. 5. Morrissey PJ, <i>et al.</i> 1989. <i>J. Exp. Med.</i> 169:707. 6. Liu W, <i>et al.</i> 2006. <i>J. Exp. Med.</i> 203:1701.
Gene ID	3575


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Spark YG™ 581 anti-human CD127 (IL-7Rα), Brilliant Violet 750™ anti-human CD127 (IL-7Rα), APC/Fire™ 810 anti-human CD127 (IL-7Rα) Antibody, PE/Cyanine7 anti-human CD127, GMP PE/Cyanine7 anti-human CD127 (IL-7Rα), PE/Fire™ 744 anti-human CD127 (IL-7Rα), TotalSeq™-D0390 anti-human CD127 (IL-7Rα), APC/Fire™ 750 anti-human CD127, PE anti-human CD127, Alexa Fluor® 700 anti-human CD127 (IL-7Rα), Biotin anti-human CD127 (IL-7Rα), APC/Fire™ 750 anti-human CD127 (IL-7Rα), APC/Cyanine7 anti-human CD127 (IL-7Rα), Purified anti-human CD127 (IL-7Rα), PE anti-human CD127 (IL-7Rα), Brilliant Violet 421™ anti-human CD127 (IL-7Rα), FITC anti-human CD127 (IL-7Rα), Alexa Fluor® 488 anti-human CD127 (IL-7Rα), Brilliant Violet 785™ anti-human CD127 (IL-7Rα), APC anti-human CD127 (IL-7Rα), PerCP/Cyanine5.5 anti-human CD127 (IL-7Rα), Brilliant Violet 711™ anti-human CD127 (IL-7Rα), PE/Dazzle™ 594 anti-human CD127 (IL-7Rα), TotalSeq™-B0390 anti-human CD127 (IL-7Rα), Spark NIR™ 685 anti-human CD127 (IL-7Rα), PE/Fire™ 640 anti-human CD127 (IL-7Rα), PE/Fire™ 700 anti-human CD127 (IL-7Rα) Antibody, Spark PLUS B550™ anti-human CD127 (IL-7Rα), GMP PerCP/Cyanine5.5 anti-human CD127 (IL-7Rα), GMP PE anti-human CD127 (IL-7Rα), GMP APC/Fire™ 750 anti-human CD127 (IL-7Rα), StarBright UltraViolet 740 anti-human CD127 (IL-7Rα), APC anti-human CD127, Spark PLUS UV395™ anti-human CD127 (IL-7Rα) (Flexi-Fluor™), PE/Cyanine5 anti-human CD127 (IL-7Rα), PerCP/Fire™ 780 anti-human CD127 (IL-7Rα), Spark YG™ 593 anti-human CD127 (IL-7Rα) (Flexi-Fluor™) Antibody, PerCP/Cyanine5.5 anti-human CD127, Spark Red™ 718 anti-human CD127 (IL-7Rα), PerCP/Fire™ 806 anti-human CD127 (IL-7Rα), Purified anti-human CD127 (IL-7Rα) (Maxpar® Ready), Pacific Blue™ anti-human CD127 (IL-7Rα), Alexa Fluor® 647 anti-human CD127 (IL-7Rα), PE/Cyanine7 anti-human CD127 (IL-7Rα), Brilliant Violet 570™ anti-human CD127 (IL-7Rα), Brilliant Violet 650™ anti-human CD127 (IL-7Rα), Brilliant Violet 510™ anti-human CD127 (IL-7Rα), Brilliant Violet 605™ anti-human CD127 (IL-7Rα), TotalSeq™-A0390 anti-human CD127 (IL-7Rα), TotalSeq™-C0390 anti-human CD127 (IL-7Rα), KIRAVIA Blue 520™ anti-human CD127 (IL-7Rα), Spark PLUS V475™ anti-human CD127 (IL-7Rα)

Product Data

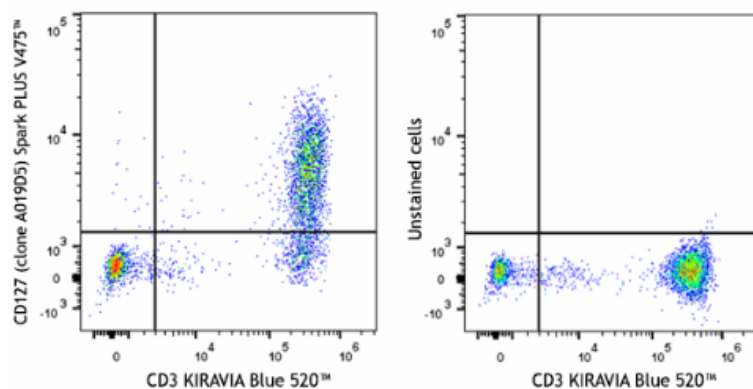

Human peripheral blood cells were stained with anti-human CD3 (clone UCHT1) KIRAVIA Blue 520™ and anti-human CD127 (IL-7Rα) (clone A019D5) Spark PLUS V475™ (left) or stained with anti-human CD3 (clone UCHT1) KIRAVIA Blue 520™ only (right). Data shown are gated on the lymphocyte population.

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587