

## Spark PLUS V475™ anti-human CD56 (NCAM) Antibody

<b>Catalog# / Size</b>	318359 / 25 tests 318360 / 100 tests
<b>Clone</b>	HCD56
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	Leu-19, NKH1
<b>Isotype</b>	Mouse IgG1, κ
<b>Description</b>	CD56 is a single transmembrane glycoprotein also known as NCAM (Neural Cell Adhesion Molecule), Leu-19, or NKH1. It is a member of the Ig superfamily. The 140 kD isoform is expressed on NK cells and NK-T cells. CD56 is also expressed in the brain (cerebellum and cortex) and at neuromuscular junctions. Certain large granular lymphocyte (LGL) leukemias, small-cell lung carcinomas, neuronal derived tumors, myelomas, and myeloid leukemias also express CD56. CD56 plays a role in homophilic and heterophilic adhesion via binding to itself or heparin sulfate.

### Product Details

---

<b>Verified Reactivity</b>	Human
<b>Reported Reactivity</b>	African Green, Baboon, Cynomolgus, Rhesus
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with Spark PLUS V475™ under optimal conditions.
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our <a href="#">Certificate of Analysis</a> online tool.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.  * Spark PLUS V475™ has a maximum excitation of 404 nm and a maximum emission of 475 nm.
<b>Excitation Laser</b>	Violet Laser (405 nm)
<b>Application Notes</b>	Clone HCD56 is not recommended for immunohistochemistry formalin-fixed paraffin-embedded tissue.
<b>Application References</b>	1. Kishimoto T, <i>et al.</i> Eds. 1997. Leucocyte Typing VI. Garland Publishing Inc. London. 2. Correia DV, <i>et al.</i> 2011. <i>Blood</i> 118:992. (FC) <a href="#">PubMed</a>
<b>(PubMed link indicates BioLegend citation)</b>	

### Antigen Details

---

<b>Structure</b>	Ig superfamily, single transmembrane or GPI-anchored glycoprotein
<b>Distribution</b>	NK cells, T subset, neural tissue, some LGL and myeloid leukemias
<b>Function</b>	Adhesion
<b>Ligand/Receptor</b>	Heparin sulfate
<b>Cell Type</b>	B cells, Leukemia, Mesenchymal Stem Cells, Neurons, NK cells, T cells
<b>Biology Area</b>	Cell Adhesion, Cell Biology, Costimulatory Molecules, Immunology, Innate Immunity, Neuroscience, Stem Cells, Synaptic Biology
<b>Molecular Family</b>	Adhesion Molecules, CD Molecules
<b>Antigen References</b>	<ol style="list-style-type: none"> <li>1. Lanier L, <i>et al.</i> 1991. <i>J. Immunol.</i> 146:4421.</li> <li>2. Hemperly J, <i>et al.</i> 1990. <i>J. Mol. Neurosci.</i> 2:71.</li> <li>3. Cremer H, <i>et al.</i> 1994. <i>Nature</i> 367:455.</li> </ol>
<b>Gene ID</b>	<a href="#">4684</a>


## Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

## Other Formats

Purified anti-human CD56 (NCAM), FITC anti-human CD56 (NCAM), PE anti-human CD56 (NCAM), PE/Cyanine5 anti-human CD56 (NCAM), APC anti-human CD56 (NCAM), Alexa Fluor® 488 anti-human CD56 (NCAM), Alexa Fluor® 647 anti-human CD56 (NCAM), Alexa Fluor® 700 anti-human CD56 (NCAM), PE/Cyanine7 anti-human CD56 (NCAM), Biotin anti-human CD56 (NCAM), PerCP/Cyanine5.5 anti-human CD56 (NCAM), Pacific Blue™ anti-human CD56 (NCAM), APC/Cyanine7 anti-human CD56 (NCAM), Brilliant Violet 421™ anti-human CD56 (NCAM), Brilliant Violet 570™ anti-human CD56 (NCAM), Brilliant Violet 605™ anti-human CD56 (NCAM), Brilliant Violet 711™ anti-human CD56 (NCAM), Brilliant Violet 510™ anti-human CD56 (NCAM), PerCP anti-human CD56 (NCAM), Brilliant Violet 650™ anti-human CD56 (NCAM), Purified anti-human CD56 (NCAM) (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD56 (NCAM), Spark Red™ 718 anti-human CD56 (NCAM), Brilliant Violet 785™ anti-human CD56 (NCAM), APC/Fire™ 810 anti-human CD56 (NCAM), Spark PLUS B550™ anti-human CD56 (NCAM), Spark PLUS V475™ anti-human CD56 (NCAM) Antibody

## Product Data

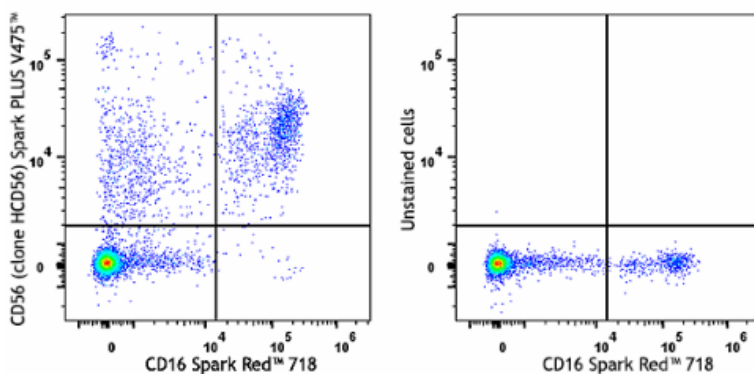

Human peripheral blood cells were stained with anti-human CD16 (clone 3G8) Spark Red™ 718 and anti-human CD56 (NCAM) (clone HCD56) Spark PLUS V475™ (left) or stained with anti-human CD16 (clone 3G8) Spark Red™ 718 only (right). Data shown are gated on the lymphocyte population.

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at [www.biolegend.com/terms](http://www.biolegend.com/terms) ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587