

Spark Red™ 718 anti-mouse CD45 Antibody

Catalog# / Size	103195 / 25 µg 103196 / 100 µg
Clone	30-F11
Regulatory Status	RUO
Other Names	T200, Ly-5, LCA
Isotype	Rat IgG2b, κ
Description	CD45 is a 180-240 kD glycoprotein also known as the leukocyte common antigen (LCA), T200, or Ly-5. It is a member of the protein tyrosine phosphatase (PTP) family, expressed on all hematopoietic cells except mature erythrocytes and platelets. There are different isoforms of CD45 that arise from variable splicing of exons 4, 5, and 6, which encode A, B, and C determinants, respectively. CD45 plays a key role in TCR and BCR signal transduction. These isoforms are very specific to the activation and maturation state of the cell as well as cell type. The primary ligands for CD45 are galectin-1, CD2, CD3, CD4, TCR, CD22, and Thy-1.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Rat
Immunogen	Mouse thymus or spleen
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Preparation	The antibody was purified by affinity chromatography and conjugated with Spark Red™ 718 under optimal conditions.
Concentration	0.5 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is ≤ 0.5 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application. * Spark Red™ 718 has a maximum excitation of 697 nm and a maximum emission of 711 nm.
Excitation Laser	Red Laser (633 nm)
Application Notes	Clone 30-F11 reacts with all isoforms and both CD45.1 and CD45.2 alloantigens of CD45. Additional reported applications (for relevant formats) include: immunoprecipitation ³ , complement-dependent cytotoxicity ^{1,5} , immunohistochemistry (acetone-fixed frozen sections, zinc-fixed paraffin-embedded sections and formalin-fixed paraffin-embedded sections) ^{4,6} , Western blotting ⁷ , and spatial biology (IBEX) ^{10,11} . The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for functional assays (Cat. No. 103163 and 103164).

Application References

(PubMed link indicates
BioLegend citation)

1. Podd BS, *et al.* 2006. *J. Immunol.* 176:6532. (FC, CMCD) [PubMed](#)
2. Haynes NM, *et al.* 2007. *J. Immunol.* 179:5099. (FC)
3. Ledbetter JA, *et al.* 1979. *Immunol. Rev.* 47:63. (IP)
4. Simon DI, *et al.* 2000. *J. Clin. Invest.* 105:293. (IHC)
5. Seaman WE. 1983. *J. Immunol.* 130:1713. (CMCD)
6. Cornet A, *et al.* 2001. *P. Natl. Acad. Sci. USA* 98:13306. (IHC)
7. Tsuboi S and Fukuda M. 1998. *J. Biol. Chem.* 273:30680. (WB) [PubMed](#)
8. Liu F, *et al.* 2012. *Blood.* 119:3295. [PubMed](#)
9. Pelletier AN, *et al.* 2012. *J. Immunol.* 188:5561. [PubMed](#)
10. Radtke AJ, *et al.* 2020. *Proc Natl Acad Sci U S A.* 117:33455-65. (SB) [PubMed](#)
11. Radtke AJ, *et al.* 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

RRID AB_3683130 (BioLegend Cat. No. 103195)
AB_3683130 (BioLegend Cat. No. 103196)

Antigen Details

Structure	Protein tyrosine phosphatase (PTP) family, 180-240 kD
Distribution	All hematopoietic cells except mature erythrocytes and platelets
Function	Phosphatase, T and B cell activation
Ligand/Receptor	Galectin-1, CD2, CD3, CD4, TCR, CD22, Thy-1
Cell Type	B cells, Dendritic cells, Mesenchymal Stem Cells, Tregs
Biology Area	Cell Biology, Immunology, Inhibitory Molecules, Innate Immunity, Neuroscience, Neuroscience Cell Markers, Stem Cells
Molecular Family	CD Molecules
Antigen References	<ol style="list-style-type: none">1. Barclay A, <i>et al.</i> 1997. <i>The Leukocyte Antigen FactsBook</i> Academic Press.2. Trowbridge IS, <i>et al.</i> 1993. <i>Annu. Rev. Immunol.</i> 12:85.3. Kishihara K, <i>et al.</i> 1993. <i>Cell</i> 74:143.4. Pulido R, <i>et al.</i> 1988. <i>J. Immunol.</i> 140:3851.
Gene ID	19264


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-mouse CD45, Biotin anti-mouse CD45, FITC anti-mouse CD45, PE anti-mouse CD45, PE/Cyanine5 anti-mouse CD45, Purified anti-mouse CD45, PE/Cyanine7 anti-mouse CD45, APC/Cyanine7 anti-mouse CD45, Alexa Fluor® 488 anti-mouse CD45, Alexa Fluor® 647 anti-mouse CD45, Pacific Blue™ anti-mouse CD45, Alexa Fluor® 700 anti-mouse CD45, PerCP/Cyanine5.5 anti-mouse CD45, PerCP anti-mouse CD45, Alexa Fluor® 594 anti-mouse CD45, Brilliant Violet 421™ anti-mouse CD45, Brilliant Violet 570™ anti-mouse CD45, Brilliant Violet 510™ anti-mouse CD45, Brilliant Violet 605™ anti-mouse CD45, Purified anti-mouse CD45 (Maxpar® Ready), PE/Dazzle™ 594 anti-mouse CD45, Brilliant Violet 711™ anti-mouse CD45, Brilliant Violet 785™ anti-mouse CD45, Brilliant Violet 650™ anti-mouse CD45, APC/Fire™ 750 anti-mouse CD45, Brilliant Violet 750™ anti-mouse CD45, TotalSeq™-A0096 anti-mouse CD45, TotalSeq™-B0096 anti-mouse CD45, Ultra-LEAF™ Purified anti-mouse CD45, Spark Blue™ 550 anti-mouse CD45, Spark NIR™ 685 anti-mouse CD45, TotalSeq™-C0096 anti-mouse CD45, Spark YG™ 570 anti-mouse CD45, PE/Fire™ 640 anti-mouse CD45, APC/Fire™ 810 anti-mouse CD45, PE/Fire™ 700 anti-mouse CD45, Spark Violet™ 538 anti-mouse CD45, Spark YG™ 593 anti-mouse CD45, Spark Blue™ 574 anti-mouse CD45 Antibody, Spark Blue™ 515 anti-mouse CD45, Spark UV™ 387 anti-mouse CD45, PE/Fire™ 810 anti-mouse CD45, Spark PLUS UV395™ anti-mouse CD45, Spark Red™ 718 anti-mouse CD45, Spark YG™ 581 anti-mouse CD45 (Flexi-Fluor™), Spark PLUS V475™ anti-mouse CD45 Antibody

Product Data

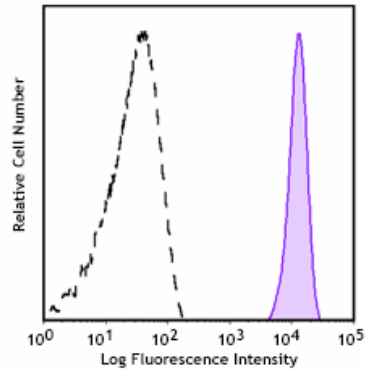

C57BL/6 mouse splenocytes were stained with anti-mouse CD45 (clone 30-F11) Spark Red™ 718 (filled histogram) or rat IgG2b, κ Spark Red™ 718 isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587