

Spark Red™ 718 anti-mouse CD80 Antibody

Catalog# / Size	104765 / 25 µg 104766 / 100 µg
Clone	16-10A1
Regulatory Status	RUO
Other Names	B7-1, B7, Ly-53
Isotype	Armenian Hamster IgG
Description	CD80 is a 60 kD highly glycosylated protein. It is a member of the Ig superfamily and is also known as B7-1, B7, and Ly-53. CD80 is constitutively expressed on dendritic cells and monocytes/macrophages, and inducibly expressed on activated B and T cells. The ligation of CD28 on T cells with CD80 and CD86 (B7-2) on antigen presenting cells (such as dendritic cells, macrophages, and B cells) elicits co-stimulation of T cells resulting in enhanced cell activation, proliferation, and cytokine production. CD80 appears to be expressed later in the immune response than CD86. CD80 can also bind to CD152, also known as CTLA-4, to deliver an inhibitory signal to T cells.

Product Details

Verified Reactivity	Mouse
Reported Reactivity	Dog
Antibody Type	Monoclonal
Host Species	Armenian Hamster
Immunogen	CHO cell line transfected with mouse B7 (CD80)
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Preparation	The antibody was purified by affinity chromatography and conjugated with Spark Red™ 718 under optimal conditions.
Concentration	0.5 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is ≤ 0.5 µg per million cells in 100 µL volume. It is recommended that the reagent be titrated for optimal performance for each application. * Spark Red™ 718 has a maximum excitation of 697 nm and a maximum emission of 711 nm.
Excitation Laser	Red Laser (633 nm)
Application Notes	Additional reported applications (for the relevant formats) include: immunoprecipitation ² , <i>in vitro</i> and <i>in vivo</i> blocking of CTLA-4 Ig to CD80 by blocking costimulation of T cells by activated B cells ²⁻⁴ , and immunohistochemical staining of acetone-fixed frozen sections ^{1,4} . The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for functional assays (Cat. Nos. 104747-104752).

Application References

(PubMed link indicates
BioLegend citation)

1. Harlan DM, *et al.* 1994. *P. Natl. Acad. Sci. USA* 91:3137. (IHC)
2. Razi-Wolf Z, *et al.* 1992. *P. Natl. Acad. Sci. USA* 89:4210. (Block, IP)
3. Hathcock KS, *et al.* 1994. *J. Exp. Med.* 180:631. (Block)
4. Herold KC, *et al.* 1997. *J. Immunol.* 158:984. (Block, IHC)
5. Ma XT, *et al.* 2006. *Cancer Res.* 66:1169.
6. Andoniou CE, *et al.* 2005. *Nature Immunology* 6:1011. (FC)
7. Lawson BR, *et al.* 2007. *J. Immunol.* 178:5366.
8. Turnquist HR, *et al.* 2007. *J. Immunol.* 178:7018.
9. Misra RS, *et al.* 2010. *J. Exp. Med.* 207:1775. [PubMed](#)
10. del Rio ML, *et al.* 2011. *Transpl. Int.* 24:501. (FC) [PubMed](#)
11. Philipsen L, *et al.* 2013. *Mol Cell Proteomics.* 12:2551. [PubMed](#)

Antigen Details

Structure	Ig superfamily, 60 kD
Distribution	Macrophages, activated B cells and T cells, dendritic cells
Function	T cell costimulation
Ligand/Receptor	CD28 (stimulatory), CD152(CTLA4) (inhibitory)
Cell Type	B cells, Dendritic cells, Macrophages, T cells, Tregs
Biology Area	Cell Biology, Costimulatory Molecules, Immunology, Neuroscience, Neuroscience Cell Markers
Molecular Family	CD Molecules, Immune Checkpoint Receptors
Antigen References	1. Barclay AN, <i>et al.</i> 1997. The Leukocyte Antigen FactsBook Academic Press. 2. Linsley PS, <i>et al.</i> 1991. <i>J. Exp. Med.</i> 174:561. 3. Salomon B, <i>et al.</i> 2001. <i>Annu. Rev. Immunol.</i> 19:225.
Gene ID	12519


Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

PE/Fire™ 640 anti-mouse CD80, Spark NIR™ 685 anti-mouse CD80, PE/Cyanine5 anti-mouse CD80, Spark Red™ 718 anti-mouse CD80 (Flexi-Fluor™), PE/Dazzle™ 594 anti-mouse CD80, Purified anti-mouse CD80 (Maxpar® Ready), APC/Fire™ 750 anti-mouse CD80, Pacific Blue™ anti-mouse CD80, Brilliant Violet 650™ anti-mouse CD80, Brilliant Violet 421™ anti-mouse CD80, Brilliant Violet 605™ anti-mouse CD80, Spark Blue™ 574 anti-mouse CD80 (Flexi-Fluor™), APC anti-mouse CD80, Alexa Fluor® 488 anti-mouse CD80, Alexa Fluor® 647 anti-mouse CD80, Biotin anti-mouse CD80, Brilliant Violet 510™ anti-mouse CD80, TotalSeq™-C0849 anti-mouse CD80, Spark PLUS B550™ anti-mouse CD80 Antibody, Spark Red™ 718 anti-mouse CD80, FITC anti-mouse CD80, Brilliant Violet 711™ anti-mouse CD80, TotalSeq™-A0849 anti-mouse CD80, Ultra-LEAF™ Purified anti-mouse CD80, Alexa Fluor® 594 anti-mouse CD80, PE/Cyanine7 anti-mouse CD80, PerCP/Cyanine5.5 anti-mouse CD80, PE anti-mouse CD80, Purified anti-mouse CD80, TotalSeq™-B0849 anti-mouse CD80 Antibody

Product Data

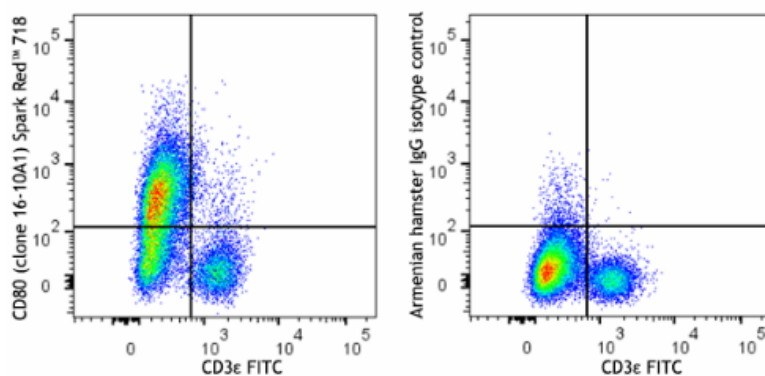

LPS stimulated (3 days) C57BL/6 mouse splenocytes were stained with anti-mouse CD3ε (clone 145-2C11) FITC and anti-mouse CD80 (clone 16-10A1) Spark Red™ 718 (left) or Armenian hamster IgG Spark Red™ 718 isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587