

## Spark Violet™ 423 anti-human CD10 Antibody

<b>Catalog# / Size</b>	312245 / 25 tests
<b>Clone</b>	HI10a
<b>Regulatory Status</b>	RUO
<b>Workshop</b>	V CD10.7
<b>Other Names</b>	Common acute lymphoblastic leukemia antigen (CALLA), Enkephalinase, Neutral endopeptidase, Neprilysin
<b>Isotype</b>	Mouse IgG1, κ
<b>Description</b>	CD10 is a 100 kD neutral endopeptidase and a member of the metalloprotease family. It is a type II transmembrane protein also known as common acute lymphoblastic leukemia antigen (CALLA), enkephalinase, and neprilysin. CD10 is expressed on B cell precursors, T cell precursors, and neutrophils. CD10 is involved in B cell development and has been shown to bind opioid enkephalins, bradykinin, angiotensins I and II, and other biologically active peptides.

### Product Details

---

<b>Verified Reactivity</b>	Human, Cynomolgus, Rhesus
<b>Reported Reactivity</b>	African Green, Baboon, Capuchin monkey, Chimpanzee
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with Spark Violet™ 423 under optimal conditions.
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our <a href="#">Certificate of Analysis</a> online tool.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.  * Spark Violet™ 423 has a maximum excitation of 400 nm and a maximum emission of 415 nm.
<b>Excitation Laser</b>	Violet Laser (405 nm)
<b>RRID</b>	AB_3662349 (BioLegend Cat. No. 312245)

### Antigen Details

---

<b>Structure</b>	Type II transmembrane metalloprotease family, 100 kD
<b>Distribution</b>	B cell precursors, T cell precursors, neutrophils
<b>Function</b>	Zinc-binding metalloproteinase, B cell development

<b>Ligand/Receptor</b>	Biologically active peptides including opioid enkephalins, bradykinin, angiotensins I & II
<b>Cell Type</b>	B cells, Neutrophils
<b>Biology Area</b>	Immunology
<b>Molecular Family</b>	CD Molecules
<b>Antigen References</b>	1. Shipp M, <i>et al.</i> 1993. <i>Blood</i> 82:1052. 2. Lu B, <i>et al.</i> 1995. <i>J. Exp. Med.</i> 181:2271.
<b>Gene ID</b>	<a href="#">4311</a>


## Related Protocols

- [Cell Surface Flow Cytometry Staining Protocol](#)

## Other Formats

Purified anti-human CD10, PE anti-human CD10, PE/Cyanine5 anti-human CD10, FITC anti-human CD10, APC anti-human CD10, APC/Cyanine7 anti-human CD10, PE/Cyanine7 anti-human CD10, PerCP/Cyanine5.5 anti-human CD10, Brilliant Violet 421™ anti-human CD10, Brilliant Violet 510™ anti-human CD10, Brilliant Violet 605™ anti-human CD10, Purified anti-human CD10 (Maxpar® Ready), Brilliant Violet 711™ anti-human CD10, PE/Dazzle™ 594 anti-human CD10, APC/Fire™ 750 anti-human CD10, TotalSeq™-A0062 anti-human CD10, TotalSeq™-C0062 anti-human CD10, TotalSeq™-B0062 anti-human CD10, Brilliant Violet 785™ anti-human CD10, TotalSeq™-D0062 anti-human CD10, GMP APC anti-human CD10, GMP FITC anti-human CD10, GMP PE/Cyanine7 anti-human CD10, GMP PE anti-human CD10, GMP PerCP/Cyanine5.5 anti-human CD10, Brilliant Violet 650™ anti-human CD10, GMP PE/Dazzle™ 594 anti-human CD10, Brilliant Violet 750™ anti-human CD10, Spark Violet™ 423 anti-human CD10, Spark YG™ 593 anti-human CD10, Spark Blue™ 574 anti-human CD10, PE/Fire™ 640 anti-human CD10, Pacific Blue™ anti-human CD10 Antibody

## Product Data

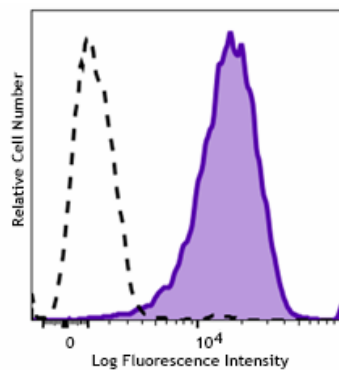

Human peripheral blood granulocytes were stained with anti-human CD10 (clone HI10a) Spark Violet™ 423 (filled histogram) or mouse IgG1, κ Spark Violet™ 423 isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at [www.biolegend.com/terms](http://www.biolegend.com/terms) ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587