

## Spark Violet™ 423 anti-human CD8 Antibody

<b>Catalog# / Size</b>	344797 / 25 tests 344798 / 100 tests
<b>Clone</b>	SK1
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	T8, Leu2
<b>Isotype</b>	Mouse IgG1, $\kappa$
<b>Description</b>	CD8a is a 32-34 kD type I glycoprotein. It forms a homodimer (CD8a/a) or heterodimer (CD8a/b) with CD8b. CD8, also known as T8 and Leu2, is a member of the immunoglobulin superfamily found on the majority of thymocytes, a subset of peripheral blood T cells, and NK cells (which express almost exclusively CD8a homodimers). CD8 acts as a co-receptor with MHC class I-restricted T cell receptors in antigen recognition and T cell activation and has been shown to play a role in thymic differentiation. Two domains in CD8a are important for function: the extracellular IgSF domain binds the $\alpha_3$ domain of MHC class I and the cytoplasmic CXCP motif binds the tyrosine kinase p56 Lck.

### Product Details

---

<b>Verified Reactivity</b>	Human, Cynomolgus, Rhesus
<b>Reported Reactivity</b>	African Green, Chimpanzee, Pigtailed Macaque, Sooty Mangabey
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with Spark Violet™ 423 under optimal conditions.
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our <a href="#">Certificate of Analysis</a> online tool.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 $\mu$ L per million cells in 100 $\mu$ L staining volume or 5 $\mu$ L per 100 $\mu$ L of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.  * Spark Violet™ 423 has a maximum excitation of 400 nm and a maximum emission of 415 nm.
<b>Excitation Laser</b>	Violet Laser (405 nm)
<b>Application Notes</b>	Clone SK1 recognizes the $\alpha$ chain of CD8. Additional reported applications (for the relevant formats) include: proteogenomics <sup>8</sup> , immunohistochemistry of acetone-fixed frozen tissue sections, and spatial biology (IBEX) <sup>9,10</sup> . This clone was tested in-house and does not demonstrate utility for formalin-fixed paraffin-embedded (FFPE) human tonsil sections.

## Application References

(PubMed link indicates  
BioLegend citation)

1. Ledbetter JA, *et al.* 1981. *J. Exp. Med.* 153:310.
2. Campanelli R, *et al.* 2002. *Intl. Immunol.* 14:39.
3. Evans RL, *et al.* 1981. *Immunol.* 78:544.
4. Wooldridge L, *et al.* 2005. *J. Bio. Chem.* 280:27491.
5. Ch'el IL, *et al.* 2011. *J Exp Med.* 208:633. [PubMed](#)
6. Carbone A, *et al.* 1999. *Blood* 93:2319. (IHC-F)
7. Ahmed A, *et al.* 2001. *J. Pathol.* 193:383. (IHC)
8. Peterson VM, *et al.* 2017. *Nat. Biotechnol.* 35:936. (PG)
9. Radtke AJ, *et al.* 2020. *Proc Natl Acad Sci USA.* 117:33455-33465. (SB) [PubMed](#)
10. Radtke AJ, *et al.* 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

## RRID

AB\_3683368 (BioLegend Cat. No. 344797)  
AB\_3683368 (BioLegend Cat. No. 344798)

## Antigen Details

---

<b>Structure</b>	Ig superfamily, homodimer or heterodimer with CD8b, 32-34 kD
<b>Distribution</b>	Majority of thymocytes, T cell subset, NK cells
<b>Function</b>	MHC class I co-receptor, thymic differentiation, T cell activation
<b>Ligand/Receptor</b>	MHC Class I molecules
<b>Cell Type</b>	NK cells, T cells, Thymocytes
<b>Biology Area</b>	Immunology
<b>Molecular Family</b>	CD Molecules
<b>Antigen References</b>	1. Barclay N, <i>et al.</i> 1993. <i>The Leucocyte Antigen FactsBook.</i> Academic Press Inc. San Diego.
<b>Gene ID</b>	<a href="#">925</a>

## Related Protocols

---

- [Cell Surface Flow Cytometry Staining Protocol](#)


## Other Formats

---

Alexa Fluor® 647 anti-human CD8, Brilliant Violet 650™ anti-human CD8, Purified anti-human CD8, FITC anti-human CD8, PE anti-human CD8, PerCP anti-human CD8, PerCP/Cyanine5.5 anti-human CD8, PE/Cyanine7 anti-human CD8, APC/Cyanine7 anti-human CD8, Alexa Fluor® 488 anti-human CD8, Pacific Blue™ anti-human CD8, Biotin anti-human CD8, APC anti-human CD8, Alexa Fluor® 700 anti-human CD8, Purified anti-human CD8 (Maxpar® Ready), Brilliant Violet 510™ anti-human CD8, Brilliant Violet 711™ anti-human CD8, Brilliant Violet 785™ anti-human CD8, Brilliant Violet 605™ anti-human CD8, PE/Dazzle™ 594 anti-human CD8, APC/Fire™ 750 anti-human CD8, Brilliant Violet 421™ anti-human CD8, TotalSeq™-A0046 anti-human CD8, TotalSeq™-C0046 anti-human CD8, Brilliant Violet 750™ anti-human CD8, TotalSeq™-B0046 anti-human CD8, Spark Blue™ 550 anti-human CD8, APC/Fire™ 810 anti-human CD8, PE/Fire™ 640 anti-human CD8, PE/Fire™ 700 anti-human CD8, TotalSeq™-D0046 anti-human CD8, GMP APC anti-human CD8, PE/Cyanine5 anti-human CD8 Antibody, Spark UV™ 387 anti-human CD8, GMP PE anti-human CD8, GMP PE/Cyanine7 anti-human CD8, Spark NIR™ 685 anti-human CD8, KIRAVIA Blue 520™ anti-human CD8, GMP FITC anti-human CD8, GMP Pacific Blue™ anti-human CD8, GMP PerCP anti-human CD8, Spark Violet™ 500 anti-human CD8, GMP APC/Fire™ 750 anti-human CD8, GMP PerCP/Cyanine 5.5 anti-human CD8, Alexa Fluor® 660 anti-human CD8a, Spark Blue™ 515 anti-human CD8, Spark Blue™ 574 anti-human CD8, Spark Violet™ 538 anti-human CD8, PE/Fire™ 810 anti-human CD8, Spark YG™ 593 anti-human CD8, Brilliant Violet 570™ anti-human CD8, PerCP/Fire™ 806 anti-human CD8, PerCP/Fire™ 780 anti-human CD8, Spark Violet™ 423 anti-human CD8 Antibody, StarBright UltraViolet 795 anti-human CD8, StarBright UltraViolet 575 anti-human CD8

## Product Data

---

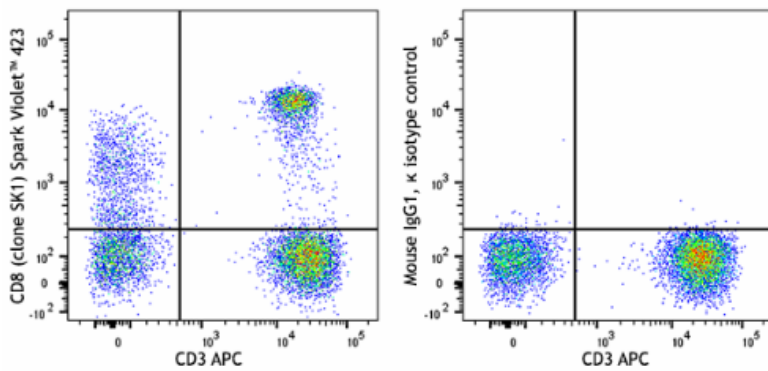

Human peripheral blood lymphocytes were stained with anti-human CD3 (clone UCHT1) APC and anti-human CD8 (clone SK1) Spark Violet™ 423 (left) or mouse IgG1, κ Spark Violet™ 423 isotype control (right).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at [www.biolegend.com/terms](http://www.biolegend.com/terms) ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587